


---

# CAS 740

## Vital Signs Monitor

---

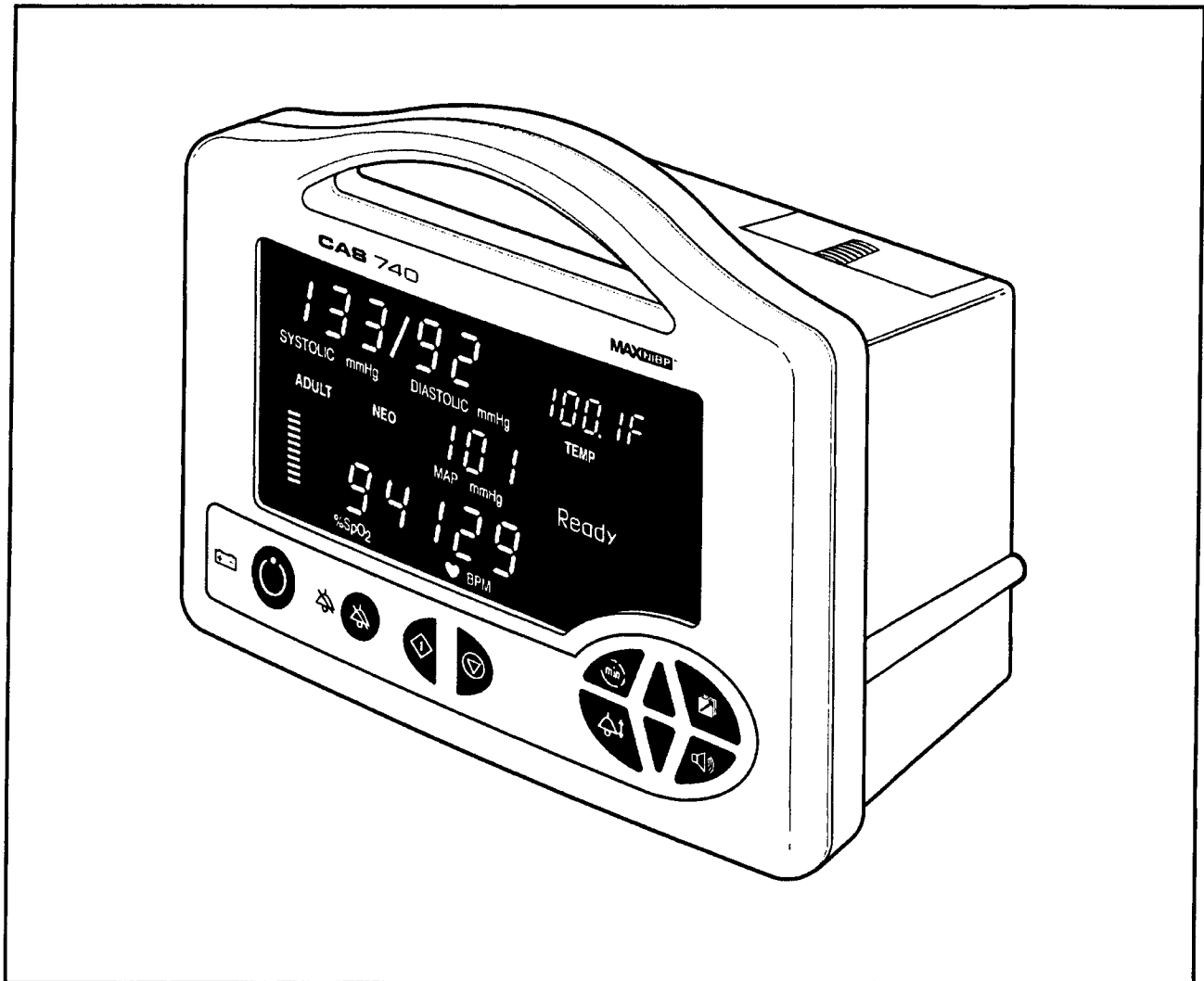

User's Manual

 **CAS** MEDICAL SYSTEMS, INC.  
TECHNOLOGY APPLIED TO MEDICINE

CE 0086





In the U.S. the following Caution applies:

**CAUTION:**  
**Federal law restricts this device to sale by or on  
the order of a physician or properly licensed  
practitioner.**

First Printing: 03/2003

Revised: 10/2004

Revised: 10/2005


## INITIAL SETUP

Before using the monitor for the first time, the following items should be performed:

- Select the operating Language
- Select the Temperature scale (if installed)
- Set the monitor's Date and Time

These items are found in the monitor's Configuration Menu.

To enter the Configuration Menu, depress and hold the ALARM LIMITS and AUDIO push button keys while turning "ON" the unit.

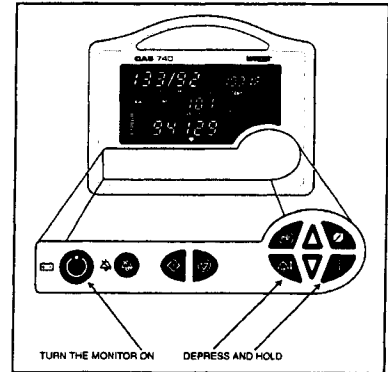

### **NOTE:**

Unit will beep once and software version will be displayed.  
Steps must be completed within 60 seconds or process must be restarted.

### Language

Press until the menu for language is displayed. Press either or to make your selection.

### Temperature Scale (if installed)

Press until the menu for temperature scale is displayed. Press either or to make your selection.

### Date

Press until the date is displayed. The day parameter is flashing. Press either or to change the Day.

Press one time. Press either or to change the Month.

Press one time. Press either or to change the Year.

### Time

Press one time. Press either or to change the Hour.

Press one time. Press either or to change the Minute.

Press one time to exit and save your changes.



## Manufacturers Declaration of Conformity

### Electronic Emissions and Immunity

The Model 740 Monitor is intended for use in the electromagnetic environment specified below. The customer or the user of the Model 740 Monitor should assure it is used in such an environment.


Emissions Test	Compliance	Electromagnetic Environment
RF emissions – CISPR 11	Group 1	The Model 740 Monitor uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions – CISPR 11	Class B	The Model 740 Monitor is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Harmonic emissions IEC 61000-3-2	Class B	
Voltage fluctuations / flicker emissions	Complies	

Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment Guidance
Electrostatic discharge (ESD) IEC 61000-4-2	+/-6 kV contact +/-8 kV air	+/-6 kV contact +/-8 kV air	Floors should be wood concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient/burst IEC 61000-4-4	+/-2 kV for power supply lines +/-1 kV for input/output lines	+/-2 kV for power supply lines +/-1 kV for input/output lines	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	+/-1 kV differential mode +/-2 kV common mode	+/-1 kV differential mode +/-2 kV common mode	Mains power quality should be that of a typical commercial or hospital environment.
Voltage Dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	< 5% $U_T$ (>95% dip in $U_T$ ) for 0.5 cycle. 40% $U_T$ (60% dip in $U_T$ ) for 5 cycles. 70% $U_T$ (30% dip in $U_T$ ) for 25 cycles. < 5% $U_T$ (> 95% dip in $U_T$ ) for 5 seconds.	< 5% $U_T$ (>95% dip in $U_T$ ) for 0.5 cycle. 40% $U_T$ (60% dip in $U_T$ ) for 5 cycles. 70% $U_T$ (30% dip in $U_T$ ) for 25 cycles. < 5% $U_T$ (> 95% dip in $U_T$ ) for 5 seconds.	Mains power quality should be that of a typical commercial or hospital environment. If user of the Model 740 Monitor requires continued operation during power mains interruptions, it is recommended that the Model 740 Monitor be powered from an uninterruptible power supply or a battery.
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.

NOTE:  $U_T$  is the A.C. mains voltage prior to application of the test level.

## Guidance and Manufacturer's Declaration – Electromagnetic Immunity

The Model 740 Monitor is intended for use in the electromagnetic environment specified below. The customer or the user of the Model 740 Monitor should insure that it is used in such an environment.

Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment - Guidance
Conducted RF IEC 61000-4-6	3 Vrms 150 kHz to 80 MHz	3 Vrms	Portable and mobile RF communications equipment should be used no closer to any part of the Model 740 Monitor, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. <b>Recommended separation distance:</b>  $d = 1.2\sqrt{P}$  $d = 1.2\sqrt{P}$ 80 MHz to 800 MHz $d = 2.3\sqrt{P}$ 800 MHz to 2.5 GHz  Where $P$ is the maximum output power rating of the transmitter in watts according to the transmitter manufacturer and $d$ is the recommended separation distance in meters.  Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey <sup>a</sup> , should be less than the compliance level in each frequency range. <sup>b</sup>  Interference may occur in the vicinity of equipment marked with the following symbol:  
Radiated RF IEC 61000-4-3	3 V/m 80 MHz to 2.5 GHz	3 V/m	

NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies.

NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is effected by absorption and reflection from structures, objects and people.

<sup>a</sup> Field strengths from fixed transmitters, such as base stations for radio (cellular / cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the Model 740 Monitor is used exceeds the applicable RF compliance level above, the Model 740 Monitor should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the Model 740 Monitor.

<sup>b</sup> Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m.

### Recommended Separation Distances Between Portable and Mobile RF Communications Equipment and the Model 740 Monitor

The Model 740 Monitor is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the Model 740 Monitor can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the Model 740 Monitor as recommended below, according to the maximum output power of the communications equipment.

Rated maximum output power of transmitter (Watts)	Separation distance according to frequency of transmitter (Meters)		
	150 kHz to 80 MHz	80 MHz to 800 MHz	800 MHz to 2.5 GHz
	$d = 1.2\sqrt{P}$	$d = 1.2\sqrt{P}$	$d = 2.3\sqrt{P}$
0.01	0.12	0.12	0.23
0.1	0.38	0.38	0.73
1	1.2	1.2	2.3
10	3.8	3.8	7.3
100	12	12	23

For transmitters operating at a maximum output power not listed above, the recommended separation distance  $d$  in meters can be estimated using the equation applicable to the frequency of the transmitter, where  $P$  is the maximum output power rating of the transmitter in watts according to the transmitter manufacturer.

NOTE 1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.



## **WARRANTY POLICY**

### **MONITOR (CAS 740)**

CAS Medical Systems, Inc. warrants the monitor, when new, to be free from defects in material and workmanship and to perform in accordance with manufacturer's specifications for a period of two (2) years from the date of original purchase from CAS or its authorized distributors or agents except as noted below.

The same warranty conditions are made for a period of one (1) year with respect to printers and battery and ninety (90) days on non-disposable accessories and certain components consisting of reusable SpO<sub>2</sub> sensors, reusable temperature probes and other accessories provided by CAS as part of the original purchase. CAS warrants blood pressure cuffs and disposable or single-patient-use products for out-of-box failure only. Where the accessory is not a CAS manufactured product, the manufacturers own warranty conditions apply.

CAS reserves the right to perform warranty service operations in its own factory, at an authorized repair facility, or at the customers' site.

Our obligation under this warranty is limited to repairing or, at our option, replacing any defective parts or our equipment, without charge, if such defects occur in normal service and with prompt notification.

Damage to any part through misuse, neglect, or accident, or by affixing any accessories or attachments other than CAS, Masimo<sup>®</sup>, Nellcor<sup>®</sup>, Nonin<sup>®</sup>, and Welch Allyn<sup>®</sup> manufactured accessories or attachments, is not covered by this warranty.

### **ACCESSORIES, BATTERIES, CUFFS, AND CERTAIN COMPONENTS**

In all cases, policy applies from date of purchase from CAS or its authorized distributors or agents.

Accessories:	(90) Days - Masimo, Nellcor and Nonin Sensors, Welch Allyn Temperature Probes.
Batteries:	(1) Year
Chargers:	(1) Year (not including power cord: see other accessories).
Cuffs (all):	Out-of-box failure only.
Other Accessories:	Out-of-box failure only.
Certain Components:	(1) Year - Printer mechanism, but not including Thermal Print Heads.
Print Heads:	Out-of-box failure only.

**THERE ARE NO WARRANTIES, WHICH EXTEND BEYOND THOSE EXPRESSLY DESCRIBED IN THIS AGREEMENT AND THE COMPANY MAKES NO WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE.**



## HOW TO CONTACT US

<p><b>CAS Medical Systems, Inc</b> 44 East Industrial Road Branford, CT 06405 U.S.A.</p> <p>Phone: (800) 227-4414 (203) 488-6056</p> <p>Fax: (203) 488-9438</p> <p>E-Mail: custsrv@camed.com sales@camed.com techsrv@camed.com</p> <p>Web: www.camed.com</p>	<p>Representative in European Union: <b>Mossa Consulting GmbH</b> Bollbrügg 22 23570 Lübeck, Germany</p> <p>Phone: +49-4502-880-557</p> <p>Fax: +49-4502-880-559</p> <p>E-Mail: mossa.rod@t-online.de</p> <p>Should service be required, contact the dealer in the country of purchase.</p>
--	---

**Copyright 2003 CAS Medical Systems, Inc.**

**All rights reserved. No part of this manual may be reproduced without the written permission of CAS Medical Systems, Inc. CAS reserves the right to make changes to this manual and improvements to the product it describes at any time without notice or obligation.**



## TABLE OF CONTENTS

<b>1. INTRODUCTION AND INTENDED USE .....</b>	<b>3</b>
INTRODUCTION .....	3
INDICATIONS FOR USE.....	3
CONTRAINDICATIONS .....	3
BRIEF DEVICE DESCRIPTION.....	3
PATIENT ENVIRONMENT .....	5
DEFINITION OF TERMS.....	6
<b>2. UNPACKING THE MONITOR .....</b>	<b>9</b>
INITIAL INSPECTION.....	9
MONITOR CHECKLIST.....	9
OPTIONAL ACCESSORIES .....	10
<b>3. SYMBOLS .....</b>	<b>13</b>
<b>4. SAFETY MEASURES AND WARNINGS .....</b>	<b>21</b>
AUTOMATIC SAFETY FEATURES .....	25
<b>5. BLOOD PRESSURE MONITORING .....</b>	<b>29</b>
CUFF SELECTION AND APPLICATION .....	29
NIBP HOSES .....	32
<b>6. PULSE OXIMETRY MONITORING .....</b>	<b>35</b>
TAKING A SpO <sub>2</sub> MEASUREMENT .....	35
FINGER CLIP SENSORS.....	35
DISPOSABLE FLEX - TYPE SENSORS .....	36
MASIMO® OXIMETER (if so equipped).....	37
ATTACHING THE CABLES .....	37
REMOVING THE INTERFACE CABLE .....	38
MASIMO MESSAGES .....	38
NELLCOR® OXIMETER (if so equipped) .....	39
ATTACHMENT PROCEDURE.....	39
REMOVING THE INTERFACE CABLE .....	39
NONIN® OXIMETER (if so equipped).....	40
ATTACHING THE SENSOR CABLE .....	40
REMOVING THE SENSOR CABLE.....	40
<b>7. TEMPERATURE MONITORING.....</b>	<b>43</b>
WELCH ALLYN TEMPERATURE.....	43
TAKING AN ORAL TEMPERATURE .....	43
TAKING A RECTAL TEMPERATURE .....	45
TAKING AN AXILLARY TEMPERATURE.....	47

<b>8. MONITOR OPERATION</b> .....	<b>53</b>
FRONT PANEL.....	53
DIGITAL DISPLAY AND INDICATORS .....	54
FRONT PANEL CONTROLS .....	56
INFRARED (Ir) DATA PORT .....	59
REAR PANEL.....	60
AC / DC CONNECTION .....	61
FUSE COMPARTMENT.....	61
BATTERY COMPARTMENT.....	61
TEMPERATURE PROBE ELECTRICAL CONNECTION.....	61
EQUIPOTENTIALITY GROUND POST .....	61
EXTERNAL DEVICE INTERFACE.....	61
LEFT SIDE VIEW .....	62
CUFF HOSE CONNECTION.....	62
SpO <sub>2</sub> SENSOR CONNECTION.....	62
RIGHT SIDE VIEW.....	63
TEMPERATURE HOLDER .....	63
MONITOR OPERATING INSTRUCTIONS .....	64
ADULT/ NEONATE OPERATING MODE .....	64
TURNING THE CAS 740 MONITOR "ON".....	64
FRONT PANEL INTENSITY CONTROL.....	66
DISPLAYING THE TIME .....	66
MANUAL MODE FOR BLOOD PRESSURE DETERMINATION .....	66
AUTOMATIC CYCLE FOR BLOOD PRESSURE DETERMINATION.....	68
STAT MODE.....	68
HISTORY MODES.....	69
EVENT HISTORY (EVENT-Ev).....	70
TREND HISTORY (TREND-Tr).....	70
PRINT HISTORY .....	71
CLEARING HISTORY .....	72
REAL TIME CLOCK (RTC) .....	72
PATIENT ALARM MODE .....	72
CHANGING ALARM LIMITS .....	73
SAVING ALARM LIMITS .....	73
RESTORE FACTORY DEFAULTS .....	73
ALARM LIMIT VALUES .....	74
AUDIBLE AND VISUAL INDICATORS .....	75
CLEARING ALARMS.....	76
NIBP PATIENT ALARMS .....	76
HIGH / LOW %SpO <sub>2</sub> ALARMS .....	76
SpO <sub>2</sub> PULSE RATE ALARM.....	76
EQUIPMENT ALARMS.....	77
ADJUSTING THE AUDIO ALARM VOLUME.....	77
ADJUSTING THE SpO <sub>2</sub> "BEEP" VOLUME .....	77
2-MINUTE AUDIO SILENCE.....	77
PERMANENT AUDIO ALARM SILENCE.....	78
ALARM LIMITS OFF.....	79

MONITOR CONFIGURATION .....	79
ENTERING THE CONFIGURATION MENU .....	79
SAVING YOUR CHANGES .....	80
SOFTWARE REVISIONS .....	80
SELECTING THE LANGUAGE .....	81
SELECTING THE PATIENT MODE .....	81
SELECTING THE POWER-OFF DELAY TIME .....	82
SELECTING THE TEMPERATURE SCALE .....	82
AUDIO ALARM SILENCE (SILENCE/RESET Pushbutton) .....	82
2-MINUTE AUDIO ALARM SILENCE .....	83
PERMANENT AUDIO ALARM SILENCE .....	83
ALARM LIMITS OFF .....	83
MAP VALUE ENABLE / DISABLE .....	84
SET THE SpO <sub>2</sub> ALARM DELAY .....	84
SETTING THE DATE .....	85
SETTING THE TIME .....	85
DAYLIGHT SAVING TIME OPTION .....	86
BATTERY POWER .....	86
CHECKING BATTERY STATUS .....	88
AUTO OFF FEATURE .....	89
POWER FAIL .....	89
USER MESSAGES .....	89
SpO <sub>2</sub> USER MESSAGES .....	90
TEMPERATURE FUNCTION USER MESSAGES .....	91
ERROR MESSAGES ON THE MESSAGE WINDOW .....	93
<b>9. EXTERNAL PRINTER .....</b>	<b>99</b>
PRINTER OVERVIEW .....	99
PRINTER CONTROLS AND INDICATORS .....	100
PRINTER OPERATION .....	101
CHARGING THE PRINTER BATTERY .....	103
INSTALLING PAPER .....	103
REPLACING THE BATTERY PACK .....	104
INSTALLING A NEW BATTERY PACK .....	105
<b>10. CLEANING .....</b>	<b>109</b>
CLEANING OVERVIEW .....	109
THE MONITOR .....	109
THE DISPLAY .....	110
CUFFS .....	110
TUFF-CUFF® REUSABLE CUFFS .....	110
PEDISPHYG® CUFFS .....	110
PNEUMATIC TUBING .....	110
PRINTER .....	111
SpO <sub>2</sub> SENSORS .....	111
TEMPERATURE PROBES .....	111

<b>11. MAINTENANCE.....</b>	<b>115</b>
MAINTENANCE INTERVALS .....	115
TEST MODE .....	115
ENTER THE TEST MODE .....	115
EXIT THE TEST MODE.....	116
LED CHECK .....	116
CALIBRATION CHECK .....	116
SYSTEM PRESSURE .....	116
OVERPRESSURE .....	117
PNEUMATIC PRESSURE CHECKS .....	117
TEMPERATURE CALIBRATION CHECK.....	119
OXIMETRY CALIBRATION CHECK .....	120
REPLACING THE MONITOR BATTERY .....	120
REMOVING THE BATTERY .....	120
INSTALLING THE BATTERY.....	121
CHANGING THE FUSES .....	121
STORAGE .....	122
<b>12. EXTERNAL DEVICE INTERFACING .....</b>	<b>125</b>
OVERVIEW.....	125
RS232 .....	125
NURSE CALL INTERFACE.....	125
<b>13. ACCESSORIES .....</b>	<b>129</b>
BLOOD PRESSURE CUFFS .....	129
Tuff-Cuff® .....	129
Safe-Cuff™ .....	129
Pedisphyg®.....	129
INFLATION HOSES .....	129
OXIMETRY .....	130
Masimo.....	130
Nellcor.....	130
Nonin.....	130
TEMPERATURE.....	131
Welch Allyn .....	131
OTHER ACCESSORIES .....	131
MONITOR CONFIGURATIONS .....	132
<b>14. GLOSSARY .....</b>	<b>135</b>
<b>15. SPECIFICATIONS .....</b>	<b>139</b>
<b>16. PURCHASING RECORD.....</b>	<b>147</b>



## FIGURES

Figure 1: Patient Environment.....	5
Figure 2: Cuff Application .....	29
Figure 3: Cuff Positioning .....	30
Figure 4: SpO <sub>2</sub> Finger Clip Sensor Application .....	35
Figure 5: Flex – Type Adult Application .....	36
Figure 6: Flex V – Type Infant Application .....	36
Figure 7: Flex – Type Neonatal Application .....	36
Figure 8: Loading the Probe Cover.....	43
Figure 9: Location of Sublingual Pockets.....	44
Figure 10: Loading the Probe Cover.....	46
Figure 11: Loading the Probe Cover.....	48
Figure 12: Front Panel Views.....	53
Figure 13: Front Panel Controls .....	56
Figure 14: 740 Rear Panel View.....	60
Figure 15: 740M Rear Panel View.....	60
Figure 16: Left Side Panel Views .....	62
Figure 17: Right Side Panel View.....	63
Figure 18: Turning the Monitor On.....	64
Figure 19: Printer Controls and Indicators .....	100
Figure 20: Sample Printouts .....	102
Figure 21: Paper Installation.....	103
Figure 22: Opening the Battery Door.....	104
Figure 23: Installing the New Battery .....	105
Figure 24: Removing the Monitor Battery Pack.....	120
Figure 25: DB9 Male Connector Pin Layout.....	126

## TABLES

Table 1: <b>Parts of the System</b> .....	5
Table 2: <b>Neonatal Cuff Size Selection</b> .....	31
Table 3: <b>Adult - Infant Cuff Size Selection</b> .....	31
Table 4: <b>Selectable Initial Inflation Pressures</b> .....	66
Table 5: <b>Default Alarm Values</b> .....	74
Table 6: <b>Audible and Visual Indicators</b> .....	75
Table 7: <b>Software Revisions</b> .....	80
Table 8: <b>Temperature Error Codes</b> .....	92
Table 9: <b>Error Messages on the Message Window</b> .....	93
Table 10: <b>DB9 Pin Out</b> .....	126
Table 11: <b>Monitor Configurations</b> .....	132

# Section 1

## Introduction and Intended Use



# 1. INTRODUCTION AND INTENDED USE

## INTRODUCTION

The CAS 740 Monitor is a multi parameter monitor measuring blood pressure, oxygen saturation and temperature. Non-invasive blood pressure is measured using the oscillometric technique determining systolic, diastolic, mean arterial pressure and pulse rate. The pulse oximeter function continuously monitors and displays values for functional arterial hemoglobin saturation and a pulse rate. Temperature is obtained in the normal (predictive) mode in as little as four (4) seconds. A monitoring mode is available for taking axillary temperatures.

## INDICATIONS FOR USE

The CAS 740 Series Vital Signs Monitor is indicated for use for non-invasive monitoring of blood pressure, oxygen saturation, pulse and temperature of adult, pediatric and neonatal patients, in the care of health professionals.

### **WARNING:**

The CAS 740 Monitor is intended only as an adjunct in patient assessment. It must be used in conjunction with clinical signs and symptoms.

## CONTRAINDICATIONS

- Oral and Rectal Temperature measurements are not intended for neonatal use.
- Reusable SpO<sub>2</sub> sensors are contraindicated for use for prolonged periods of use. It is not intended for long term monitoring. It must be removed and repositioned every four (4) hours and if indicated by circulatory condition or skin integrity, reapplied to a different monitoring site.
- Disposable SpO<sub>2</sub> sensors are contraindicated for patients that exhibit allergic reactions to adhesive tape. The sensors must be removed and repositioned every eight (8) hours and if indicated by circulatory condition or skin integrity, reapplied to a different monitoring site.
- No other contraindications are known at this time.

## BRIEF DEVICE DESCRIPTION

The CAS 740 Monitor is compact, lightweight and portable, allowing it to be easily carried and used in a variety of clinical settings. The monitor is powered by AC Line Power, +12 VDC or by a Nickel Metal Hydride (NiMH) rechargeable battery pack. The internal battery pack charges when the monitor is plugged into a power source (AC Line Power or +12 VDC). The CAS 740 Monitor can be set to operate in one of nine (9) different languages: English, German, French, Italian, Spanish, Dutch, Swedish, Portuguese or Norwegian. The message window can display various system alarm messages. These messages direct the user to check conditions such as the battery state, air leaks and measurement problems. The message window also displays the operational mode of the monitor (automatic or manual).

## CAS 740 Monitors

The non-invasive blood pressure (NIBP) parameter automatically inflates an occluding cuff and, using the oscillometric measurement technique, determines systolic, diastolic and mean arterial pressure and pulse rate. Measurement results along with operator prompts and error messages are displayed on the front panel. The frequency of NIBP determination can be selected by the operator in varied times between one and ninety minutes. The auto and manual operating modes cover a variety of clinical uses.

The pulse oximeter parameter (%SpO<sub>2</sub>) determines arterial oxyhemoglobin saturation by measuring the absorption of red and infrared light passing through the tissue. Changes in absorption caused by pulsations of blood in the vascular bed are used to determine arterial saturation and pulse rate. The oximeter requires no routine calibration or maintenance. Oxygen saturation and heart rate are displayed on light emitting diode (LED) digital displays. On each detected pulse, the perfusion LED does indicate patient perfusion signals. This bar graph gives the user a pulse-by-pulse visual indication of waveform signal quality. An audio "beep" can be enabled that is generated each time the SpO<sub>2</sub> module detects a pulse.

**NOTE:**

The bar graph is not proportional to the pulse volume.

The temperature parameter has the capability of taking temperature in either normal (predictive) or monitor mode. In the normal mode, the thermometer's microprocessor "predicts" body temperature in about four (4) seconds for oral temperatures, about ten (10) seconds for axillary temperatures and in about fifteen (15) seconds for rectal temperatures.

Monitor mode is normally used for longer term monitoring and when difficult situations prevent accurate temperature from being taken in the predictive mode. In monitoring mode, the probe must be in contact with tissue for at least three (3) minutes for accurate oral / rectal temperature measurement and five (5) minutes for accurate axillary temperature measurement.

The default setting used by the CAS 740 Monitor for temperature determinations is the normal (predictive) mode.

**NOTE:**

Axillary temperature readings may only be taken in the Neonate monitoring mode.

**PATIENT ENVIRONMENT**

The CAS 740 Monitor has been tested with specific parts of the "system" used within the Patient Environment. Figure 1, defines the Patient Environment.

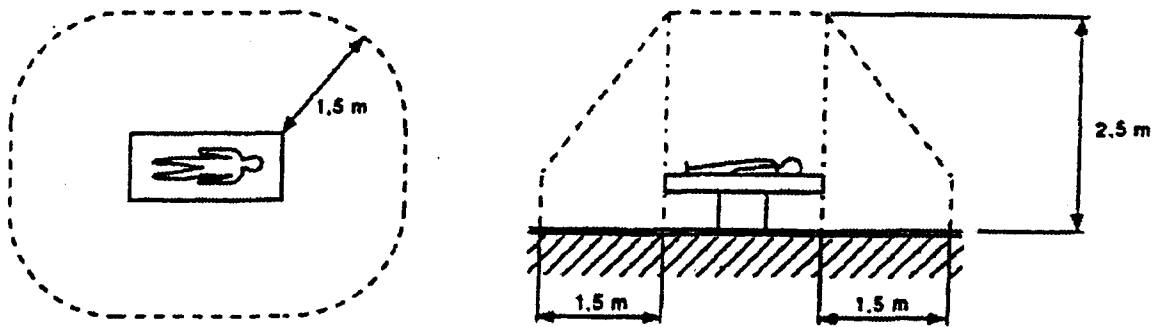


Figure 1: Patient Environment

The parts of the CAS 740 Monitor "system" that can be used in the Patient Environment are defined as;

The CAS 740 Monitor
Appropriate Accessories, listed in the ACCESSORIES section of this User's Manual
Line Cord
RS232 Interface
Citizen CMP-10 Mobile Printer
RS232 Interconnect Cable (supplied with printer)
AC Adapter / Charger, Model TRC-09-1100-M from Group West or equivalent (supplied with printer)

Table 1: Parts of the System

## DEFINITION OF TERMS

In this manual, "WARNING", "CAUTION", "IMPORTANT" and "NOTE" mean the following:

**WARNING:**

Directions that warn of conditions that put the patient or caregiver at risk.

**CAUTION:**

Directions that help you avoid damaging your monitor or losing data.

**IMPORTANT:**

Directions you should be particularly aware of; something not readily apparent.

**NOTE:**

Directions that make it easier to use your monitor.



# Section 2

## Unpacking the Monitor



## 2. UNPACKING THE MONITOR

### INITIAL INSPECTION

Before unpacking the monitor, inspect the packaging for damage. If there are any signs of damage to the package, a claim should be filed immediately with the shipping agent. It is the receiver's responsibility to notify the carrier's local office to arrange for the pickup of the damaged items. Save the damaged shipping carton as evidence.

Contact your distributor, CAS sales representative, or call CAS Medical Systems, Inc. to report external damage and to arrange for repair or replacement of damaged equipment.

The shipping carton should contain the items listed below. Unpack the monitor and account for each item. Inspect each item for signs of external damage, dents, cracks, scratches, etc. If an item is missing or damaged, contact your distributor, CAS sales representative, or CAS Medical Systems, Inc.

Record the monitor model, serial number and date of purchase at the back of this manual.

### MONITOR CHECKLIST

Qty	Description
1	CAS 740 Monitor
1	AC Power Cord or DC Power Cord – depending on model purchased
1	Ten (10) Foot Coiled Inflation Hose
1	Tuff-Cuff® Blood Pressure Cuffs, Adult
1	Tuff-Cuff® Blood Pressure Cuffs, Child
1	SpO <sub>2</sub> Interconnect Cable - For models with SpO <sub>2</sub> installed.(*)
1	SpO <sub>2</sub> Finger Sensor - For models with SpO <sub>2</sub> installed.(*)
1	SureTemp® Oral Probe and a box of Probe Covers - For models with Temperature installed.
1	P9 Calibration Kit (includes T - connector with tubing and male luer plug)
1	CAS 740 Monitor User's Manual

(\*) CAS 740 Monitors configured with Nonin SpO<sub>2</sub>, ship with a 2-meter sensor/cable assembly.

#### NOTE:

The monitor is shipped with the appropriate line cord for the country and or voltage being used.

### OPTIONAL ACCESSORIES

The CAS 740 Monitor is available with several mounting configurations to fit your needs. They consist of:

Swiveled Hard Mount  
Carrying Case

GCX Roll Stand and Basket  
Universal Mount

Refer to Section 13, ACCESSORIES for part number information. Contact CAS Medical Systems' Customer Service Department or your local distributor for more information.

# Section 3

## Symbols



### 3. SYMBOLS

Units may display the following symbols:



Alternating Current



**CAUTION:** Before using, read instructions included.



The CE Mark and Notified Body Registration Number signify the device has met all essential requirements of European Medical Device Directive 93/42/EEC.

**IPX1**

Protection against ingress of water.



Indicates this monitor is subject to the Waste Electrical and Electronic Equipment Directive in the European Union.



Medical Electrical Equipment Classification  
Class II device (if applicable)

The CAS 740 Monitor is normally a Class I device.

The CAS 740 Monitor becomes a Class II device when it is mounted and connected to a DC power source (740M).

## SYMBOLS (CONT.)



Symbol used on the rear panel of the CAS 740M, to indicate the polarity of the DC power input.



Direct Current



Indicates protection against the effects of the discharge of a cardiac defibrillator. Patient connections are Type BF and protected against defibrillation.




Equipotentiality Ground Post



NIBP Hose and Cuff Connector

SpO<sub>2</sub>

Pulse Oximeter Probe Input Connector

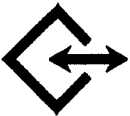

Two-way Communication Port  
RS232 Interface Connector



Temperature Probe Input Connector



## SYMBOLS (CONT.)


These symbols appear on the front panel in place of text.




ON/STANDBY – Turns “ON” the Monitor’s display.



SILENCE/RESET

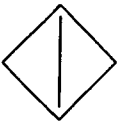

START/STAT

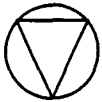


CANCEL




CYCLE TIME



HISTORY

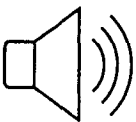

AUDIO / VOLUME

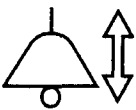

ALARM LIMITS

**SYMBOLS (CONT.)**



ARROW UP



ARROW DOWN



Bar graph display of SpO<sub>2</sub> signal strength.



Pulse Rate Display

**ADULT**


A lighted LED used to indicate NIBP operating in Adult Mode.

**NEO**

A lighted LED used to indicate NIBP operating in Neonatal Mode.

**TEMP**

A lighted LED used to indicate the Temperature Option is installed.

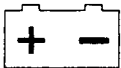


A tri-colored LED used to indicate the status of the monitors power source.

**SYMBOLS (CONT.)**

**These symbols appear on the battery pack in place of text.**



Recycling suggested (see General Notes).

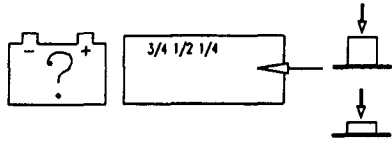

Located on the Smart Pack batteries, a set of four (4) LEDs used to indicate the approximate amount of charge remaining in the battery pack. See Page 88, CHECKING BATTERY STATUS for more information.

**These symbols appear on the packaging in place of text.**



Symbol used to indicate where Relative Humidity information concerning storage and transport can be located.



Symbol used to indicate the minimum and maximum storage and transport Temperatures.

**This symbol appears on the printer in place of text.**



**WARNING:** Before removing, read instructions in Section 9.



# Section 4

## **Safety Measures, Warnings and Precautions**



## 4. SAFETY MEASURES AND WARNINGS

### **WARNING:**

Do not use this instrument for any purpose other than specified in this manual. Doing so will invalidate the monitor's warranty.

Do not connect more than one (1) patient to the monitor.

Do not plug the monitor into an outlet controlled by a wall switch.

Before each use, verify that the alarm limits are appropriate for the patient being monitored.

The position of subject, physiological condition, and other factors affect the readings.

Blood pressure and pulse can fluctuate greatly between measurements; the monitor cannot alert the user to changes in vital signs occurring between measurement cycles.

Occasionally, electrical signals at the heart do not produce a peripheral pulse. If a patient's beat-to-beat pulse amplitude varies significantly (for example, pulsus alternans, atrial fibrillation, rapid-cycling artificial ventilator), blood pressure and pulse rate readings can be erratic and an alternate measuring method should be used for confirmation.

Where the integrity of the external protective conductor in the installation or its arrangement is in doubt, EQUIPMENT shall be operated from its INTERNAL ELECTRICAL POWER SOURCE.

Isolation of product from mains can only be achieved by removal of external power cord.

Do not, under any circumstances, perform any testing or maintenance on the monitor or power cord while the unit is being used to monitor a patient. Unplug the power cord before cleaning or servicing the monitor. The operator should not perform any servicing except as specifically stated in this manual.

Do not touch part of non-medical electrical equipment in the patient environment after removal of covers, connectors etc... without the use of a tool which operate at voltages not exceeding 25 VAC or 60 VDC and the patient at the same time.

Do not use a frayed or damaged power cord, or any accessory if you notice any sign of damage. Contact CAS Medical Systems for assistance.

Equipment not suitable for use in the presence of FLAMMABLE ANESTHETICS.

Equipment is not intended to be used in Oxygen Enriched Atmospheres.

Do not gas sterilize or autoclave the monitor.

Do not use the monitor in the presence of Magnetic Resonance Imaging (MRI) equipment.

### **WARNING:**

Do not apply the blood pressure cuff on an extremity being used for an intravenous infusion.

Do not place liquids on top of the monitor. Do not immerse the monitor or power cord in water or any liquid. If unit is accidentally wetted it should be thoroughly dried. The rear cover can be removed by a qualified service technician to verify absence of water.

During use and testing, single-use disposable temperature probe covers will limit patient cross-contamination and ensure the safety of the patient, user and device. The use of any other probe covers or failure to use a probe cover may produce temperature errors and will invalidate the monitor's warranty.

A pulse oximeter should be considered an early warning device. As a trend toward patient deoxygenation is indicated, blood samples should be analyzed by a laboratory co-oximeter to completely understand the patient's condition.

Accurate oxygen saturation measurements cannot be obtained when the oximeter is not measuring the pulse properly. If the perfusion LED is erratic or the Pulse Rate display is erratic or inaccurate, first examine the patient for any sign of distress and only then re-examine sensor placement.

**ACCURACY** – If the accuracy of any measurement does not seem reasonable, first check the patient's vital signs by alternate means and then check the CAS 740 Monitor for proper functioning.

**CABLES** – Route all cables away from patient's throat to avoid possible strangulation.

**DEFIBRILLATION** – Do not come in contact with patients during defibrillation. Serious injury or death could result.

**DISPOSAL** – Dispose of the packaging material, observing the applicable waste control regulations.

**SITE REQUIREMENTS** – For safety reasons, all connectors for patient cables and sensor leads are designed to prevent inadvertent disconnection, should someone pull on them. Do not route cables in a way that they may present a stumbling hazard. For devices installed above the patient, adequate precautions must be taken to prevent them from dropping on the patient.

### **CAUTION:**

Before each use, make sure that the monitor default alarm settings are appropriate for the specific patient being monitored.

Pressing the front panel keyswitch with a sharp or pointed instrument may permanently damage the keyswitch. Press the keyswitch using only your finger.



**CAUTION:**

Do not operate the monitor unless it has been properly calibrated. Inaccurate blood pressure readings may result. A calibration check is recommended once every year. A pneumatic check is recommended once every six (6) months.

As with any non-invasive oscillometric blood pressure monitor, the accuracy of the measurements obtained may be adversely affected by the presence of agents that alter the patient's cardiovascular system.

Do not alter the monitor's air hose. CAS Medical Systems cannot ensure proper monitor performance if the tubing is altered. Modification of the air hose will void the warranty. Avoid compression or restriction of pressure tubes.

A NIBP monitor does not operate effectively if a patient is having seizure activity, convulsions or tremors or is connected to a heart/lung machine.

In shock conditions, the low amplitude of the blood pressure waveform may make it difficult for the monitor to accurately determine the systolic and diastolic pressures.

When a patient is experiencing arrhythmias during a NIBP measurement, the accuracy of the pulse determination may be affected or the time needed to complete a measurement may be extended. The monitor will not make a determination beyond 120 seconds.

If the cuff is applied on a limb being used for oxygen saturation monitoring %SpO<sub>2</sub> results will be altered during each blood pressure measurement due to the occlusion of blood flow.

Inspect the monitor, air hose and sensors for any damage prior to operation. If any damage is noted, the monitor should not be used until it has been serviced. The monitor should be repaired only by personnel authorized to do so by CAS Medical Systems, Inc.

Use only CAS Medical Systems approved accessories and sensors to preserve the integrity, accuracy and the electromagnetic compatibility of the monitor.

Consult a physician for interpretation of blood pressure measurements.

The oximeter is factory calibrated to determine the percentage of arterial oxygen saturation of functional hemoglobin.

Significant levels of dysfunctional hemoglobins such as carboxyhemoglobin or methemoglobin may affect the accuracy of the measurement.

Cardiogreen and other intravascular dyes, depending on the concentration, may affect the accuracy of the oximeter measurement.

### **CAUTION:**

Some sensors may not be appropriate for a particular patient. If at least ten (10) seconds of one bar pulses cannot be observed for a given sensor, change sensor location or sensor type until this condition is achieved.

If the monitor fails to respond, do not use it until the situation has been corrected by qualified personnel.

**ACCIDENTAL SPILLS** – In the event that fluids are accidentally spilled on the monitor, take the monitor out of operation and inspect for damage.

**BATTERY POWER** – If the monitor will not be used or not connected to AC line power for a period over six (6) months, remove the battery.

**ELECTRICAL SHOCK** – To reduce the risk of electrical shock, do not remove the back cover. Refer all servicing to qualified personnel.

**ELECTROMAGNETIC COMPATIBILITY (EMC)** – The equipment needs special precautions regarding EMC. Be aware that strong electromagnetic fields may interfere with monitor operation. Interference prevents the clear reception of signals by the monitor. If the hospital is close to a strong transmitter such as TV, AM, or FM radio, police or fire stations, a HAM radio operator, an airport, or cellular phone, their signals could be picked up as signals by the monitor.

**ELECTROSURGERY** – Measurements may be affected in the presence of strong electromagnetic sources such as electro surgery equipment.

**GROUNDING** – Do not defeat the three-wire grounding feature of the power cord by means of adaptors, plug modifications, or other methods. Do not use extension cords of any type. Do not connect the monitor to an electrical outlet controlled by a wall switch or dimmer.

**INTERFACING OTHER EQUIPMENT** – Monitoring equipment must be interfaced with other types of medical equipment by qualified biomedical engineering personnel. Be certain to consult manufacturers' specifications to maintain safe operation.

**LEAKAGE CURRENT TEST** – The interconnection of auxiliary equipment with this device may increase the total leakage current. When interfacing with other equipment, a test for leakage current must be performed by a qualified biomedical engineering personnel before using with patients. Serious injury or death could result if the leakage current exceeds applicable standards. The use of accessory equipment not complying with the equivalent safety requirements of this equipment may lead to a reduced level of safety of the resulting system. Consideration relating to the choice shall include: use of the accessory in the patient vicinity; and evidence that the safety certification of the accessory has been performed in accordance with the appropriate IEC 601.1 and/or IEC 601.1.1 harmonized national standard.

**STACKING** – Where monitor is used adjacent to or stacked with other equipment, the monitor should be observed to verify normal operation in the configuration in which it will be used.

**GENERAL NOTES:**

There are no known risks with common disposal of equipment or accessories; however, the disposing of accessories should follow in accordance with local hospital policies. The user should ensure these policies do not conflict with any local, state or federal guidelines.

The monitor is suitable for use in the presence of electro surgery.

The monitor is suitable to be connected to public AC mains power.

The CAS 740 Monitor is not "Category AP or APG Equipment".

The CAS 740 Monitor is for "Continuous Operation".

The CAS 740 Monitor applied parts are "Type BF Defibrillation Proof".

The CAS 740 Monitor provides "DRIP-PROOF" level of protection from ingress to moisture. Do not expose the CAS 740 Monitor to extreme moisture levels such as direct exposure to rain. Exposure to extreme moisture levels may cause incorrect or inaccurate performance, or device failure during or after exposure.

**AUTOMATIC SAFETY FEATURES**

The monitor has been designed to promote patient safety. The maximum amount of time allowed to complete a blood pressure measurement is 120 seconds in adult mode and 90 seconds in neonate mode. If the measurement has not been completed within that time, the cuff is deflated automatically and a message is displayed indicating the problem.

To prevent exposure of the extremity to an inordinately high pressure, the cuff is deflated automatically when the pressure in the system is greater than 290 mmHg in the adult mode or 145 mmHg in the neonatal mode.

The cuffs used by the CAS 740 Monitor are designed without transducers for patient safety. The transducers used for NIBP measurement are located inside the monitor on the NIBP board and are isolated from the patient.

In the event of a microprocessor failure, the cuff will be deflated automatically within ten (10) seconds.

All equipment parts are protected against the effects of the discharge of a defibrillator. No separate actions are required when using this equipment with a defibrillator.

Should the AC or DC power be interrupted coming into the monitor, the monitor automatically runs off battery power. An indication of this would be a change in color of the Battery Power Visual Indicator LED from Green to either Orange or Red.

Whenever the power is disconnected from the monitor and the monitor is not allowed to shut down in an orderly fashion, the monitor, when re-powered alerts the user. Refer to Page 89, POWER FAIL for more information.

**CAUTION:**

Regardless of these safety features, always be sure to check that there are no signs of prolonged impairment of patient circulation and that the monitor is functioning properly.

# Section 5

## Blood Pressure Monitoring

## CAS 740 Monitors

Verify the hook and loop sections of the cuff are fully engaged when it is wrapped around the limb. Measurements made above the level of the heart will give reduced blood pressure readings while measurements made below the heart level will give increased readings. These errors are mainly due to the weight of the blood.

The CAS cuff is very easy to use because it does not require exact placement over the brachial artery. The fully encircling bladder ensures that the artery is properly compressed every time.

Do not compress the cuff or the cuff hose. The hose must not be kinked or pinched. It can be placed in any position.

For best results, a cuff should be wrapped for a snug fit and the limb should be positioned to be at heart level.

Do not wrap a cuff over the patient's clothing; inaccuracies could be introduced into the measurement.

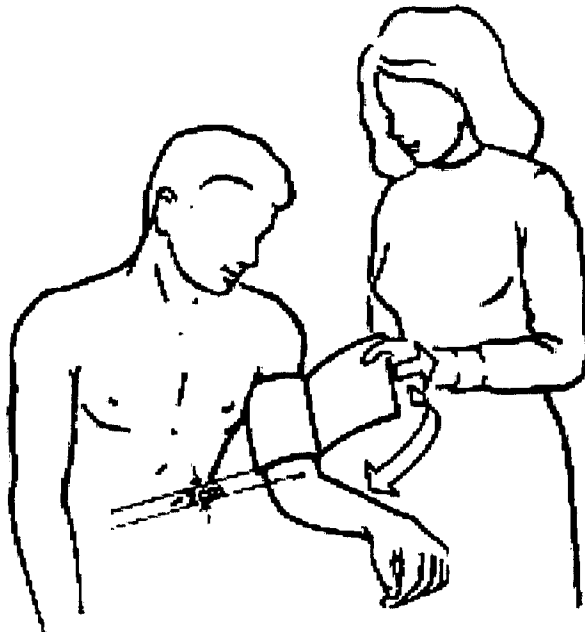


Figure 3: Cuff Positioning

**NOTE:**

Remember that there may be a marked difference between readings taken from the left and the right arms. Be consistent with each patient.

**WARNING:**

The cuff should not be applied on a limb being used for an intravenous infusion. Do not place the cuff on any extremity being used for SpO<sub>2</sub> monitoring.

**WARNING:**

When monitoring over an extended period of time, or at frequent intervals, periodically observe the patient's limb to make sure that the circulation is not impaired for a prolonged period of time.

Use the following tables as a guide to select the correct size cuff.

**Pedisphyg Neonatal Cuffs**

<b>Cuff Model Number</b>	<b>Bladder Size (Width)</b>	<b>Limb Circumference Range</b>
C26	2.5 cm (1.0 in.)	Up to 9.0 cm (3.5 in.)
C39	3.0 cm (1.2 in.)	6.0 – 11.5 cm (2.4 – 4.5 in.)
C412	4.0 cm (1.6 in.)	8.0 – 14.5 cm (3.1 – 5.7 in.)
C515	5.0 cm (2.0 in.)	10.0 – 17.75 cm (3.9 – 7.0 in.)

**Table 2: Neonatal Cuff Size Selection**

**Adult and Infant Cuffs**

<b>Cuff Model Number</b>	<b>Bladder Size (Width)</b>	<b>Limb Circumference Range</b>
CR5206	6 cm (2.4 in.)	11 – 18 cm (4.3 – 7.1 in.)
CR5207	7 cm (2.8 in.)	14 – 21 cm (5.5 – 8.3 in.)
CR5209	9 cm (3.5 in.)	18 – 27 cm (7.1 – 10.6 in.)
CR5212	12 cm (4.7 in.)	21 – 30 cm (8.3 – 11.8 in.)
CR5214	14 cm (5.5 in.)	24 – 37 cm (9.4 – 14.6 in.)
CR5216	16 cm (6.3 in.)	28 – 42 cm (11.0 – 16.5 in.)

**Table 3: Adult - Infant Cuff Size Selection**

### NIBP HOSES

CAS Medical offers two (2) NIBP Inflation Hoses.

For Adult and Pediatric patient monitoring using the ADULT NIBP mode and Cuffs, the coiled ten (10) foot NIBP hose (CAS p/n 01-02-0131) is recommended.

For Neonatal and Infant patient monitoring using the NEO NIBP mode and Neonatal cuffs, the Neonatal six (6) foot NIBP hose (CAS p/n 01-02-0185) is recommended.

See Section 13, ACCESSORIES for CAS Medical Systems cuff size and part number information.

**NOTE:**

Avoid compression or restriction of NIBP hoses.



# Section 6

## Pulse Oximetry Monitoring



## 6. PULSE OXIMETRY MONITORING

The CAS 740 Monitor features a wide variety of SpO<sub>2</sub> technology that is ideal for every application.

### TAKING A SpO<sub>2</sub> MEASUREMENT

**IMPORTANT:**

Prior to patient monitoring, ensure the monitor is configured to the appropriate patient mode – Neonate or Adult.

The following is a general procedure for taking a SpO<sub>2</sub> measurement:

- 1) Select a sensor based on the patient size and monitoring conditions and properly attach the sensor to the patient.

### FINGER CLIP SENSORS

The finger clip sensor is designed for spot check monitoring of pediatric and adult patients or continuous monitoring where patient movement is not expected, and the patient's finger is large enough for the sensor to be attached securely.

**NOTE:**

If patient movement is occurring or the finger size is inappropriate, select a different sensor that is appropriate for the patient and the monitoring environment.

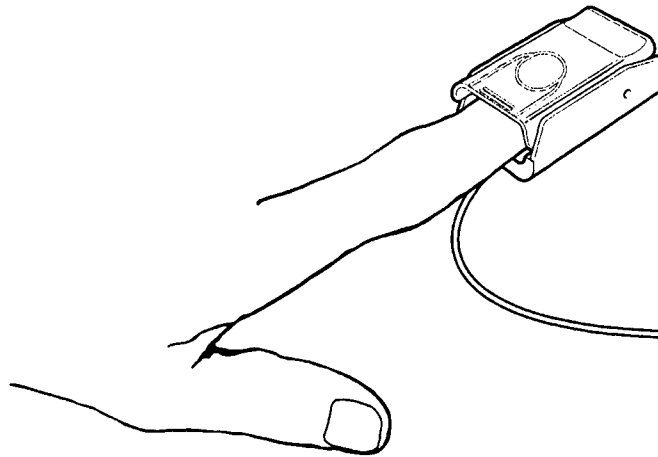


Figure 4: SpO<sub>2</sub> Finger Clip Sensor Application

Insert the finger (preferably left or right index finger) completely into the sensor. Place the sensor with the LEDs positioned on the nail side. The thumb is specifically not recommended for use with the finger clip sensor.


**DISPOSABLE FLEX - TYPE SENSORS**

These types of sensors are designed for patients as a single patient use sensor and is intended for use where moderate patient movement is expected or cross contamination is possible.




**Figure 5: Flex – Type Adult Application**

Adult and Pediatric: The preferred application site is the index finger, however, other fingers or toes may be used where the tissue thickness is between 5 and 17 millimeters.

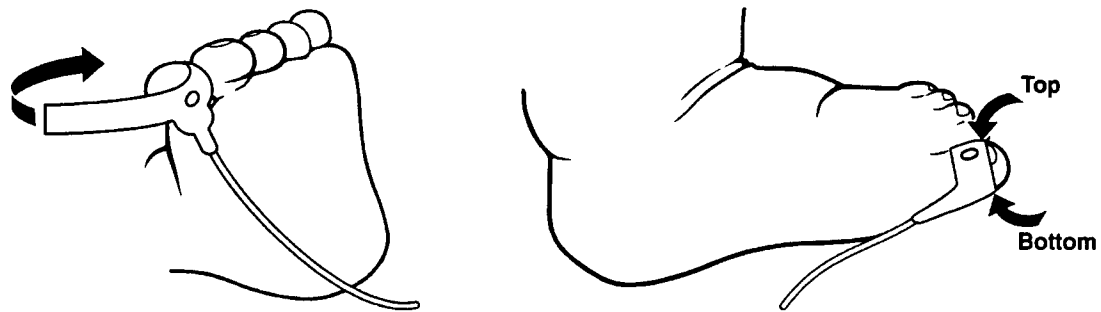

**Figure 6: Flex V – Type Infant Application**

Infant: The preferred application site is the large toe of infants greater than 2 kilograms in weight.

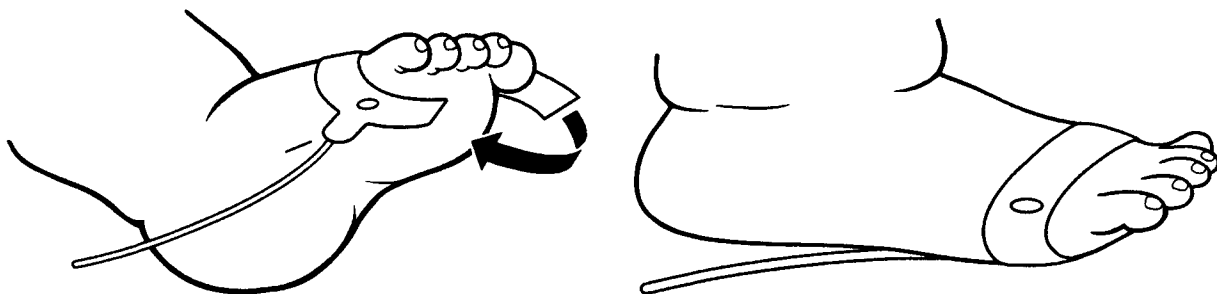

**Figure 7: Flex – Type Neonatal Application**

Neonatal: The preferred application site is on the foot close to the toes for infants less than 2 kilograms in weight.

**NOTE:**

For best results, secure the sensor cable independently from the sensor, preferably around the base of the finger. Tape may be used to secure the cable to the patient. Make sure that the tape being used does not restrict the blood flow.

- 2) Once the sensor has been attached to the patient and to the monitor, the SpO<sub>2</sub> signal bar graph will illuminate and indicate the relative signal strength and signal quality at the sensor site. Wait for the monitor to determine the initial %SpO<sub>2</sub> and Pulse Rate values. When the values have been determined, they will be displayed in their respective %SpO<sub>2</sub> and ♥BPM front panel LED display windows.

**NOTE:**

Inspect the pulse oximeter site every 2 to 4 hours or per hospital protocol. If there is any skin irritation caused by the sensor, remove the sensor and apply it to a different location.

- 3) When SpO<sub>2</sub> monitoring is completed, remove the sensor from the patient.

When the probe is removed from the patient, the message "Prb OFF" appears in the Message Window and an audible alarm sounds, indicating a connection has been lost.

Press the SILENCE/RESET pushbutton. The monitor silences the audible alarm tone, but the message remains.

**NOTE:**

If either the 2-Minute Audio Off or Permanent Audio Off is enabled, no audio will be heard but a visual message will appear in the Message Window.

**MASIMO® OXIMETER (if so equipped)**

The CAS 740 Monitor can be equipped to use SpO<sub>2</sub> sensors manufactured by Masimo. No other manufacturer's sensors should be used.

**NOTE:**

CAS 740 Monitors equipped with Masimo oximetry will have the Masimo SET® logo next to the SpO<sub>2</sub> connector.

**ATTACHING THE CABLES**

- 1) Select and apply a sensor that is the appropriate size for the patient's digit or extremity, according to the instructions provided by Masimo.
- 2) Orient the connecting tab of the sensor so that the shiny contacts are pointed up. Mate the logo on the sensor to the logo on the patient cable.
- 3) Insert the tab of the sensor into the patient cable connector until there is a tactile or audible "click" connection. Verify a secure connection and gently tug on the patient cable connector.
- 4) Plug the Interface Cable into the SpO<sub>2</sub> connector on the side panel of the monitor. The connector is shaped like a "D". Line up the "D" on the Interface Cable with the "D" on the receptacle. Push the cable in until you hear an audible "click".

- 5) Press the ON/STANDBY pushbutton to turn "ON" the monitor.
- 6) Check the Alarm Limits and configure them appropriately for the patient. Refer to Page 72, PATIENT ALARM MODE for more information.

### REMOVING THE INTERFACE CABLE

When SpO<sub>2</sub> monitoring is not required, disconnect the Interface Cable by squeezing the grey tabs with your thumb and index finger while pulling the connector away from the monitor.

**NOTE:**

To avoid damage to the Interface Cable, always hold it by the connector rather than the cable when connecting or disconnecting either end.

When the probe is disconnected from the monitor, the message "No Probe" appears in the Message Window and an audible alarm sounds indicating a connection has been broken. Press the SILENCE/RESET pushbutton to silence the visual and audible alarm.

**NOTE:**

If either the 2-Minute Audio Off or Permanent Audio Off is enabled, no audio will be heard but a visual message will appear in the Message Window.

Refer to Section 13, ACCESSORIES for Masimo oximeter probe types and part number information. Consult instructions enclosed with each sensor for proper application.

### MASIMO MESSAGES

When the message "Low Perf" or "Low Sig" appear in the Message Window, and the monitor is displaying valid %SpO<sub>2</sub> numerics, no audible alarm will be heard.

When the message "Low Perf" or "Low Sig" appear in the Message Window, and the monitor is not displaying valid %SpO<sub>2</sub> numerics, an audible alarm will be heard. Press the SILENCE/RESET pushbutton. The monitor silences the audio alarm tone, but the message remains.

**NELLCOR® OXIMETER (if so equipped)**

The CAS 740 Monitor can be equipped to use Nellcor OxiMax® sensors. No other manufacturer's sensors should be used.

**NOTE:**

CAS 740 Monitors equipped with Nellcor oximetry will have the Nellcor OxiMax logo next to the SpO<sub>2</sub> connector.

**ATTACHMENT PROCEDURE**

- 1) Select and apply a sensor that is the appropriate size for the patient's digit or extremity, according to the instructions provided by Nellcor.
- 2) Connect the sensor assembly to the Interface Cable:
  - a) Place the plastic hinged cover in the unlocked position (perpendicular to the connector).
  - b) Connect the sensor assembly to the Interface Cable.
  - c) Lock the plastic hinged cover to prevent accidental cable disconnection.
- 3) Plug the Interface Cable into the SpO<sub>2</sub> connector on the side panel of the monitor. The connector is shaped like a "D". Line up the "D" on the Interface Cable with the "D" on the receptacle. Push the cable in until you hear an audible "click".
- 4) Press the ON/STANDBY pushbutton to turn "ON" the monitor.
- 5) Check the Alarm Limits and configure them appropriately for the patient. Refer to Page 72, PATIENT ALARM MODE for more information.

**REMOVING THE INTERFACE CABLE**

When SpO<sub>2</sub> monitoring is not required, disconnect the Interface Cable by squeezing the grey tabs with your thumb and index finger while pulling the connector away from the monitor.

**NOTE:**

To avoid damage to the Interface Cable, always hold it by the connector rather than the cable when connecting or disconnecting either end.

When the probe is disconnected from the monitor, the message "No Probe" appears in the Message Window and an audible alarm sounds indicating a connection has been broken. Press the SILENCE/RESET pushbutton to silence the visual and audible alarm.

**NOTE:**

If either the 2-Minute Audio Off or Permanent Audio Off is enabled, no audio will be heard but a visual message will appear in the Message Window.

Refer to Section 13, ACCESSORIES for Nellcor oximeter probe types and part number information. Consult instructions enclosed with each sensor for proper application.

### NONIN® OXIMETER (if so equipped)

The CAS 740 Monitor can be equipped to use SpO<sub>2</sub> sensors manufactured by Nonin. No other manufacturer's sensors should be used.

**NOTE:**

CAS 740 Monitors equipped with Nonin oximetry will have the Nonin logo next to the SpO<sub>2</sub> connector.

### ATTACHING THE SENSOR CABLE

- 1) Select and apply a sensor that is the appropriate size for the patient's digit or extremity, according to the instructions provided by Nonin.
- 2) Plug the sensor assembly into the SpO<sub>2</sub> connector on the side panel of the monitor. The connector is shaped like a "D". Line up the "D" on the sensor cable with the "D" on the receptacle. Push the connector in completely.
- 3) Press the ON/STANDBY pushbutton to turn "ON" the monitor.
- 4) Check the Alarm Limits and configure them appropriately for the patient. Refer to Page 72, PATIENT ALARM MODE for more information.

### REMOVING THE SENSOR CABLE

When SpO<sub>2</sub> monitoring is not required, disconnect the Sensor Cable by carefully removing the connector from the SpO<sub>2</sub> input receptacle.

**NOTE:**

To avoid damage to the Sensor Cable, always hold it by the connector rather than the cable when connecting or disconnecting either end.

When the probe is disconnected from the monitor, the message "No Probe" appears in the Message Window and an audible alarm sounds indicating a connection has been broken. Press the SILENCE/RESET pushbutton to silence the visual and audible alarm.

**NOTE:**

If either the 2-Minute Audio Off or Permanent Audio Off is enabled, no audio will be heard but a visual message will appear in the Message Window.

Refer to Section 13, ACCESSORIES for Nonin oximeter probe types and part number information. Consult instructions enclosed with each sensor for proper application.



# Section 7

## Temperature Monitoring



## 7. TEMPERATURE MONITORING

### WELCH ALLYN TEMPERATURE

**WARNING:**

During use, single-use disposable temperature probe covers supplied by CAS Medical Systems or Welch Allyn will limit patient cross-contamination. The use of any other probe cover or the failure to use a probe cover may produce temperature errors and will invalidate the monitor's warranty. Temperature probe covers are required to ensure the safety of the patient and user.

### TAKING AN ORAL TEMPERATURE

**CAUTION:**

Never use the Temperature Probe without a probe cover. Accurate oral temperatures can only be obtained by using the *blue* temperature probe. The use of the wrong probe will produce temperature errors.

To take an oral temperature (in either Predictive or Monitor mode) follow this procedure:

- 1) Insure that the oral probe is connected to the monitor and that the probe is secured into its holder. An audible "click" should be heard when the probe is completely placed into its holder. The oral probe has a BLUE probe cover ejection button.
- 2) Remove the probe from the probe holder. A short self-test mode will be initiated where every LED segment in the TEMP display window are illuminated briefly. Following this self-test the display will show "OrL". The "OrL" in the TEMP display window indicates that the oral algorithm will be used to take a predictive temperature measurement.
- 3) Load a probe cover onto the probe by holding the probe collar with the thumb and forefinger, being careful not to hold or press the ejection button. Refer to Figure 8.

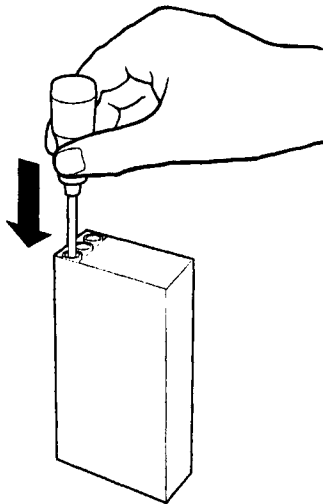


Figure 8: Loading the Probe Cover


- 4) When a double dash “--” appears in the TEMP display window, the temperature sensor is ready to take a measurement. To take a normal predictive temperature, simply insert the probe tip gently into the patient’s slightly opened mouth. Carefully slide the probe under the tongue on either side of the mouth to reach the sublingual pocket. Refer to Figure 9.

To take a direct read temperature (monitor mode) wait 1 minute until a temperature value appears in the TEMP display. The decimal point associated with the temperature value will flash at a rate of one (1) second “ON” and one (1) second “OFF” as a reminder that you are in monitor mode. Insert the probe tip into the patient’s mouth as described above.

**NOTE:**

Accurate temperatures can only be obtained in this location. Temperatures in other mouth locations can vary by as much as 2°F or 1°C.

**Sublingual Pockets**

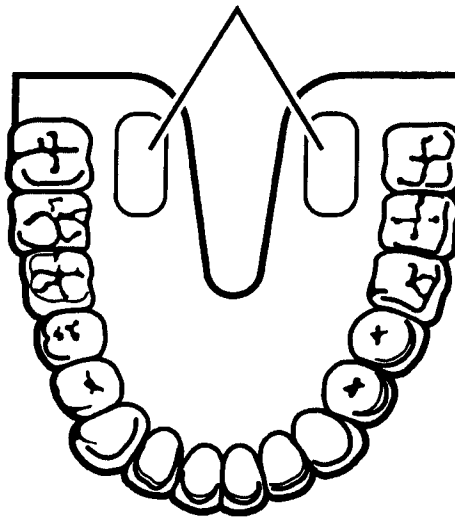

**Figure 9: Location of Sublingual Pockets**

- 5) The probe should be held by the clinician during the entire temperature measurement process to insure the probe tip maintains tissue contact.
- 6) During the predictive temperature measurement cycle, a walking  $\cong$  segments will appear in the TEMP display window indicating that a predictive measurement is in process.
- 7) When the final temperature has been reached, the temperature value will be displayed in the TEMP display window and an audible tone will be generated. The time of the measurement and the measured temperature will be stored in the monitor’s history memory.

**NOTE:**

To change the Temperature scale, refer to Page 82, SELECTING THE TEMPERATURE SCALE.

**NOTE:**

If the probe tip did not maintain tissue contact during the entire predictive measurement, the final temperature displayed in the TEMP display window will flash. If this occurs it is recommended that a new temperature be taken. Depress the SILENCE/RESET pushbutton or after waiting for two (2) minutes the monitor will automatically blank the display.

- 8) After the temperature measurement is complete, remove the probe from the patient's mouth and eject the probe cover by firmly pressing the ejection button on the probe. Properly dispose of the used probe cover per protocol.
- 9) Insert the probe into the probe holder before attempting to take another temperature measurement. An audible "click" should be heard when the probe is completely placed into its holder.

**NOTE:**

Failure to correctly place the probe back into its holder may result in failure of the next predictive temperature measurement.

- 10) Following a completed temperature measurement, the current temperature measurement is displayed for five (5) minutes after which time the TEMP display will go blank.

**TAKING A RECTAL TEMPERATURE****CAUTION:**

Never use the Temperature Probe without a probe cover. Accurate rectal temperatures can only be obtained by using the *red* temperature probe. The use of the wrong probe will produce temperature errors.

To take a rectal temperature, follow this procedure:

- 1) Insure that the rectal probe is connected to the monitor and that the probe is secured into its holder. An audible "click" should be heard when the probe is completely placed into its holder. The rectal probe has a RED probe cover ejection button.
- 2) Remove the probe from the probe holder. A short self-test mode will be initiated where every LED segment in the TEMP display window are illuminated briefly. Following this self-test the display will show "rEC". The "rEC" in the TEMP display window indicates that the rectal algorithm will be used to take a predictive temperature measurement.
- 3) Load a probe cover onto the probe by holding the probe collar with the thumb and forefinger, being careful not to hold or press the ejection button. Refer to Figure 10.

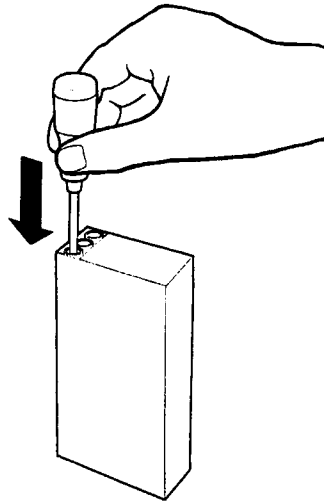


Figure 10: Loading the Probe Cover

- 4) When a double dash “--” appears in the TEMP display window, the temperature sensor is ready to take a measurement. To take a normal predictive temperature, separate the buttocks with one hand. Apply a thin coat of water-based lubricant when necessary. Using the other hand, gently insert the probe ONLY 1cm (3/8 inch ONLY) inside the rectal sphincter.

To take a direct read temperature (monitor mode) wait 1 minute until a temperature value appears in the TEMP display. The decimal point associated with the temperature value will flash at a rate of one (1) second “ON” and one (1) second “OFF” as a reminder that you are in monitor mode. Insert the probe tip into the patient’s rectal sphincter as described above.

**WARNING:**

Extreme caution should be used to avoid risk of bowel perforation in children.

- 5) Tilt the probe to insure good tissue contact and continue to keep the buttocks separated while the measurement is in process.
- 6) During the predictive temperature measurement cycle, “walking” segments will appear in the TEMP display window indicating that a predictive measurement is in process.
- 7) When the final temperature has been reached, the temperature value will be displayed in the TEMP display window and an audible tone will be generated. The time of the measurement and the measured temperature will be stored in the monitor’s history memory.

**NOTE:**

To change the Temperature scale, refer to Page 82, SELECTING THE TEMPERATURE SCALE.

**NOTE:**

If the probe tip did not maintain tissue contact during the entire predictive measurement, the final temperature displayed in the TEMP display window will flash. If this occurs it is recommended that a new temperature be taken. Depress the SILENCE/RESET pushbutton or after waiting for two (2) minutes the monitor will automatically blank the display.

- 8) After the temperature measurement is complete, remove the probe from the patient's rectum and eject the probe cover by firmly pressing the ejection button on the probe. Properly dispose of the used probe cover per protocol.
- 9) Insert the probe into the probe holder before attempting to take another temperature measurement. An audible "click" should be heard when the probe is completely placed into its holder.

**NOTE:**

Failure to correctly place the probe back into its holder may result in failure of the next predictive temperature measurement.

- 10) Following a completed temperature measurement, the current temperature measurement is displayed for five (5) minutes after which time the TEMP display will go blank.

**TAKING AN AXILLARY TEMPERATURE****CAUTION:**

Never use the Temperature Probe without a probe cover. Accurate axillary temperatures can only be obtained by using the *blue* temperature probe. The use of the wrong probe will produce temperature errors.

To take an axillary temperature in Normal (Predictive) mode follow this procedure:

- 1) Insure that the oral probe is connected to the monitor, secured into its holder and that the monitor is in NEO mode (Axillary temperatures can only be taken in NEO mode). An audible "click" should be heard when the probe is completely placed into its holder. The oral probe has a BLUE probe cover ejection button.
- 2) Remove the probe from the probe holder. A short self-test mode will be initiated where every LED segment in the TEMP display window are illuminated briefly. Following this self-test the display will show "ALy". The "ALy" in the TEMP display window indicates that the axillary algorithm will be used to take a predictive temperature measurement.
- 3) Load a probe cover onto the probe by holding the probe collar with the thumb and forefinger, being careful not to hold or press the ejection button. Refer to Figure 11.

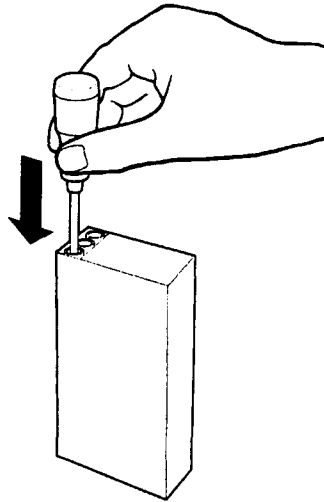


Figure 11: Loading the Probe Cover

When a double dash “--” appears in the TEMP display window the temperature sensor is ready to take a measurement. To take a normal predictive temperature, simply place the probe as high as possible in the axilla. Do not allow the probe tip to come into contact with the patient until it is deliberately placed in the measurement site.

- 4) Be sure that the probe tip is completely surrounded by axillary tissue. Clothing or any other material touching the probe tip may cause inaccurate readings.
- 5) Place the arm close to the patient’s side. Hold the arm in this position without movement of the arm or probe during the measurement cycle.
- 6) During the predictive temperature measurement cycle, “walking” segments will appear in the TEMP display window, indicating that a predictive measurement is in process.
- 7) When the final temperature has been reached, the temperature value will be displayed in TEMP display window and an audible tone will be generated. The time of the measurement and the measured temperature will be stored in the monitor’s history memory.

**NOTE:**

To change the Temperature scale, refer to Page 82, SELECTING THE TEMPERATURE SCALE.

**NOTE:**

If the probe tip did not maintain tissue contact during the entire predictive measurement, the final temperature displayed in the TEMP display window will flash. If this occurs it is recommended that a new temperature be taken. Depress the SILENCE/RESET pushbutton or after waiting for two (2) minutes the monitor will automatically blank the display.



- 8) After the temperature measurement is complete, remove the probe from the patient's axilla and eject the probe cover by firmly pressing the ejection button on the probe. Properly dispose of the used probe cover per protocol.
- 9) Insert the probe into the probe holder before attempting to take another temperature measurement. An audible "click" should be heard when the probe is completely placed into its holder.

**NOTE:**

Failure to correctly place the probe back into its holder may result in failure of the next predictive temperature measurement.

- 10) Following a completed temperature measurement, the current temperature measurement is displayed for five (5) minutes after which time the TEMP display will go blank.

Refer to Section 13, ACCESSORIES for Welch Allyn temperature probe types and part number information. Consult instructions enclosed with each probe for proper application.



# Section 8

## Monitor Operation




## 8. MONITOR OPERATION

### FRONT PANEL



Front Panel Keyswitch panel with English text

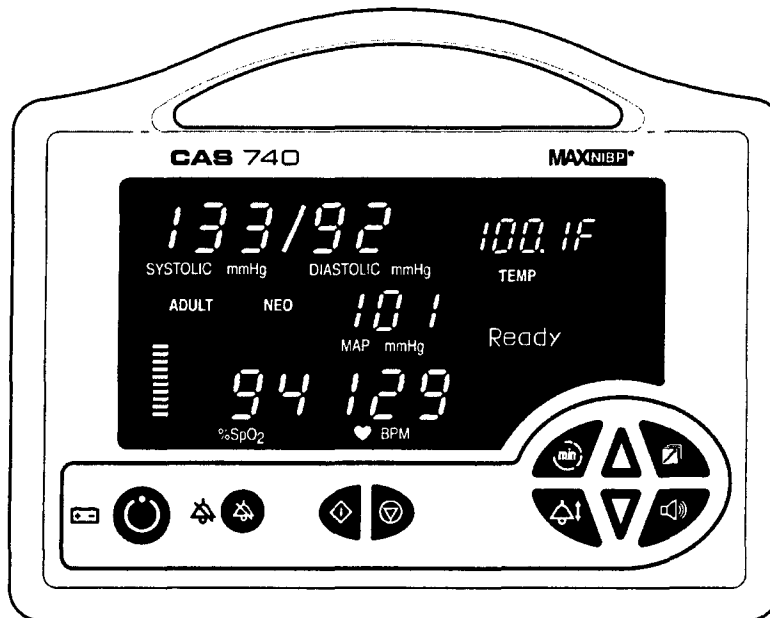

Front Panel Keyswitch panel, Symbols only

Figure 12: Front Panel Views

**DIGITAL DISPLAY AND INDICATORS**


**SYSTOLIC mmHg** Red colored LEDs indicate the Systolic pressure measurement in mmHg.


**DIASTOLIC mmHg** Red colored LEDs indicate the Diastolic pressure measurement in mmHg.

**MAP mmHg** Red colored LEDs indicate the Mean Arterial Pressure in mmHg (if enabled).

**TEMP** A Yellow LED indicator with Red colored LEDs indicates the temperature value (if installed).

**% SpO<sub>2</sub>** Green colored LEDs indicates the %SpO<sub>2</sub> value (if installed).

 **BPM** Red colored LEDs indicate the Pulse Rate in BPM (beats per minute).

 Green colored LEDs provide a visual indication of the SpO<sub>2</sub> signal strength in a bar graph form (if installed).


**ADULT** A yellow LED indicator used to inform the user that the NIBP is operating in the Adult Mode.

**NEO** A yellow LED indicator used to inform the user that the NIBP is operating in the Neonatal Mode.

**TEMP** A yellow LED indicator used to inform the user that the Temperature Option is installed.

## Ready

Message Window area used to display various messages that aid the user in monitor operation.

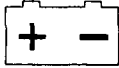

A tri-colored visual indicator used to display the status of the power source and battery condition.

The status of the LED is:

- GREEN = Monitor is connected to a main power source
- ORANGE = In Use on Battery
- RED = Battery Low or Dead Battery



A Yellow LED visual indicator used along with the SILENCE/RESET pushbutton to display the status of the Audio Alarm Silence feature. Refer to Page 77 for more information.

The status of the LED is:

- "ON" continuously = 2 Minute Audio Silence
- Flash one second "ON"/one second "OFF" = Permanent Audio Alarm Silence

## FRONT PANEL CONTROLS

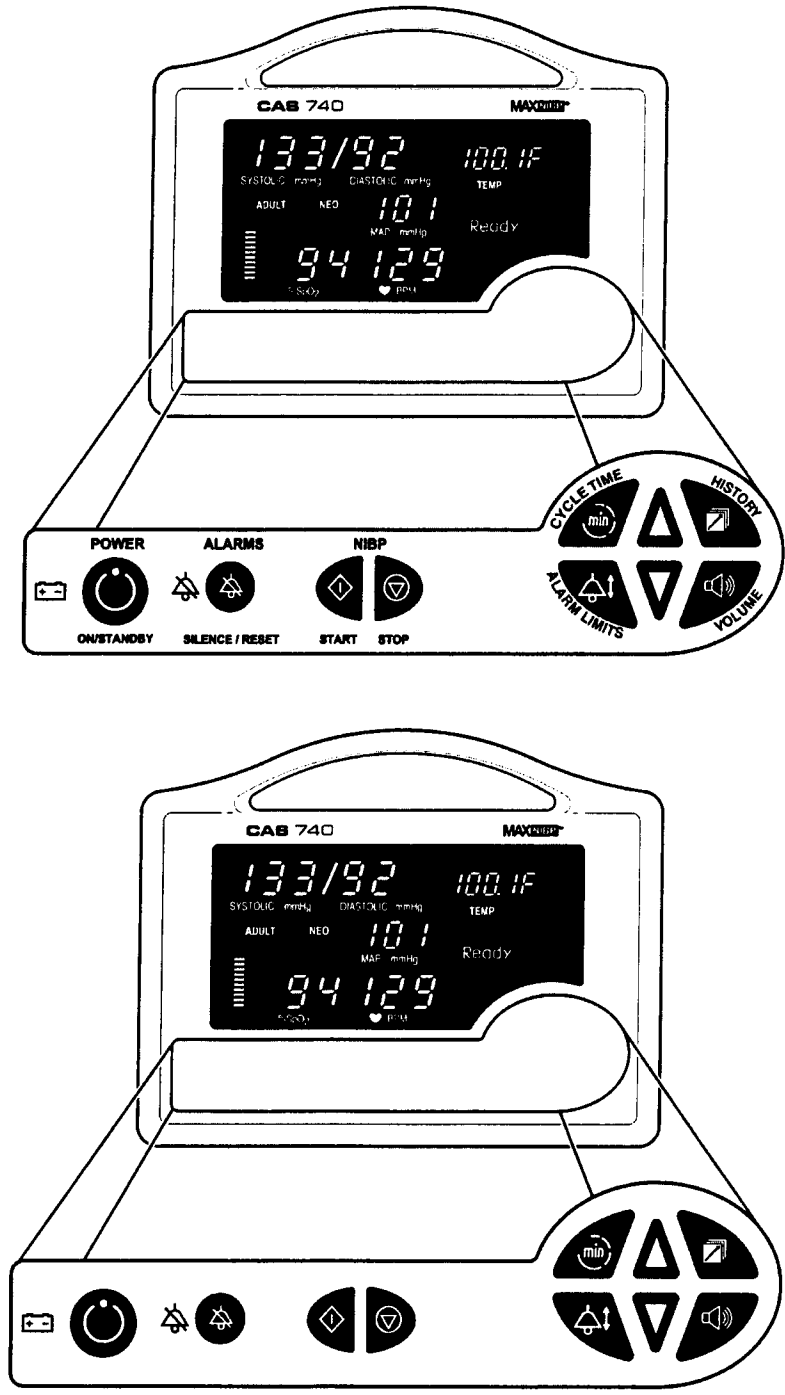


Figure 13: Front Panel Controls



**ON/STANDBY:**

Depress once, turns "ON" the CAS 740 Monitor's display (if it was OFF).

To turn the monitor's display "OFF", either depress once or depress and hold for two (2) seconds. Based on the selection for Power-Off Delay Time in the Configuration menu. Refer to Page 82 for more information.


**SILENCE/RESET:**

When depressed during an active patient alarm, silences the audio portion of that alarm for fifteen (15) seconds.

When depressed during an active equipment alarm, the alarm condition shall be acknowledged along with the audio and visual shall be removed.

Used to enable and disable the two (2) Minute Audio Alarm Silence or Permanent Audio Alarm Silence feature. Refer to Page 82, AUDIO ALARM SILENCE (SILENCE/RESET Pushbutton) for more information.


Allows the user to clear NIBP, SpO<sub>2</sub> and Temperature messages from the front panel display.

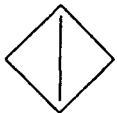**START/STAT:****START:**

Initiates a blood pressure measurement in the Manual Mode or begins the selected Automatic Cycle.

**STAT:**

Starts a series of NIBP measurements (depress and hold for two (2) seconds). Continues for five minutes.

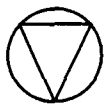**CANCEL:**

Cancels any active blood pressure function and immediately deflates the cuff.

Also used to cancel out of a menu and return to "Ready".



**CYCLE TIME:**

Allows the user to select a time interval for Automatic blood pressure measurement.

Automatic measurement cycles of 1, 2, 3, 4, 5, 10, 15, 30, 60, or 90 minutes may be chosen.

Also used to read the monitor's current time (depress and hold for two (2) seconds).



**HISTORY:**


Allows the user to review stored patient readings. Refer to Page 69, HISTORY MODES for more information.



**AUDIO / VOLUME:**

Used to set the volume level of the Alarms and the SpO<sub>2</sub> Beep (if SpO<sub>2</sub> installed). Use the Up and Down Arrows to change.

Also used to adjust the brightness of the front panel displays (depress and hold for two (2) seconds).

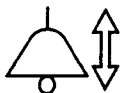

**ALARM LIMITS:**

Allows the user to enter and set the monitor's Alarm Limits.



**ARROW UP:**

Allows forward Adjustment (Auto Cycle, History, Inflation Pressure, Limits and Monitor Configuration).

Depress to cycle through menu selections or depress and hold for quicker advance.



**ARROW DOWN:**

Allows backwards Adjustment (Auto Cycle, History, Inflation Pressure, Limits and Monitor Configuration).

Depress to cycle through menu selections or depress and hold for quicker advance.

**NEXT**

The HISTORY and AUDIO pushbutton keys have been programmed to allow the user to advance forward to the next selection in the Monitor Configuration menu.

**PREVIOUS**

The CYCLE TIME and ALARM LIMITS pushbutton keys have been programmed to allow the user to advance backwards to the previous selection in the Monitor Configuration menu.

**INFRARED (Ir) DATA PORT**

An Infrared (Ir) output port, located on the bottom panel of the monitor's front cover, is available to print the NIBP, %SpO<sub>2</sub> and Temperature History data to the optional external printer or other data collection device(s).

REAR PANEL

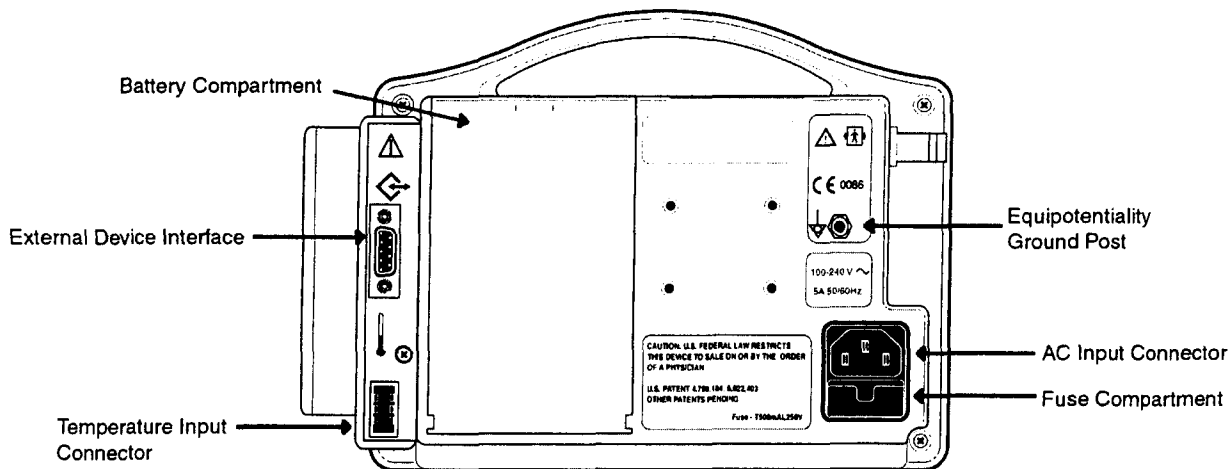


Figure 14: 740 Rear Panel View

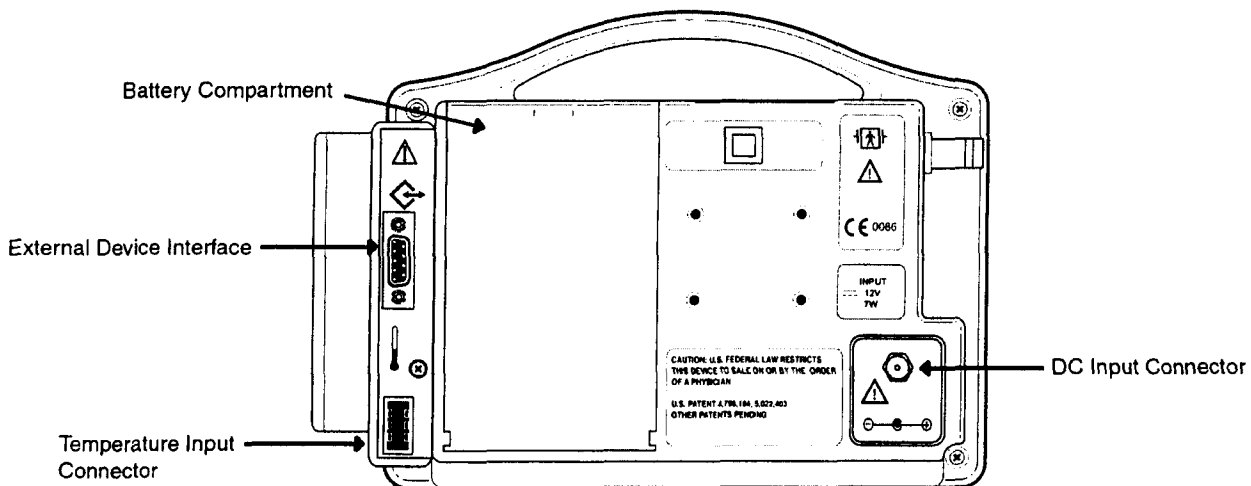


Figure 15: 740M Rear Panel View

## AC / DC CONNECTION

Receptacle for the AC power cord when the CAS 740 Monitor is purchased with an AC Line Power option and built-in battery.

When the CAS 740 Monitor is purchased as an EMS monitor equipped with a Swivel Mount, the ambulance DC power cord is attached here.

## FUSE COMPARTMENT

When the CAS 740 Monitor is purchased for AC Line power, the power input receptacle incorporates dual fuses located in the hot and neutral lines.

## BATTERY COMPARTMENT

The CAS 740 Monitor is equipped with a 7.2 Volt, 3700 mAh battery pack that, when fully charged, is capable of taking 100 NIBP readings when the monitor is set in the 5-minute Automatic Mode.

### NOTE:

The serial number label is located on the bottom of the monitor.

## TEMPERATURE PROBE ELECTRICAL CONNECTION

(if equipped)

Connect the Temperature cable in this receptacle for Temperature monitoring. Leave the probe connection in place. The Temperature function is OFF until the probe is removed from the holder. Refer to Section 7, TEMPERATURE MONITORING for more information.

## EQUIPOTENTIALITY GROUND POST

This terminal can be used to provide an auxiliary ground for the monitor.

## EXTERNAL DEVICE INTERFACE

The CAS 740 Monitor is available with a combined optional DB9 RS232 output and Nurse Call interface. The RS232 output may be used to interface to the Citizen CMP-10 Mobile Printer or another serial printing device.

Refer to Section 12, EXTERNAL DEVICE INTERFACING for more information.

LEFT SIDE VIEW

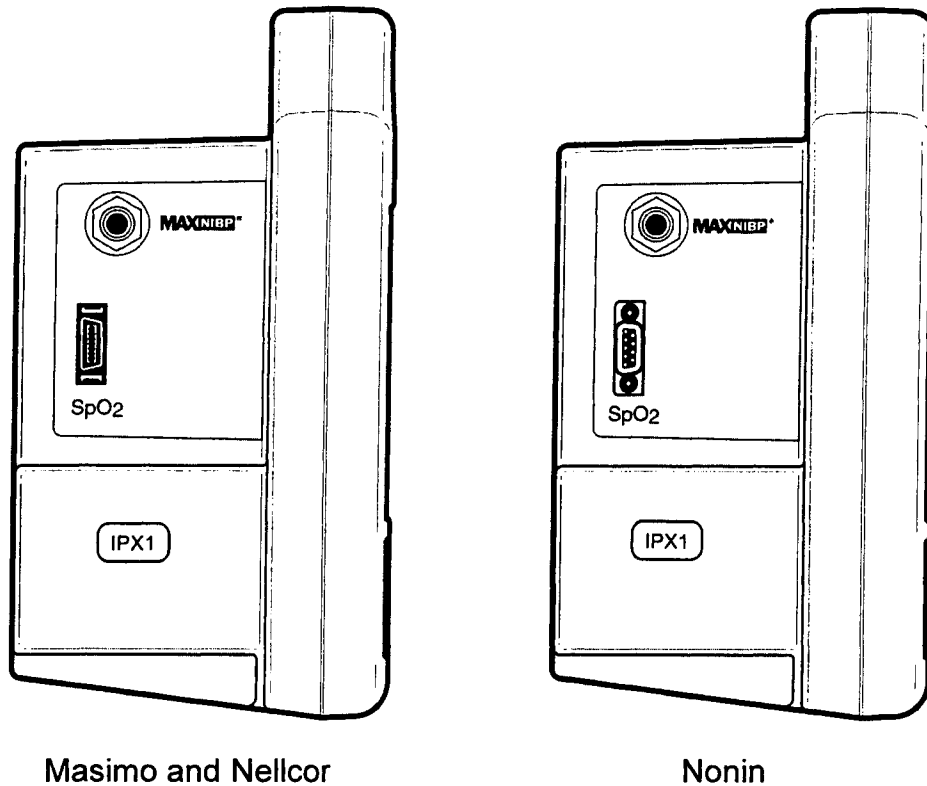


Figure 16: Left Side Panel Views

**MAXNIBP®**

**CUFF HOSE CONNECTION**

The inflation hose is connected to the monitor where the MAXNIBP logo is located as shown in Figure 16. The hose must be connected to the cuff prior to use.

**NOTE:**

An optional six (6) foot inflation hose is available when monitoring in the NEO mode. Refer to Section 13, ACCESSORIES for part number information.

**SpO<sub>2</sub> SENSOR CONNECTION**  
(if equipped)

Connect the sensor cable in this receptacle for SpO<sub>2</sub> monitoring.

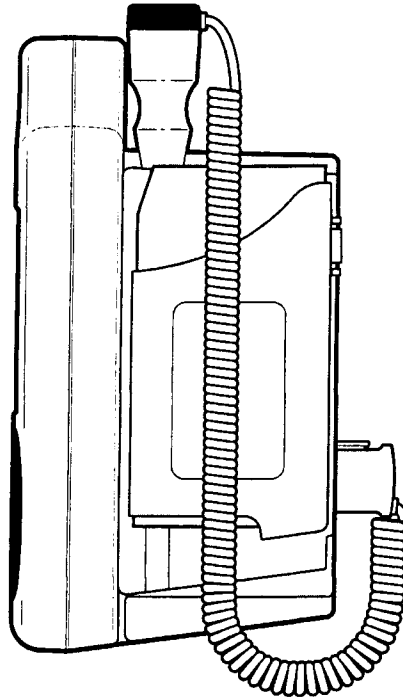
**RIGHT SIDE VIEW**

Figure 17: Right Side Panel View

**TEMPERATURE HOLDER**

(if equipped)

Store the Temperature Probe and Temperature Probe Covers in their holder locations when it is not in use.

**NOTE:**

Insure that the Temperature probe is secured in its holder. An audible "click" should be heard when the probe is completely placed into its holder.

## MONITOR OPERATING INSTRUCTIONS

### ADULT/ NEONATE OPERATING MODE

NIBP and Temperature functions are affected by changing between ADULT and NEONATE operating modes.

**IMPORTANT:**

Prior to patient monitoring, ensure the monitor is configured to the appropriate patient mode – Neonate or Adult. Refer to Page 81, SELECTING THE PATIENT MODE for more information.

Once power has been applied, a visual indicator, located on the front panel of the monitor, indicates the current operating mode.

The monitor's operating mode may also be set when the monitor power is being turned "ON".

The factory default is the Adult mode. To set the CAS 740 Monitor to the Neonate mode, depress and hold the ARROW DOWN pushbutton while turning monitor power "ON". The Neonate visual indicator will be displayed on the Main display.

To set the CAS 740 Monitor to the Adult mode, depress and hold the ARROW UP pushbutton while turning monitor power "ON". The Adult visual indicator will be displayed on the Main display.

The monitor will operate in the mode selected until it is changed.

### TURNING THE CAS 740 MONITOR "ON"

Press the ON/STANDBY pushbutton on the front panel to turn the monitor "ON".

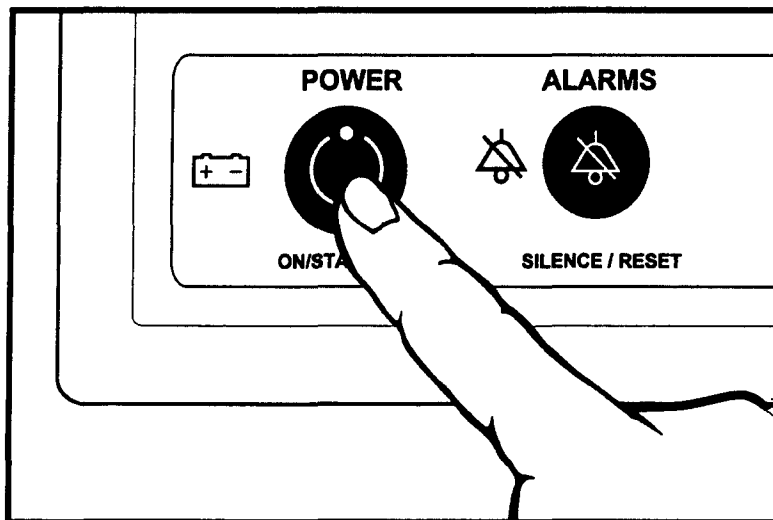


Figure 18: Turning the Monitor On



Each time the monitor is turned "ON", a one (1) second Configuration Setup Test and a four (4) second electronic Power On Self-Test (POST) is conducted to ensure that its internal circuits are functioning properly.

**NOTE:**

The user should use the Power On Self Test as a verification tool that all front panel visual indicators and the audio are functioning properly.

The one (1) second Configuration Setup Test is a visual indication of the CAS 740 Monitor's current configuration. It consists of the monitor's Model number and a description of its power source (740 or 740M), an Installed Parameter Code (1, 2, 3) and a one (1) or two (2) character module configuration code.

As an example, upon power-up the CAS 740 Monitor displays: 740-2MS.

- The first set of characters indicates the Model number and the source of its power supply (740 = Internal/AC Line Power or 740M = DC connection).
- The second character (2) describes how many parameters are installed in the monitor. The monitor in our example is configured for NIBP and SpO<sub>2</sub>. Installed Parameter Codes include;
  - (1) = NIBP; (2) = NIBP and SpO<sub>2</sub> or Temperature; (3) = NIBP, SpO<sub>2</sub> and Temperature
- The third set of characters (MS) describes the type of module installed. The monitor in our example is configured for Masimo oximeter. Configuration Codes are;
  - MS = Masimo; NL = Nellcor; NN = Nonin; T = Temperature

The four (4) second Power On Self-Test consists of:

- All equipped parameter segments are lit for one (1) second.
- All High Alarm Values are displayed for one (1) second in their corresponding numeric display window and an audible tone is emitted from the monitor's internal speaker. The Message Window indicates "HI LIMS".
- All Low Alarm Values are displayed for one (1) second in their corresponding numeric display window and an audible tone is emitted from the monitor's internal speaker. The Message Window indicates "LO LIMS".
- The monitor's current time is displayed for 1 second.

**NOTE:**

If the printer was powered "ON" during the monitor's power-up sequence, the CAS Medical Systems logo is printed.



Once the Power On Self Test is completed, the monitor indicates that it is "Ready" for use.

## FRONT PANEL INTENSITY CONTROL

The CAS 740 Monitor utilizes Super Bright LED displays in its front panel for better viewing in a variety of environments. If the displays are too bright or should the need come to operate the monitor in a dimly lit environment, the monitor incorporates a user selectable Intensity Control.

Depress and hold the front panel AUDIO pushbutton for two (2) seconds. The monitor's front panel LED displays will toggle from their bright to their dim intensity level.

Continue to depress the AUDIO pushbutton and every two (2) seconds the monitor's front panel LED displays will toggle; dim to bright, bright to dim etc. Release the AUDIO pushbutton and select the intensity level required. The CAS 740 Monitor remembers this setting and will power-up the next time to the intensity level selected.

## DISPLAYING THE TIME

The CAS 740 Monitor can be used to display the current Time of day (as set by the user) at any time. Depress and hold the CYCLE TIME pushbutton for two (2) seconds.

The monitor displays the Time, in 24 Hr. format in the Message Window for as long as the pushbutton is depressed.

## MANUAL MODE FOR BLOOD PRESSURE DETERMINATION

Select and apply the appropriate sized cuff for the patient being monitored to the extremity. Refer to Section 5, BLOOD PRESSURE MONITORING for more information.

Connect the cuff to the end of the monitor tubing. Connect the monitor tubing to the NIBP connector, located on the left side of the monitor.

To select one of the Initial Inflation Pressures shown in TABLE 4, press either the ARROW UP or the ARROW DOWN pushbutton or accept the default Initial Inflation Pressure of 150 mmHg (Adult) or 85 mmHg (Neonate).

### CAUTION:

When measuring blood pressure on a Pediatric patient, using the ADULT mode, it is recommended that the Initial Inflation Pressure be set to a lower value of 120 mmHg.

ADULT MODE	NEONATE MODE
200 mmHg	120 mmHg
180 mmHg	100 mmHg
160 mmHg	85 mmHg
150 mmHg	80 mmHg
140 mmHg	60 mmHg
120 mmHg	
100 mmHg	

Table 4: Selectable Initial Inflation Pressures

**IMPORTANT:**

Excessive patient motion can contribute to inaccurate measurements. It is important that the patient be kept still during a measurement. Make every attempt to alleviate fear, anxiety and pain.

Press the START pushbutton to begin a measurement. For the first measurement the monitor will inflate to the default setting or the Initial Inflation Pressure selected.

**NOTE:**

When the monitor is configured in the Alarms Off mode, the message "Alrm Off" is displayed in the Message Window, the cuff inflation pressure will always be set to 150 mmHg (Adult) and 85 mmHg (Neonate). Also, as an example, if the Initial Inflation Pressure was set to 100 mmHg, the first NIBP inflation pressure would be 100 mmHg followed by pressures at 150 mmHg.

The monitor's front panel NIBP displays will indicate all dashes "- -" while the measurement is in progress and the Message Window displays the Inflation Pressure in the format "Meas XXX", where "XXX" is the pressure value. For subsequent measurements, the monitor will inflate approximately 30 mmHg higher than the previously determined Systolic pressure.

The measurement typically takes less than 30 seconds to complete. In no case will the cuff remain pressurized for more than 120 seconds for Adult/Pediatric patients and no more than 90 seconds for Neonates.

When the measurement is completed, the cuff will automatically deflate and the monitor will display the Systolic, Diastolic, MAP (if enabled), and Pulse Rate values on the front panel displays. The Message Window displays a time stamp of the NIBP measurement in the format "BP HH:MM". Where "HH" is the hour and "MM" is the minute of the measurement taken.

**NOTE:**

If SpO<sub>2</sub> is active, the Pulse Rate value will be determined from the SpO<sub>2</sub> measurement.

Press the CANCEL pushbutton, at any time, to stop a measurement and deflate the cuff during the measurement process. The monitor's front panel NIBP displays will indicate all dashes "- -" and the Message Window returns to "Ready".

The CANCEL pushbutton can be depressed after a measurement is taken to clear the current reading from the displays. The monitor's front panel NIBP displays will indicate all dashes "- -" and the Message Window returns to "Ready".

**NOTE:**

If any displayed NIBP measurement were to be left on the display for up to twenty-four (24) hours, the monitor will automatically blank the displays to all dashes "- -" and the Message Window returns to "Ready".

### AUTOMATIC CYCLE FOR BLOOD PRESSURE DETERMINATION

The CAS 740 Monitor can automatically take blood pressure measurements at pre-selected time intervals.

To choose a time interval, press the **CYCLE TIME** pushbutton. The Message Window displays the format "Auto XX". Where "XX" is the time interval in minutes.

Choose the time interval using one of the following methods:

- 1) Depress the **CYCLE TIME** pushbutton to advance through the available time intervals.
- 2) Depress the **ARROW UP** or **ARROW DOWN** pushbutton keys to advance forward or backwards through the available time intervals.

Once a time interval has been selected, press the **START** pushbutton to begin the first measurement.

Between each measurement, the Message Window will display the time remaining until the next measurement, as well as the cycle time chosen in the format "MM:SS CC" where "MM" is the minutes and "SS" is the seconds until the next measurement and "CC" is the cycle time selected.

The measurement results are displayed on the front panel until the start of the next measurement cycle.

#### **NOTE:**

If a measurement is desired between measurement cycles, press the **START** pushbutton. After this measurement, the monitor will re-enter the Automatic Cycle mode and countdown to the next measurement based on the Cycle time selected.

In the Automatic Mode, while during a measurement or in the idle time between measurements, pressing the **CANCEL** pushbutton will cause the monitor to exit the Automatic mode and return to the Manual mode. The monitor's front panel NIBP numeric displays will indicate the reading of the previous measurement taken and the Message Window will display the time stamp of that reading.

### STAT MODE

#### **WARNING:**

Readings obtained during STAT mode may not meet the stated accuracy of this monitor.

#### **WARNING:**

In some cases, rapid, prolonged cycling of an oscillometric, noninvasive blood pressure monitor cuff has been associated with any or all of the following: ischemia, purpura, or neuropathy. Apply the oscillometric cuff appropriately, according to instructions, and check the cuff site and cuffed extremity regularly when blood pressure is measured at frequent intervals or over extended periods of time.

The CAS 740 Monitor can automatically take a series of blood pressure measurements for a five (5) minute interval with a brief (approx. ten (10) second) pause between determinations to allow venous blood return.

Press and Hold the START/STAT pushbutton for two (2) seconds to activate this feature. The Message Window will display "STAT XXX", where "XXX" is the real time cuff pressure. Between readings the Message Window will display "0:XX S" where "XX" is the 10-second count down until the next measurement.

After five (5) minutes of determinations the monitor will stop taking measurements and exit the STAT mode. The monitor will return to the Manual mode and the Message Window will return to the "Ready" message.

In the STAT Mode, while during a measurement or in the idle time between measurements, pressing the CANCEL pushbutton will cause the monitor to exit the STAT mode and return to the Manual mode. The monitor's front panel NIBP numeric displays will indicate the reading of the previous measurement taken and the Message Window will display the time stamp of that reading.

### NOTE:

Pressing the CYCLE pushbutton during the Manual, Automatic or STAT Modes will terminate the current NIBP measurement (if active) and place the monitor in the Automatic Cycle mode. The user will then be allowed to alter the Cycle time, if necessary, and begin a new NIBP measurement.

### WARNING:

When monitoring over an extended period of time, or at frequent intervals, periodically observe the patient's limb to make sure that the circulation is not impaired for a prolonged period of time.

## HISTORY MODES

The History Modes allow the user to recall previously taken measurements. The measurement results and the time the measurement was taken are also displayed.

Patient History is organized in two visually displayed lists.

Event History: Pressing the HISTORY pushbutton the first time reveals the first viewable list that contains up to 480 entries of Event driven data (NIBP and Predictive Temperature (if installed)) along with the corresponding %SpO<sub>2</sub> results (if installed).

Trend History: Pressing the HISTORY pushbutton again reveals the second viewable list that contains up to 480 entries of Trend History at one (1) minute intervals. The Trend History contains all the events from the Event History log combined with readings that are saved every minute consisting of: a one (1) minute average of the %SpO<sub>2</sub> values (if installed) and Monitor Mode Temperature values (if installed).

**NOTE:**

Both Event and Trend History data are only available for twenty-four (24) hours. Data older than twenty-four (24) hours is automatically removed from the History lists.

Once the HISTORY pushbutton has been depressed, use the HISTORY pushbutton to toggle between the two History menus.

When completed, press the CANCEL pushbutton to return to "Ready" or after thirty (30) seconds of button inactivity the monitor will automatically return to "Ready".

**NOTE:**

History data values for NIBP and Predictive Temperature are saved as they occur. History data values for %SpO<sub>2</sub> (one (1) minute averages) and Monitor Temperature are saved once a minute.

**NOTE:**

Turning the power "OFF" does not clear the History memory. Measurements will remain in memory for up to twenty-four (24) hours. Any measurement older than twenty-four hours is deleted. It is suggested to manually clear History between patients (See CLEARING HISTORY on Page 72).

### EVENT HISTORY (EVENT-Ev)

To review History based on Event readings, press the HISTORY pushbutton. The most current measurement (NIBP or Predictive Temperature) and the time of measurement (HH:MM-Ev) are displayed on the front panel along with their corresponding %SpO<sub>2</sub> values. If there are no readings in the History, the message "No Evnts" is displayed in the Message Window.

Press the ARROW DOWN pushbutton to review preceding measurements. The word "Oldest" will be displayed briefly on the Message Window when there are no more readings left in the memory. The oldest measurement will then be redisplayed.

Press the ARROW UP pushbutton to advance the measurements towards the most current measurement taken. The word "Newest" will be displayed briefly on the Message Window when there are no more current readings. The most recent measurement will then be redisplayed.

To exit the History Mode and return to "Ready" press the CANCEL pushbutton. If a measurement is desired immediately while in the History Mode, simply press START. This will exit the History Mode and begin a measurement.

### TREND HISTORY (TREND-Tr)

**NOTE:**

If the CAS 740 Monitor is not equipped for SpO<sub>2</sub> or Temperature monitoring, no Trend History menu will be available.

To review Trend History, press the HISTORY pushbutton twice. The most current measurement (NIBP, Predictive or Monitor Temperature or a one (1) minute averaged %SpO<sub>2</sub>) values and the time of measurement (HH:MM-Tr) are displayed on the front panel. If there are no readings in the History, the message "No Trend" is displayed in the Message Window.

Press the ARROW DOWN pushbutton to review preceding measurements. The word "Oldest" will be displayed briefly on the Message Window when there are no more readings left in the memory. The oldest measurement will then be redisplayed.

Press the ARROW UP pushbutton to advance the measurements towards the most current measurement taken. The word "Newest" will be displayed briefly on the Message Window when there are no more current readings. The most recent measurement will then be redisplayed.

To exit the History Mode and return to "Ready" press the CANCEL pushbutton. If a NIBP measurement is desired immediately while in the History Mode, simply press START. This will exit the History Mode and begin a measurement.

**NOTE:**

Turning the power "OFF" does not clear the History memory. Measurements will remain in memory for up to twenty-four (24) hours. Any measurement older than twenty-four hours is deleted. It is suggested to manually clear History between patients (See CLEARING HISTORY on Page 72).

**PRINT HISTORY**

The stored Event History and Trend History data are printed individually. Press and Hold the HISTORY pushbutton for two (2) seconds while in either the Event or Trend History screens.

The message "Print? N" appears in the Message Window.

Use either the ARROW UP or ARROW DOWN pushbutton until the message "Print? Y" appears in the Message Window.

Depress the HISTORY pushbutton again to start printing the History data. The message "Printing" appears in the Message Window as a status indicator.

**NOTE:**

Prior to accessing the Print History mode, the printer should have paper installed and the power turned "ON". Refer to Section 9, EXTERNAL PRINTER for more information about the printer and sample printouts of both History screens.

### CLEARING HISTORY

The History data can be manually cleared by pressing and holding for two (2) seconds the HISTORY pushbutton while in either of the History modes. The message "Print? N" appears in the Message Window.

Depress the HISTORY pushbutton again and the message "Erase? N" appears in the Message Window.

Use either the ARROW UP or ARROW DOWN pushbutton until the message "Erase? Y" appears in the Message Window.

Depressing the HISTORY pushbutton a final time and the message "Erasing" appears in the Message Window. The monitor erases the memory, exits the History mode and returns to the "Ready" state.

### REAL TIME CLOCK (RTC)

The CAS 740 Monitor uses an internal Real Time Clock to time stamp all entries that are stored in History. Changes made to either the time or date settings, should be performed in-between patients being monitored. Refer to Page 66, for displaying the monitors time and Page 85, for setting the time.

### PATIENT ALARM MODE

**WARNING:**

Configuring the CAS 740 Monitor's alarm settings to "OFF" will disable all audible and visual alarms. This mode should only be selected for spot check applications where the patient is receiving bedside surveillance by a trained clinician. Refer to Page 82, AUDIO ALARM SILENCE (SILENCE/RESET Pushbutton) for more information.

The CAS 740 Monitor is equipped with patient alarms to warn the user if any measurement parameter is outside the range of a user set value. This feature will allow the user to set values for:

- Systolic High and Low
- Diastolic High and Low
- MAP High and Low
- Pulse High and Low
- %SpO<sub>2</sub> High and Low



## **CHANGING ALARM LIMITS**

The CAS 740 Monitor incorporates a two (2) step process for reviewing and changing Alarm Limit values.

Press the ALARM LIMITS pushbutton once and the monitor's front panel displays indicate the current values of all High Alarm Limit values, for two (2) seconds then toggles to indicate all Low Alarm Limit values for two (2) seconds. This reviewing of all High and Low Alarm Limit values continues for thirty (30) seconds if no other pushbutton is depressed.

To change the value of an Alarm Limit, press the ALARM LIMITS pushbutton again. The Alarm Limit value is displayed in its corresponding display window and the Message Window will display the limit's name.

Use the ARROW UP or ARROW DOWN pushbuttons to change the limit value. Use the ALARM LIMITS pushbutton key to advance onto the next parameter value.

Press the CANCEL pushbutton to return to "Ready" or after thirty (30) seconds of button inactivity the monitor will automatically return to "Ready". The alarm value(s) set will now be used until power is turned "OFF".

## **SAVING ALARM LIMITS**

To save the Alarm Limit values to non-volatile memory, press and hold the ALARM LIMITS pushbutton for two (2) seconds while in the Alarm Limits set value mode. The Message Window will display "Save? N".

Press either the ARROW UP or ARROW DOWN pushbuttons will display "Save? Y" in the Message Window. Depress the ALARM LIMIT pushbutton a final time.

The Message Window will display "Saving" and the current value(s) will be saved in memory. When finished, press the CANCEL pushbutton to return to the "Ready" state. These values will be retained in memory even after the monitor is turned "OFF".

## **RESTORE FACTORY DEFAULTS**

To restore the factory default settings, press the ALARM LIMITS pushbutton to enter the Alarm Limits mode.

Press and hold the ALARM LIMITS pushbutton for 2 seconds until the message "Saving? N" appears in the Message Window.

Press the ALARM LIMITS pushbutton to bring up the Recall menu. The message "Recall? N" appears in the Message Window. Pressing either the ARROW UP or ARROW DOWN pushbutton will display the message "Recall? Y" in the Message Window.

Press the ALARM LIMITS pushbutton. The Message Window will display "Restore" and return to the Alarm Limit menu. When finished, press the CANCEL pushbutton to return to the "Ready" state.

## ALARM LIMIT VALUES

There are two sets of alarm limit parameters, one for Adult/Pediatric and one for Neonatal. The Alarm Limits pushbutton will operate on the parameters for the mode the monitor is currently in (ADULT / NEO). Table 5, lists the Alarm Limit Default Values used by the CAS 740 Monitor.

<b>ADULT</b>	<b>Systolic</b>	<b>Diastolic</b>	<b>MAP</b>	<b>Pulse</b>	<b>%SpO<sub>2</sub></b>
High	240	130	180	220	OFF
Low	OFF	OFF	OFF	OFF	88
<b>NEONATAL</b>	<b>Systolic</b>	<b>Diastolic</b>	<b>MAP</b>	<b>Pulse</b>	<b>%SpO<sub>2</sub></b>
High	120	80	100	220	96
Low	OFF	OFF	OFF	OFF	88

Table 5: Default Alarm Values

**NOTE:**

Switching modes from Adult to Neonate or Neonate to Adult, shall recall the last stored values for patient alarm limit.

**NOTE:**

Patient alarms for Systolic, Diastolic and MAP values are produced at the time the measurement is taken. Alarms for %SpO<sub>2</sub> have a ten (10) second delay. Alarm for Pulse Rate will be immediate if taken from the NIBP or have a ten (10) second delay if taken from the SpO<sub>2</sub> signal.

**AUDIBLE AND VISUAL INDICATORS**

The CAS 740 Monitor is capable of producing both an audible and a visual indicator for a variety of monitor conditions. The following table provides a cross reference for audible and visual indications.

Alarm Condition	Priority Level	Audible Indication	Visual Indication
Dead Battery	High	3 Beeps followed by 2 Beeps every 10 seconds	The message is displayed continuously in the Message Window. The Battery Power Visual Indicator is Red.
Patient Alarm Limit Violations	High	3 Beeps followed by 2 Beeps every 10 seconds	The associated 7-segment LED flashes for one second on / one second off for the parameter limit in violation.
SpO <sub>2</sub> % Alarm Limit Violations			The message is displayed for one second every two to six seconds in the Message Window *.
Power Fail	High	3 Beeps followed by 2 Beeps every 10 seconds	The message is displayed for one second every two to four seconds in the Message Window *.
Low Battery Alarm	Medium	3 Beeps every 25 seconds	The message is displayed for one second every two to six seconds in the Message Window *. The Battery Power Visual Indicator is Red.
NIBP Application Error	Medium	3 Beeps every 25 seconds	The message is displayed for one second every two to six seconds in the Message Window *.
SpO <sub>2</sub> Probe Alarms	Medium	3 Beeps every 25 seconds	The message is displayed for one second every two to six seconds in the Message Window *.
NIBP Complete	Low	1 Beep when blood pressure is completed	Numerics are updated with each blood pressure taken.
SpO <sub>2</sub> Pulse Beat	Low	1 Beep coincides with each received SpO <sub>2</sub> pulse rate	Numerics are updated with each received SpO <sub>2</sub> pulse rate value.
Temperature Complete	Low	1 Beep when final temperature is completed	Numerics are updated when a final temperature is completed.
Key Click	Low	1 Beep associated with each button action	None

**Table 6: Audible and Visual Indicators**

\* - The message interval time will vary based on what monitor conditions are present.

**NOTE:**

Refer to Table 9, on Pages 93 for a listing of messages that may be displayed in the Message Window.

### **CLEARING ALARMS**

The CAS 740 Monitor provides to the user an audible and visual indication for both patient and equipment alarm conditions.

### **NIBP PATIENT ALARMS**

During an active NIBP Patient Alarm Limit Violation (Low/High Systolic or Low/High Diastolic), the monitor flashes the LED display of the parameter in alarm and provides a Patient Alarm Limit Violation audible tone.

To acknowledge the alarm, depress the SILENCE/RESET pushbutton.

The monitor silences the audible tone and the parameter in alarm stops flashing.

### **HIGH / LOW %SpO<sub>2</sub> ALARMS**

During either a High or Low %SpO<sub>2</sub> Alarm Limit Violation, the monitor flashes the %SpO<sub>2</sub> LED display, provides a Patient Alarm Limit Violation audible tone and displays the appropriate message (SpO<sub>2</sub> Hi or SpO<sub>2</sub> Lo) in the Message Window.

To acknowledge the alarm and temporarily mute the audio, depress the SILENCE/RESET pushbutton.

The monitor silences the audible tone for a maximum of fifteen (15) seconds, the parameter in alarm continues to flash on the front panel and the appropriate message will be displayed in the Message Window. If the current patient alarm condition becomes inactive during the 15 seconds and then recurs OR a different patient alarm condition occurs during the 15 seconds, the audio alarm is reactivated. If the alarm condition continues uninterrupted for the 15 seconds, the audio alarm is reactivated.

### **SpO<sub>2</sub> PULSE RATE ALARM**

During a SpO<sub>2</sub> Pulse Rate Alarm Limit Violation, the monitor flashes the ♥BPM LED display and provides a Patient Alarm Limit Violation audible tone.

To acknowledge the alarm and temporarily mute the audio, depress the SILENCE/RESET pushbutton.

The monitor silences the audible tone for a maximum of fifteen (15) seconds and the parameter in alarm continues to flash on the front panel. If the current patient alarm condition becomes inactive during the 15 seconds and then recurs OR a different patient alarm condition occurs during the 15 seconds, the audio alarm is reactivated. If the alarm condition continues uninterrupted for the 15 seconds, the audio alarm is reactivated.


## EQUIPMENT ALARMS

During an Equipment Alarm, the monitor displays the alarm parameter in either the Message Window or flashes the LED display of the alarmed parameter and provides an audible Alarm Tone. To clear the alarm, depress the SILENCE/RESET pushbutton. The monitor silences the audible tone and clears the parameter.

**NOTE:**

The Low Battery and Dead Battery alarms cannot be silenced.

## ADJUSTING THE AUDIO ALARM VOLUME

The Alarm Volume can be adjusted to one (1) of five (5) volume levels. Depress the AUDIO pushbutton. The Message Window displays the current value of the Audio Alarm Volume "ALARM .

Use the ARROW UP and ARROW DOWN pushbuttons to adjust the volume level desired. Altering the Alarm Volume will produce a single tone at the selected volume level.

Press the CANCEL pushbutton when completed.

**NOTE:**


The Volume Level cannot be set to "OFF".

**NOTE:**

No "Key Click" will be heard when the Audio Volume is set to a MINIMUM setting.

## ADJUSTING THE SpO<sub>2</sub> "BEEP" VOLUME

(available if SpO<sub>2</sub> is installed)

The Volume for the SpO<sub>2</sub> "Beep" can be adjusted to one (1) of five (5) levels and "OFF". This setup menu follows directly after the Audio Alarm Volume. Depress the AUDIO pushbutton until the Message Window displays the SpO<sub>2</sub> Volume message "SpO<sub>2</sub> .

Use the ARROW UP and ARROW DOWN pushbuttons to adjust the volume level desired. Altering the SpO<sub>2</sub> Volume will produce a single tone at a selected "non-off" volume level.

Press the CANCEL pushbutton when completed.

## 2-MINUTE AUDIO SILENCE

**NOTE:**

Enabling the 2-Minute Audio Alarm Off mode can only be accomplished after all active alarm conditions have been addressed. To clear an alarm, refer to Page 76, CLEARING ALARMS for more information.

If the 2-Minute Audio Off mode is selected in the Monitor Configuration menu, pressing the SILENCE/RESET pushbutton will temporarily deactivate the audio alarms for two (2) minutes. The visual indicator next to the SILENCE/RESET pushbutton will illuminate continuously and the message "2Min Aud" will be displayed along with the "Ready" message in the Message Window as a reminder.

**NOTE:**

If the Message Window is currently displaying a time stamp from a previously taken NIBP measurement, the "2Min Aud" message will alternately be displayed with the time stamp message.

When the 2-Minute Audio Off is enabled, the monitor will not alarm for patient related High and Low alarms, NIBP Application Error and SpO<sub>2</sub> Probe Alarm. The monitor will alarm for Low Battery and Dead Battery.

To re-activate the audio alarm, press the SILENCE/RESET pushbutton again. The visual indicators will go out. If the button is not pressed, audible alarms will be re-armed automatically after two (2) minutes. Refer to Page 82, AUDIO ALARM SILENCE (SILENCE/RESET Pushbutton) for more information.

### PERMANENT AUDIO ALARM SILENCE

**NOTE:**

Enabling the Permanent Audio Alarm Off mode can only be accomplished after all active alarm conditions have been addressed. To clear an alarm, refer to Page 76, CLEARING ALARMS for more information.

If the Permanent Audio Off mode is selected in the Monitor Configuration menu, pressing the SILENCE/RESET pushbutton will deactivate the audio alarms. The visual indicator next to the SILENCE/RESET pushbutton will flash at a rate of one (1) second "ON" and one (1) second "OFF" and the message "Perm Aud" will be displayed along with the "Ready" message in the Message Window as a reminder.

**NOTE:**

If the Message Window is currently displaying a time stamp from a previously taken NIBP measurement, the "Perm Aud" message will alternately be displayed with the time stamp message.

When enabled, the monitor will not alarm for patient related High and Low alarms, NIBP Application Error and SpO<sub>2</sub> Probe Alarm. The monitor will alarm for Low Battery and Dead Battery.

To re-activate the audio alarm, press the SILENCE/RESET pushbutton again. The visual indicators will go out. Refer to Page 82 AUDIO ALARM SILENCE (SILENCE/RESET Pushbutton) for more information.

## ALARM LIMITS OFF

If the Alarm Limits Off mode is selected in the Monitor Configuration menu, all alarms associated with patient alarms are "OFF" except for Low Battery and Dead Battery. The visual indicator next to the SILENCE/RESET pushbutton will flash at a rate of two (2) second "ON" and two (2) seconds "OFF" and the message "Alrm Off" will be displayed along with the "Ready" message in the Message Window as a reminder.

When the Alarm Limits Off is selected, the ALARM LIMITS pushbutton is inactive. The Alarm Limits Off mode cannot be disabled with the SILENCE/RESET pushbutton.

### NOTE:

If the Message Window is currently displaying a time stamp from a previously taken NIBP measurement, the "Alrm Off" message will alternately be displayed with the time stamp message.




To re-activate the audio alarms, refer to Page 82, AUDIO ALARM SILENCE (SILENCE/RESET Pushbutton) for more information.

## MONITOR CONFIGURATION

The Monitor Configuration section allows the user to configure the CAS 740 Monitor to your individual needs. Once entered, the user can:

- Review the monitor's internal Software Revisions
- Select the Operating Language
- Select the Patient Mode
- Select the Power-Off Delay Time
- Select the Temperature Scale
- Make selections for Audio Alarms
- Choose to display the MAP value
- Set the SpO<sub>2</sub> Alarm Delay
- Set the Date
- Set the Time
- Set Daylight Saving Time Option
- Perform System Checks (Refer to Page 116, CALIBRATION CHECK)
  - Manometer Mode
  - Leak Check

## ENTERING THE CONFIGURATION MENU

To enter the monitor's Configuration Menu, depress and hold the AUDIO  and ALARM LIMITS  pushbutton keys while the monitor is being turned "ON" .

Once in the menu, use one of the NEXT (HISTORY / AUDIO) or PREVIOUS (CYCLE TIME / ALARM LIMITS) programmed pushbutton keys to advance onto the next or go back to the previous parameter in the Configuration Menu.

**NOTE:**

While in the Configuration Menu, if no pushbutton is depressed within 60 seconds, the monitor will automatically save all changes made and exit the Monitor Configuration menu. The Message Window will briefly display "Saving" and return to the "Ready" mode.

## SAVING YOUR CHANGES

When you have completed configuring the monitor, press the CANCEL pushbutton to exit and lock in your selection(s). The Message Window will briefly display "Saving" and return to the "Ready" mode.

## SOFTWARE REVISIONS

The CAS 740 Monitor displays the current software revision of its operating system and that of the internal modules being used inside. The software versions are displayed in the following order:

Software Module	Message Window
CAS 740 Control Board	Ver X.X
Boot Loader	Boot X.XX
Power Supply PIC Processor	PIC X.X
CAS NIBP Module	ND X.X
Nellcor SpO <sub>2</sub> Module <sup>(1)</sup>	NEL X.X
Masimo SpO <sub>2</sub> Module <sup>(1)</sup>	MAS X.X
Nonin SpO <sub>2</sub> Module <sup>(1)</sup>	NON XX
Welch Allyn Temperature Module <sup>(2)</sup>	WA X.X

Table 7: Software Revisions

Use the ARROW UP or ARROW DOWN pushbuttons to view the messages.

Press one of the NEXT programmed pushbutton keys to advance onto the next menu or use the CANCEL pushbutton to Save your changes and exit to the "Ready" mode.

(1) The SpO<sub>2</sub> module is optional, in the case when it is not installed the Version text advances to the next Module.

(2) The Temperature module is optional, in the case when it is not installed the Version text advances to the next Module.



**SELECTING THE LANGUAGE**

The CAS 740 Monitor can operate in one (1) of nine (9) languages: English, German, French, Italian, Spanish, Dutch, Swedish, Portuguese or Norwegian.

To configure the monitor's operating language, first enter the Configuration menu. Refer to Page 79, ENTERING THE CONFIGURATION MENU.

Once in the menu, use one of the NEXT programmed pushbutton keys until the Message Window displays the current language being used.

Use the ARROW UP or ARROW DOWN pushbuttons to make your selection.

Press one of the NEXT programmed pushbutton keys to advance onto the next menu or use the CANCEL pushbutton to exit to the "Ready" mode.

**SELECTING THE PATIENT MODE**

The CAS 740 Monitor can be used on patients from Neonates to Adults.

To configure the monitor's operating mode, first enter the Configuration Menu. Refer to Page 79, ENTERING THE CONFIGURATION MENU.

Once in the menu, use one of the NEXT programmed pushbutton keys until the Message Window displays "Patient".

Use the ARROW UP or ARROW DOWN pushbuttons to select the patient mode. The front panel display will illuminate with the patient mode selected (ADULT = Adult / Pediatric; NEO = Neonate).

**NOTE:**

Switching modes from Adult to Neonate or Neonate to Adult, shall recall the last stored values for patient alarm limit.

Press one of the NEXT programmed pushbutton keys to advance onto the next menu or use the CANCEL pushbutton to exit to the "Ready" mode.

**NOTE:**

The Patient Mode may also be changed at monitor power-up. Refer to Page 64, ADULT/NEONATE OPERATING MODE for more information.

### SELECTING THE POWER-OFF DELAY TIME

The CAS 740 Monitor incorporates a user selectable Power-Off Delay time feature. The amount of time the ON/STANDBY pushbutton is depressed to turn the display "OFF" can be configured to be either "Off 0S" (no delay) or "Off 2S" (depress and hold for 2 seconds). The default setting used by the monitor is "Off 0S".

To configure the Power-Off Delay Time, first enter the Configuration menu. Refer to Page 79, ENTERING THE CONFIGURATION MENU.

Once in the menu, use one of the NEXT programmed pushbutton keys until the Message Window displays the monitor's current setting for the Power-Off Delay Time.

Use the ARROW UP or ARROW DOWN pushbuttons to make your selection.

Press one of the NEXT programmed pushbutton keys to advance onto the next menu or use the CANCEL pushbutton to exit to the "Ready" mode.

### SELECTING THE TEMPERATURE SCALE

(available if Temperature is installed)

The CAS 740 Monitor can display Temperature readings in either the Celsius or Fahrenheit scales.

To select the operating Temperature scale, first enter the Configuration menu. Refer to Page 79, ENTERING THE CONFIGURATION MENU.

Once in the menu, use one of the NEXT programmed pushbutton keys until the Message Window displays the Temperature setup menu "°F" or "°C".

Use the ARROW UP or ARROW DOWN pushbuttons to make your selection.

Press one of the NEXT programmed pushbutton keys to advance onto the next menu or use the CANCEL pushbutton to exit to the "Ready" mode.

### AUDIO ALARM SILENCE (SILENCE/RESET Pushbutton)

The CAS 740 Monitor's SILENCE/RESET pushbutton can be configured to have the audio associated with patient alarms set to one of three selections. The selections are:

- 2-Minute Audio Alarm Silence (Default)
- Permanent Audio Alarm Silence
- Alarm Limits Off

To configure the alarms, first enter the Configuration menu. Refer to Page 79, ENTERING THE CONFIGURATION MENU.

Once in the menu, use one of the NEXT programmed pushbutton keys until the Message Window displays the monitor's current value of the Audio Alarms menus.

Use the ARROW UP or ARROW DOWN pushbuttons to make your selection.

Press one of the NEXT programmed pushbutton keys to advance onto the next menu or use the CANCEL pushbutton to exit to the "Ready" mode.

### **2-MINUTE AUDIO ALARM SILENCE**

When the monitor is configured for the 2-Minute Audio Alarm Silence setting, use the SILENCE/RESET pushbutton to "enable or disable" audio alarms for a two (2) minute period. The SILENCE visual indicator, located on the front panel of the monitor will be illuminated constantly and the message "2Min Aud" will be displayed on the Message Window as a reminder when enabled. At the end of two (2) minutes, the monitor will automatically exit the 2-Minute Audio Alarm Silence setting and return to normal operation.

During a two-minute silence period, if an alarm (patient or equipment) occurs, except for Low Battery and Dead Battery, the audio alarm remains silenced for the remainder of the two-minutes and only a visual indicator is provided.

### **PERMANENT AUDIO ALARM SILENCE**

When the monitor is configured to the Permanent Audio Alarm Silence setting, use the SILENCE/RESET pushbutton to "enable or disable" audio alarms. The SILENCE visual indicator, located on the front panel of the monitor will flash at a rate of one (1) second "ON" and one (1) second "OFF" and the message "Perm Aud" is displayed on the Message Window as a reminder when enabled.

During a permanent audio alarm off period, if an alarm (patient or equipment) occurs, except for Low Battery and Dead Battery, the audio alarm remains silenced and only a visual indicator is provided.

### **ALARM LIMITS OFF**

When the monitor is configured to the Alarm Limit Off setting, all alarms associated with patient alarms are "OFF". Also, when a NIBP measurement is started using the Manual mode, the inflation pressure will always be set to 150 mmHg (Adult mode) and 85 mmHg (Neonate mode).

This mode is useful for spot check applications or if the monitor is being moved from patient to patient and the user may not want to be disturbed by any audible alarms or have to worry about what pressure to set the monitor to.

The SILENCE visual indicator, located on the front panel of the monitor will flash at a rate of two (2) seconds "ON" and two (2) seconds "OFF" and the message "Alrm Off" will be displayed on the Message Window as a reminder when enabled.

During an alarm limit off period, if an equipment alarm occurs, except for Low Battery and Dead Battery, the audio alarm remains silenced and only a visual indicator is provided.

### MAP VALUE ENABLE / DISABLE

During a blood pressure reading, the user can elect to display or not to display the MAP value.

Should it need to be changed, first enter the Configuration menu. Refer to Page 79, ENTERING THE CONFIGURATION MENU.

Once in the menu, use one of the NEXT programmed pushbutton keys until the Message Window displays the monitor's current setting for the MAP value "MAP On" or "MAP Off".

Use the ARROW UP or ARROW DOWN pushbuttons to make your selection.

**NOTE:**

When "MAP Off" is selected, MAP values are omitted from History Display and Printing as well. All alarms associated with MAP values are also disabled.

Press one of the NEXT programmed pushbutton keys to advance onto the next menu or use the CANCEL pushbutton to exit to the "Ready" mode.

### SET THE SpO<sub>2</sub> ALARM DELAY

The delay time until an alarm is generated for %SpO<sub>2</sub> and Pulse Rate can be configured to be either zero (0) seconds "Delay 0S" or ten (10) seconds "Delay10S". The default value used by the Model 740 Monitor is ten (10) seconds.

Should it need to be changed, first enter the Configuration menu. Refer to Page 79, ENTERING THE CONFIGURATION MENU.

Once in the menu, use one of the NEXT programmed pushbutton keys until the Message Window displays the monitor's current setting for the SpO<sub>2</sub> Alarm Delay.

Use the ARROW UP or ARROW DOWN pushbuttons to make your selection.

Press one of the NEXT programmed pushbutton keys to advance onto the next menu or use the CANCEL pushbutton to exit to the "Ready" mode.

**SETTING THE DATE**

The CAS 740 Monitor's Date value is set at the factory. Should it need to be changed, first enter the Configuration menu. Refer to Page 79, ENTERING THE CONFIGURATION MENU.

Once in the menu, use one of the NEXT programmed pushbutton keys until the Message Window displays the monitor's date using the following format: "DDMMMYY". Where DD = Day of the Month, MMM = Month of the Year (JAN, FEB, etc.) and YY = Last 2 digits of the year (2002 is displayed as 02). The flashing parameter indicates the parameter that can be changed.

Use the ARROW UP or ARROW DOWN pushbuttons to make your selection. Press one of the NEXT programmed pushbutton keys to advance to the next parameter to set within the Date menu.

Press one of the NEXT programmed pushbutton keys to advance onto the next menu or use the CANCEL pushbutton to exit to the "Ready" mode.

**SETTING THE TIME**

The CAS 740 Monitor's Time value is set for Eastern Time and is set at the factory. Should it need to be changed, first enter the Configuration menu. Refer to Page 79, ENTERING THE CONFIGURATION MENU.

Once in the menu, use one of the NEXT programmed pushbutton keys until the Message Window displays the monitor's time using the following format: "TM HH:MM". Where HH = Hour of the Day (0 – 23) and MM = Minute of the Hour (0 – 59). The flashing parameter indicated the parameter that can be changed.

Use the ARROW UP or ARROW DOWN pushbuttons to make your selection. Press one of the NEXT programmed pushbutton keys to advance to the next parameter to set within the Time menu.

**NOTE:**

Altering the Date and Time will affect the History readings, but not erase them.

Press one of the NEXT programmed pushbutton keys to advance onto the next menu or use the CANCEL pushbutton to exit to the "Ready" mode.

### DAYLIGHT SAVING TIME OPTION

The CAS 740 Monitor can be configured to automatically respond to time changes associated with Daylight Saving Time. The monitor can be configured to one of five Daylight Saving Time Option settings. They are:

- **DST OFF** Daylight Saving Time is "OFF". The user is responsible for changing the time if needed. This is the default setting for the CAS 740 Monitor.
- **DST N AM** Daylight Saving Time "North America". Use this setting and the monitor will automatically *add* one (1) hour the first Sunday in April at 2 a.m. and *subtract* (1) hour the last Sunday in October at 2 a.m.
- **DST EU 1, 2, 3** Daylight Saving Time "European Union".

In the European Union, Daylight Saving Time begins and ends at 1 a.m. Universal Time (Greenwich Mean Time). It starts the last Sunday in March, and ends the last Sunday in October. In the EU, all time zones change at the same moment.

Select EU 1 if the monitor will be located in Ireland, Portugal or the United Kingdom. Select EU 3 for Finland. EU 2 can be used for all remaining countries within the European Union.

**NOTE:**

Enabling Daylight Saving Time will affect the History readings, but not erase them.

Should it need to be changed, first enter the Configuration menu. Refer to Page 79, ENTERING THE CONFIGURATION MENU.

Once in the menu, use one of the NEXT programmed pushbutton keys until the Message Window displays the monitor's current setting for Daylight Savings Time.

Press one of the NEXT programmed pushbutton keys to advance onto the next menu or use the CANCEL pushbutton to exit to the "Ready" mode.

### BATTERY POWER

The CAS 740 Monitor is equipped with an internal rechargeable battery. The battery is charging whenever the monitor is plugged into a power source (AC Line Power or +12 VDC). A Battery Power Visual Indicator, located on the front panel, indicates the status condition of the monitor's battery.

Batteries will self-discharge when they are not used. It is recommended that the battery be maintained at full charge by leaving the monitor connected to a power source whenever possible.

The standard 7.2 Volt 3700 mAhr battery pack, when fully charged, is capable of taking 100 NIBP readings when the monitor is set in the 5-minute Automatic Mode.

When the message "Low Batt" appears in the Message Window, at least thirty (30) minutes of battery operation remain. The "Low Batt" message will alternate continuously with the "Ready" message to indicate that the battery should be charged as soon as possible. Also, the monitor's front panel Battery Power Visual Indicator will change from Orange to Red and three (3) audio "beeps" are heard every twenty-five (25) seconds.

**WARNING:**

Upon the detection of a Low Battery condition and if the battery is not charged by the user, the monitor may no longer function as intended. The monitor should be plugged into a power source as soon as possible and the battery allowed to charge for four (4) hours.

When the "Dead Bat" message appears, the battery is no longer able to power a measurement. The message "Dead Bat" is displayed continuously in the Message Window, the Battery Power Visual Indicator is colored Red and three (3) audio "beeps" followed by two (2) audio "beeps" once every ten (10) seconds are heard until the power is turned off.

**WARNING:**

Upon the detection of a Dead Battery condition and if the monitor is not turned off by the user, the monitor shuts down and turns "OFF" after three (3) minutes of operations.

When either of these messages appears, it is necessary to recharge the battery. A depleted battery may be fully recharged in four (4) hours. The monitor can be used to obtain measurements while the battery is charging.

**NOTE:**

Using the monitor while charging may lengthen the time to restore battery charge.

**NOTE:**

During charging of the battery, the case may feel warm to the touch.

## CHECKING BATTERY STATUS

### NOTE:

If your CAS 740 Monitor is not equipped with a Smart Pack battery, the monitor, when connected to a power source (AC Line Power or +12 VDC) and is powered "OFF", will display "Charging" in the Message Window until the battery reaches a full charge condition. Once the battery reaches full charge, the monitor is capable of supplying battery charge information as described in items 1, 2 and 3 below.

If your CAS 740 Monitor is equipped with a Smart Pack battery, the following applies:

The monitor's Main Board receives information from electronics enclosed within the battery pack. The status of the battery pack can be verified using one of the following methods:

- 1) When the monitor is connected to a power source (AC Line Power or +12 VDC) and is powered "ON", depress and hold the SILENCE/RESET pushbutton for two (2) seconds. The monitor will display in the Message Window, for as long as the button is depressed, "XXX % CHG", where "XXX" represents the percent (0 to 100) of the full charge condition of the battery.
- 2) When the monitor is connected to a power source (AC Line Power or +12 VDC) and is powered "OFF", the Message Window will display "XXX % CHG", where "XXX" represents the percent (0 to 100) of the full charge condition of the battery.
- 3) When the monitor is running on battery power and is powered "ON", depress and hold the SILENCE/RESET pushbutton for two (2) seconds. Based on the revision of the Control Board software installed in the monitor, one of the following will apply;
  - a) If your CAS 740 Monitor has Control Board software Version 1.1, 1.2, or 1.3 installed, the Message Window will display the "Estimated Battery Run Time" remaining for as long as the button is depressed in one of formats listed below;
    - If less than one-half hour is available, the Message Window will display "<.5 hrs".
    - If greater than one-half hour is available, the Message Window will display time in hours (.5, 1, 1.5, 2, 2.5 etc).
    - If more than ten (10) hours is available, the Message Window will display ">10 hrs".
  - b) If your CAS 740 Monitor has Control Board software Version 1.4 or greater installed, the Message Window will display "XXX % CHG", where "XXX" represents the percent (0 to 100) of the full charge condition of the battery for as long as the button is depressed.
- 4) The CAS 740 Monitor's Battery Pack may contain on-board electronics, which can be used to check the status of battery charge remaining when not connected to the monitor. If available, depress the pushbutton located on the inside panel of the battery pack. Once depressed, four (4) LED indicators will illuminate to display the status of battery charge remaining (each LED lit is equivalent to approximately 25%).



### **CAUTION:**

This product contains a rechargeable battery that is recyclable. Under various state and local laws, it may be illegal to dispose of this battery into the municipal waste stream. Check with your local authorities for instructions on recycling options in your area.

### **AUTO OFF FEATURE**

The CAS 740 Monitor incorporates an Auto "OFF" feature to preserve battery life. The Auto Off feature is enabled automatically when the user selects the "Alrm Off" option in the Configuration Menu.

When the monitor is in this mode and if accidentally left powered "ON" for ten (10) minutes without any key presses or any parameters active, the monitor will beep for five (5) seconds prior to automatically shutting itself down to preserve battery life.

For the Auto "OFF" feature to work properly, the following conditions must apply;

- "Alrm Off" is selected in the Configuration Menu.
- The CAS 740 Monitor is operating on battery power.
- No NIBP, SpO<sub>2</sub> or Temperature readings were taken in the past ten (10) minutes.
- No key presses were made in the past ten (10) minutes.
- The monitor cannot be in the Automatic NIBP cycle mode.

### **POWER FAIL**

The CAS 740 Monitor incorporates a Power Fail feature. Whenever the power is disconnected from the monitor and the monitor is not allowed to shut down in an orderly fashion, the monitor, when re-powered alerts the user. The message "Pwr Fail" is displayed in the Message Window and three (3) audio "beeps" followed by two (2) audio "beeps" are heard every ten (10) seconds. During this condition, all other pushbuttons are inactive except for the ON/STANDBY and SILENCE/RESET.

To clear the Power Fail condition, and return to "Ready", depress the SILENCE/RESET pushbutton or recycle the monitor's power.

### **USER MESSAGES**

The CAS 740 Monitor displays a variety of messages to aid the user in monitor operation. If a troubleshooting message is displayed during a measurement, follow the actions listed to correct the situation.

If the monitor does not turn on, or exhibits a flashing display and failure to operate, the battery is most likely below the Dead Battery point. Connect the monitor to a power source (AC Line Power or +12 VDC) and allow it to charge for four (4) hours.

## CAS 740 Monitors

Refer to Table 9: Error Messages on the Message Window for more information.

If the monitor is in need of repair, it must be referred to the appropriate service personnel. Service performed by unauthorized personnel could be detrimental to the monitor and will void the warranty. For service, contact your dealer or CAS Medical Systems, Inc.

### SpO<sub>2</sub> USER MESSAGES

(available if SpO<sub>2</sub> is installed)

If the accuracy of any measurement does not seem reasonable, first check the patient's vital signs by an alternate method.

#### NOTE:

The SpO<sub>2</sub> probe must be kept as motionless as possible to make a proper determination. Use the SpO<sub>2</sub> strength bar graph to determine if a strong rhythmic pulse signal is present.

When no oximeter probe is attached to the monitor, the %SpO<sub>2</sub> window and signal strength window will be blank. When no SpO<sub>2</sub> pulse data is available, the monitor will display the last NIBP pulse.

When the probe is connected to the monitor, but is off of the patient, the message “-” is displayed in the %SpO<sub>2</sub> and Pulse Rate windows. The Message Window flashes the message “Prb OFF” and three (3) audio “beeps” are heard every twenty-five (25) seconds.

Depress the SILENCE/RESET pushbutton. The monitor silences the audible alarm tone, but the message remains.

If the message “Prb” should appear in the %SpO<sub>2</sub> window, verify that the probe being used is the correct one for the monitor's SpO<sub>2</sub> configuration (Masimo, Nellcor, Nonin) or that the probe is not defective.

Depress the SILENCE/RESET pushbutton. The monitor silences the audio alarm tone, but the message remains. Remove the defective probe and replace it with a working probe.

Each sensor is designed for a specific clinical application.

#### NOTE:

Inaccurate measurements may be caused by:

- anemia or low hemoglobin concentrations
- electro-surgical interference
- excessive ambient light
- excessive patient movement
- incorrect sensor application or use
- intravascular dyes such as indocyanine green or methylene blue
- moisture in the sensor
- placement of a sensor on an extremity with a blood pressure cuff, arterial catheter, or intravascular line
- venous pulsations

**NOTE:**

The loss of a pulse signal can occur in any of the following situations:

- a blood pressure cuff is inflated on the same extremity as the one with the SpO<sub>2</sub> sensor attached
- excessive ambient light such as from a surgical lamp, a bilirubin lamp, or sunlight
- the patient has hypotension, severe vasoconstriction, severe anemia, or hypothermia
- the patient is in cardiac arrest or is in shock
- the sensor is too tight
- there is arterial occlusion proximal to the sensor

If the SpO<sub>2</sub> Module located inside the CAS 740 Monitor should fail, the message "Err" will appear in the %SpO<sub>2</sub> display window.

Depress the SILENCE/RESET pushbutton. The monitor silences the audio alarm tone, but the message remains.

Should any of the above problems persist, contact your dealer or CAS Medical Systems, Inc.

**TEMPERATURE FUNCTION USER MESSAGES**

If the probe becomes unattached to the CAS 740 Monitor's rear panel connector, the message "Prb" is shown in the TEMP display window.

**NOTE:**

Depress the SILENCE/RESET pushbutton to clear the TEMP display.

When an attached probe is removed from the probe holder, the message OrL, ALy or rEC is shown briefly in the TEMP display window indicating the predictive algorithm is being used by the Temperature function.

When the message "Err XX", where "XX" is an Error Number, is displayed in the TEMP display window an error condition has occurred and a reliable temperature reading could not be obtained.

**NOTE:**  
Depress the SILENCE/RESET pushbutton to clear the TEMP display.

Refer to Table 8: Temperature Error Codes for more information.

Verify the monitor's operating environment are within its limits and start the procedure from the beginning.

ERROR NUMBER	ERROR DESCRIPTION
00	Transmit buffer overflow.
01	Probe heater energy accumulation too high.
02	Probe a/d pulse width out of range.
03	Adaptive probe gain too high or too low.
11	Ambient temperature above 104 °F.
12	Ambient temperature below 60.8 °F.
21	Battery voltage below error threshold of 3.0 volts.
31	RAM read/write error.
32	ROM checksum error.
33	CPU instruction error.
40	PTB resistor a/d pulse width out of range.
41	RatioCal resistor a/d pulse width out of range.
42	External ambient thermistor a/d pulse width out of range.
50	Heater circuit failure.
51	Probe heated above 112 °F.
52	Heater watchdog timeout failure.
60	PTB resistor "temperature" out of range.
99	Temperature option no longer recognized.

**Table 8: Temperature Error Codes**

Should any of these problems persist, contact your dealer or CAS Medical Systems, Inc.

## ERROR MESSAGES ON THE MESSAGE WINDOW

ERROR MESSAGE	POSSIBLE CAUSE	POSSIBLE SOLUTION
"Air Leak"	Air leak in cuff/hose/monitor pneumatic system.	Check that the cuff/hose/monitor connection is secure.  Check cuff for leaks. Do not use a known leaky cuff.
"Appl Err"	Neonate cuff is detected in Adult Mode.	Check cuff. Replace cuff or change operating mode
"Chk Prb" (Masimo)	The monitor is questioning the quality of the signal being received by the SpO <sub>2</sub> sensor.  The sensor is receiving too much ambient light.	Verify that the sensor is being used according to the manufacturer's recommendations.  Verify that the sensor emitter and detector are parallel to and directly opposing each other.
"ChksumEr"	An electronic failure has occurred within the monitors' Main Control Board.	Contact CAS Medical Systems to have the monitor serviced.
"Dead Bat"	The battery is fully discharged.	Recharge the battery for at least 4 hours.
"Flow Err"	Stable cuff pressure cannot be maintained by the pneumatic system.	Check the external tube for kinks.  Perform a Pneumatic Check as detailed in the Maintenance section of this manual.  Replace cuff.
"LooseCuf"	Cuff applied too loosely.	Check cuff for proper fit on patient.
"Low Batt"	The battery is almost discharged.	At least 30 minutes of operation is available from when the message first appears.  Recharge the battery as soon as possible.
"Low Perf" (Masimo)	The perfusion level being received by the SpO <sub>2</sub> sensor is low.	Determine if an extreme change in the patient's physiology and blood flow at the monitoring site occurred (e.g. an inflated blood pressure cuff, a squeezing motion).  Try to warm the patient or sensor site.  Move sensor to a site with better perfusion.

Table 9: Error Messages on the Message Window

**ERROR MESSAGES ON THE MESSAGE WINDOW (cont.)**

<b>ERROR MESSAGE</b>	<b>POSSIBLE CAUSE</b>	<b>POSSIBLE SOLUTION</b>
"Motion"	There was too much extremity motion for the monitor to accurately complete the NIBP measurement in 120 seconds.	Measurements can be obtained when there is limited extremity movement, but the measurement time may be extended.  Measurement time is limited to 120 seconds.  Restrain patient extremity motion.
"NBP Cal"	Pressure calibration data corrupted within NIBP module.	Pressure module needs recalibration.  Contact CAS Medical Systems to have the monitor serviced.
"NIBP Err"	An electronic failure has occurred within the NIBP module.	Contact CAS Medical Systems to have the monitor serviced.
"No Probe"	The monitor is not detecting the SpO <sub>2</sub> probe.	The probe was disconnected from either the Interface Cable or from the monitor.
"OverPres"	Cuff pressure exceeded 290 mmHg in the Adult mode or 145 mmHg in the Neonatal mode.	Very rapid squeezing of the cuff can cause this error.  Repeat the measurement.  If this message repeatedly occurs during normal use, the monitor must be serviced.
"P Search"	The monitor is searching for a Pulse signal.	Normal at power-up as the monitor searches for a pulse.  The probe position may have changed.  Check the probe site.
"Prb OFF"	The monitor is no longer receiving a patient signal from the SpO <sub>2</sub> probe.	The probe is no longer in contact with the patient.  Check the probe site.
"Pwr Fail"	Power was disconnected from the monitor.	Depress the SILENCE/RESET pushbutton to clear the message.  OR Re-cycle the monitor's power.

**Table 9: Error Messages on the Message Window**

**ERROR MESSAGES ON THE MESSAGE WINDOW (cont.)**

<b>ERROR MESSAGE</b>	<b>POSSIBLE CAUSE</b>	<b>POSSIBLE SOLUTION</b>
"RangeErr"	The systolic reading exceeds the measurement range of 255 mmHg in the Adult mode or 135 mmHg in the Neonatal mode.	Repeat measurement. If the message is displayed again, use another method to measure the patient's blood pressure.
"SetClock"	The monitor's clock needs to be set.	The monitor's time and date values are incorrect. Refer to Page 85, for information to set the Time and Date. The monitor's internal clock battery needs to be replaced. Contact CAS Medical Systems.
"Signal ?" (Masimo)	The quality of the signal level being received by the SpO <sub>2</sub> sensor is in question.	Ensure proper sensor type and application. Verify that the sensor emitter and detector are parallel to and directly opposing each other. Clean or replace the sensor.
"Sig Sat"	Motion pulses too strong.	Limit patient activity; the arm must be still and/or relaxed. Repeat measurement.
"Sys Err"	An electronic failure has occurred within the monitor.	Contact CAS Medical Systems to have the monitor serviced.
"Time Out"	The monitor was unable to complete a measurement within 120 seconds in the Adult mode or 90 seconds in the Neonatal mode.	An extremely long measurement can be due to a loose cuff, high blood pressure, or monitor re-pumps. Try measurement again. Try higher initial pressure. If message consistently reappears try using another means to obtain patient's blood pressure.
"Weak Sig"	The monitor did not detect any pulses during a NIBP measurement.	Check the fit of the cuff. Repeat measurement.

**Table 9: Error Messages on the Message Window**





# Section 9

## External Printer



## 9. EXTERNAL PRINTER

The following section is provided as an overview of the Citizen CMP-10 Mobile Printer as it is used with the CAS 740 Monitor.

**NOTE:**

For more detailed information on the Citizen Model CMP-10 Mobile Printer, refer to the User's Manual that was supplied with the printer.

**WARNING:**

The CAS 740 Monitor has been tested with the Citizen CMP-10 Mobile printer to comply with IEC 60601-1-1 and is the only printer that is recommended to be used with the monitor. If another printer is to be used, the user must read the Caution on Page 24 under LEAKAGE CURRENT TEST and follow the guidance given.

### PRINTER OVERVIEW

The Citizen CMP-10 Mobile Printer interfaces to the CAS 740 Monitor via an infrared IrDA port or by using the direct connect RS232 cable (supplied with printer). The IrDA ports are located on the top of the printer and on the bottom front panel of the CAS 740 Monitor.

When using the IrDA port, it is important to keep the two devices close together (less than three (3) feet/one (1) one meter) and in-line to maintain proper communications.

**CAUTION:**

For safe and proper usage of the external printer, please observe the following:

When using the printer:

- Avoid placing monitor in areas where fluid may enter the printer opening accidentally.
- Do not drop or bump the monitor.
- Avoid places subject to high or low temperature extremes.
- Avoid direct sunlight.
- Avoid dusty places and where corrosive gasses are generated.
- Never attempt to dismantle or repair the printer mechanism.

When handling the thermal paper:

- Store in a dark, cool and dry place.
- Do not place near organic solvents.
- Avoid contact with vinyl chloride films erasers or adhesive tapes for extended periods.
- Avoid exposure to high temperature, humidity, liquid, or sunlight.
- Always use specified thermal paper (CAS P/N 28-02-0077).



## PRINTER OPERATION

Position the printer's IrDA port window in-line with the IrDA window of the CAS 740 Monitor OR connect the RS232 direct cable from the printer to the 9-pin RS 232 connector, located on the rear panel of the CAS 740 Monitor.

**NOTE:**

When using the IrDA port for printing, it is recommended to keep all items clear of the communications path between the monitor and the printer.

Turn the printer "ON". To turn the printer "ON", press and hold the ON/OFF pushbutton for one (1) second. The Power LED illuminates Green.

To turn the printer "OFF", press and hold the ON/OFF pushbutton. The Power LED will illuminate Red and change back to Green. When the LED illuminates Green, remove your finger from the pushbutton.

The CMP-10 Mobile Printer contains an Auto Power Off feature. If the monitor and printer become separated and after ten (10) minutes during which no data has been sent to the printer and the Line Feed (LF) pushbutton has not been depressed, the printer will automatically shutoff.

The CAS 740 Monitor will periodically send to the printer a "wake-up" message that will disable the Auto Power Off feature.

**NOTE:**

If the printer was powered "ON" during the monitor's power-up sequence, the CAS Medical Systems logo is printed.



# CAS 740 Monitors

Sample printouts of both History Modes are shown below.

Header	<p>740X Series Monitor 04-Oct-02 13:56</p> <p>Patient: _____</p> <p>Notes: _____</p> <hr/> <p>Events:</p> <table border="0"> <thead> <tr> <th>Hr:Mn</th> <th>Sys</th> <th>Dia</th> <th>MAP</th> <th>Pls</th> <th>%O2</th> <th>Temp</th> </tr> </thead> <tbody> <tr> <td>13:57</td> <td>112</td> <td>78</td> <td>89</td> <td>72</td> <td>***</td> <td>****</td> </tr> <tr> <td>13:59</td> <td>115</td> <td>77</td> <td>88</td> <td>60</td> <td>98</td> <td>****</td> </tr> <tr> <td>14:01</td> <td>117</td> <td>78</td> <td>90</td> <td>70</td> <td>***</td> <td>****</td> </tr> <tr> <td>14:01</td> <td>120</td> <td>80</td> <td>100</td> <td>65</td> <td>***</td> <td>****</td> </tr> <tr> <td>14:03</td> <td>118</td> <td>78</td> <td>91</td> <td>77</td> <td>98</td> <td>****</td> </tr> <tr> <td>14:04</td> <td>***</td> <td>***</td> <td>***</td> <td>***</td> <td>***</td> <td>37.0°C</td> </tr> </tbody> </table>	Hr:Mn	Sys	Dia	MAP	Pls	%O2	Temp	13:57	112	78	89	72	***	****	13:59	115	77	88	60	98	****	14:01	117	78	90	70	***	****	14:01	120	80	100	65	***	****	14:03	118	78	91	77	98	****	14:04	***	***	***	***	***	37.0°C	<p>740X Series Monitor 04-Oct-02 13:56</p> <p>Patient: _____</p> <p>Notes: _____</p> <hr/> <p>Trends:</p> <table border="0"> <thead> <tr> <th>Hr:Mn</th> <th>Sys</th> <th>Dia</th> <th>MAP</th> <th>Pls</th> <th>%O2</th> <th>Temp</th> </tr> </thead> <tbody> <tr> <td>13:57</td> <td>112</td> <td>78</td> <td>89</td> <td>72</td> <td>***</td> <td>****</td> </tr> <tr> <td>13:59</td> <td>115</td> <td>77</td> <td>88</td> <td>60</td> <td>98</td> <td>****</td> </tr> <tr> <td>14:00</td> <td>***</td> <td>***</td> <td>***</td> <td>64</td> <td>99</td> <td>****</td> </tr> <tr> <td>14:01</td> <td>117</td> <td>78</td> <td>90</td> <td>70</td> <td>***</td> <td>****</td> </tr> <tr> <td>14:01</td> <td>120</td> <td>80</td> <td>100</td> <td>65</td> <td>***</td> <td>****</td> </tr> <tr> <td>14:03</td> <td>118</td> <td>78</td> <td>91</td> <td>77</td> <td>98</td> <td>****</td> </tr> <tr> <td>14:04</td> <td>***</td> <td>***</td> <td>***</td> <td>***</td> <td>***</td> <td>37.0°C</td> </tr> </tbody> </table>	Hr:Mn	Sys	Dia	MAP	Pls	%O2	Temp	13:57	112	78	89	72	***	****	13:59	115	77	88	60	98	****	14:00	***	***	***	64	99	****	14:01	117	78	90	70	***	****	14:01	120	80	100	65	***	****	14:03	118	78	91	77	98	****	14:04	***	***	***	***	***	37.0°C
	Hr:Mn	Sys	Dia	MAP	Pls	%O2	Temp																																																																																																				
13:57	112	78	89	72	***	****																																																																																																					
13:59	115	77	88	60	98	****																																																																																																					
14:01	117	78	90	70	***	****																																																																																																					
14:01	120	80	100	65	***	****																																																																																																					
14:03	118	78	91	77	98	****																																																																																																					
14:04	***	***	***	***	***	37.0°C																																																																																																					
Hr:Mn	Sys	Dia	MAP	Pls	%O2	Temp																																																																																																					
13:57	112	78	89	72	***	****																																																																																																					
13:59	115	77	88	60	98	****																																																																																																					
14:00	***	***	***	64	99	****																																																																																																					
14:01	117	78	90	70	***	****																																																																																																					
14:01	120	80	100	65	***	****																																																																																																					
14:03	118	78	91	77	98	****																																																																																																					
14:04	***	***	***	***	***	37.0°C																																																																																																					
	<p>740X Series Monitor 04-Oct-02 15:02</p> <p>Patient: _____</p> <p>Notes: _____</p> <hr/> <p>Events:</p> <table border="0"> <thead> <tr> <th>Hr:Mn</th> <th>Sys</th> <th>Dia</th> <th>MAP</th> <th>Pls</th> <th>%O2</th> <th>Temp</th> </tr> </thead> <tbody> <tr> <td>15:04</td> <td>112</td> <td>78</td> <td>89</td> <td>72</td> <td>***</td> <td>****</td> </tr> <tr> <td>15:09</td> <td>115</td> <td>77</td> <td>88</td> <td>60</td> <td>98</td> <td>****</td> </tr> </tbody> </table>	Hr:Mn	Sys	Dia	MAP	Pls	%O2	Temp	15:04	112	78	89	72	***	****	15:09	115	77	88	60	98	****	<p>740X Series Monitor 04-Oct-02 15:02</p> <p>Patient: _____</p> <p>Notes: _____</p> <hr/> <p>Trends:</p> <table border="0"> <thead> <tr> <th>Hr:Mn</th> <th>Sys</th> <th>Dia</th> <th>MAP</th> <th>Pls</th> <th>%O2</th> <th>Temp</th> </tr> </thead> <tbody> <tr> <td>15:04</td> <td>112</td> <td>78</td> <td>89</td> <td>72</td> <td>***</td> <td>****</td> </tr> <tr> <td>15:09</td> <td>115</td> <td>77</td> <td>88</td> <td>60</td> <td>98</td> <td>****</td> </tr> </tbody> </table>	Hr:Mn	Sys	Dia	MAP	Pls	%O2	Temp	15:04	112	78	89	72	***	****	15:09	115	77	88	60	98	****																																																															
Hr:Mn	Sys	Dia	MAP	Pls	%O2	Temp																																																																																																					
15:04	112	78	89	72	***	****																																																																																																					
15:09	115	77	88	60	98	****																																																																																																					
Hr:Mn	Sys	Dia	MAP	Pls	%O2	Temp																																																																																																					
15:04	112	78	89	72	***	****																																																																																																					
15:09	115	77	88	60	98	****																																																																																																					

Figure 20: Sample Printouts

**NOTE:**

When no patient information is stored in the Trend History, for consecutive minute storage, a break appears in the Trends printout. If the %O2 reading is "\*\*\*\*", then the Pulse value (Pls) shown is associated with the NIBP reading.

## CHARGING THE PRINTER BATTERY

The CMP-10 Mobile Printer is equipped with a rechargeable Lithium Ion (LiION) battery pack.

When the printer detects a Low Battery condition within itself, the message "Low Battery" is printed and an audio indicator, located inside the printer sounds three (3) times.

### **WARNING:**

Charge the printer battery using the AC Adapter, Model TRC-09-1100-M from GROUP WEST, or equivalent, included with the printer.


Plug the battery charger's cord into the printer battery charger jack, located on the rear panel. Plug the charger into an AC wall outlet of the appropriate voltage. Verify the CHARGE LED indicator is lit Red. Battery charge time is approximately three (3) hours. Once the battery is fully charged, the CHARGE LED indicator switches to Green.

## INSTALLING PAPER

### **NOTE:**

A red line appears when the remaining supply of thermal paper becomes low.

- 1) Switch the printer "OFF".
- 2) Press the Cover Open button to access the paper compartment. Remove any remaining paper before installing the new roll.
- 3) Place the new paper roll as shown on the illustration and pull out enough paper to reach out over the control panel of the printer.
- 4) Close the paper door.

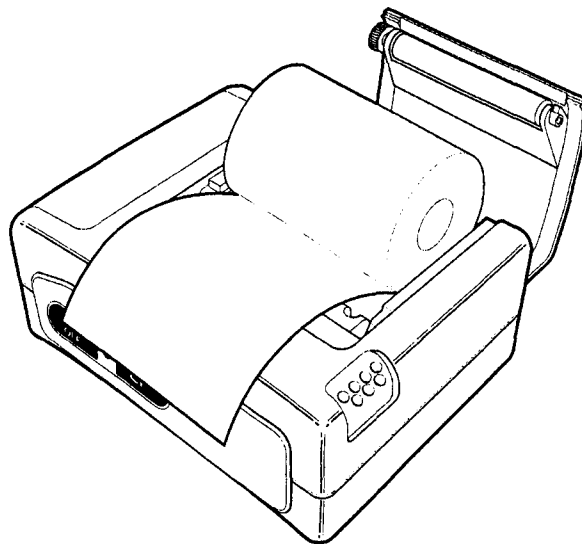

**Figure 21: Paper Installation**

**NOTE:**

Make sure that the paper is correctly placed. If it is tilted in one or another direction and does not come out straight from under the cover, open the door and reposition the roll again.

**WARNING:**

Do not touch the print head or paper cutter while replacing the printer paper.

### REPLACING THE BATTERY PACK

**WARNING:**

Do not operate the printer or connect the printer to the CAS 740 Monitor with the battery pack removed.

**WARNING:**

Never change the battery pack while the battery charger is plugged in and/or the CAS 740 Monitor is being operated.

- 1) Switch the printer "OFF".
- 2) Disconnect the printer from the CAS 740 Monitor and unplug the wall charger cord.
- 3) Open the battery door by pressing in on the battery cover and pushing upward.

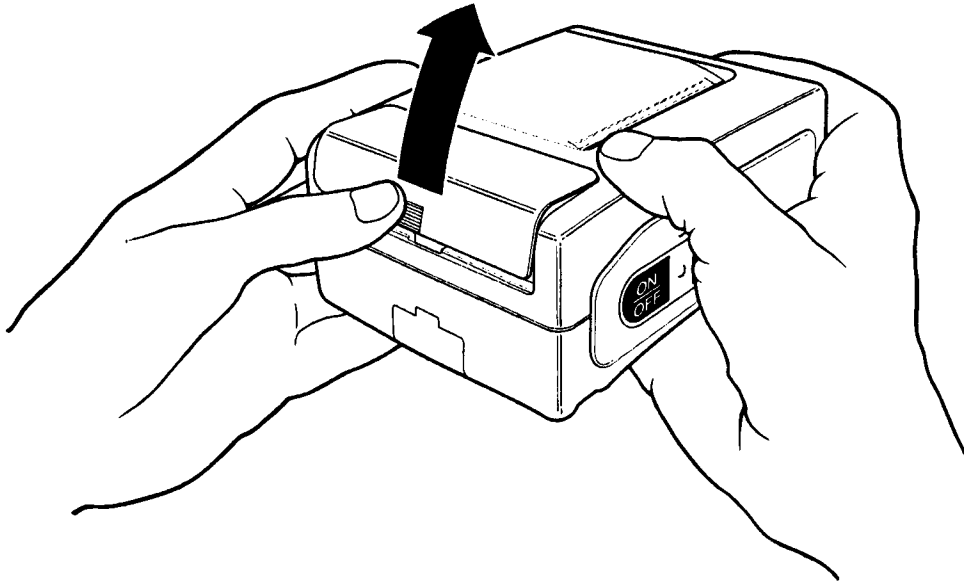


Figure 22: Opening the Battery Door

- 4) Remove the battery cover.
- 5) Remove the battery pack from the compartment and disconnect its connecting cable.



## INSTALLING A NEW BATTERY PACK

- 1) Connect the battery cable into the battery connector.
- 2) Insert the battery and its connecting cable into the battery compartment.

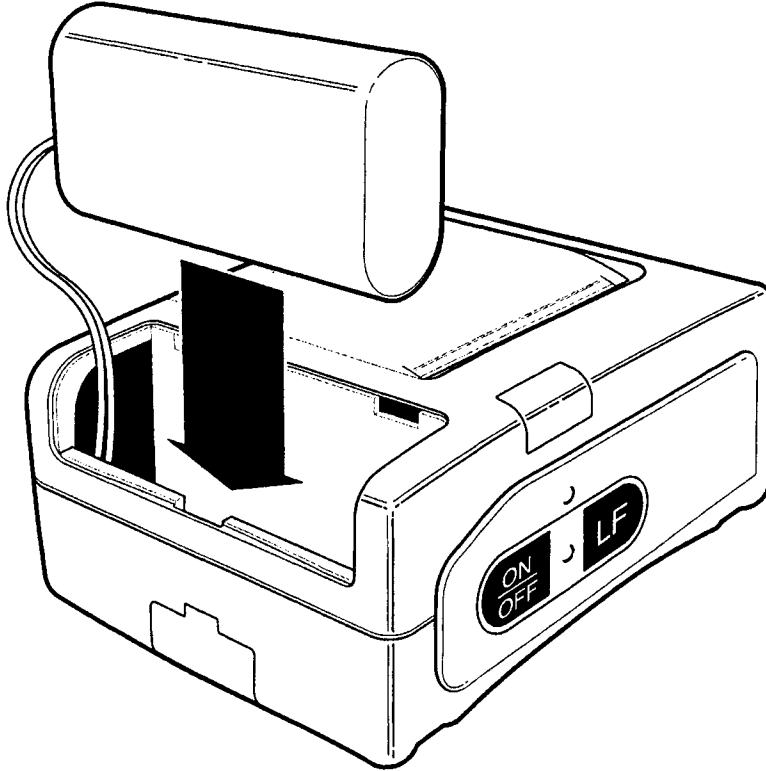


Figure 23: Installing the New Battery

- 3) Replace the battery cover by sliding it in from the back of the printer and pushing down to lock it in place.

**CAUTION:**

Be sure to place the battery cover firmly in its position after installing the new battery pack.

**WARNING:**

Do not disassemble the battery pack or batteries. The batteries contain electrolytes, which can cause injury to eyes, skin and clothing.

**NOTE:**

This product contains a rechargeable battery that is recyclable. Under various state and local laws, it may be illegal to dispose of this battery into the municipal waste stream. Check with your local authorities for instructions on recycling options in your area.



# Section 10

## Cleaning



## 10. CLEANING

### CLEANING OVERVIEW

**WARNING:**

Do not, under any circumstances, perform any testing or maintenance on the monitor while the monitor is being used to monitor a patient. The monitor must be turned "OFF". Unplug the monitor from the AC power source and remove the internal battery.

**CAUTION:**

Do not open the monitor to clean or repair it. Contact CAS Medical System for service needs.

**CAUTION:**

Disconnect all accessories from the monitor before cleaning. Do not immerse any part of the electrical connector of the cable or accessories in the cleaning or disinfection solution at any time. Do not use an abrasive cloth or cleaner on the accessories.

### THE MONITOR

On a daily basis, examine the monitor's case for any damages and check the AC power cord for bent or broken prongs, cracks or fraying. Neither the monitor nor the power cord should be used if damaged. If any damage is noted, contact the appropriate service personnel.

**CAUTION:**

Do not spray any water or cleaning solution directly onto the monitor.

As needed, clean the monitor using a soft cloth dampened with a mild dishwashing detergent solution and gently rub the soiled area until clean. Use a clean soft cloth to dry the monitor. Do not use abrasive cleaners on the monitor. Do not use either isopropyl alcohol or solvent to clean the monitor. Use of these cleaners can cause damage to the monitors' surface. Do not immerse the monitor or power cord in the cleaning solution.

When necessary, the monitor surfaces may be disinfected using a soft cloth saturated with a 10% (1:10) solution of chlorine bleach in tap water. When all of the surfaces have been disinfected, wipe the entire surface of the monitor using a soft cloth dampened with fresh water to remove any trace amounts of residue and/or fumes.

**NOTE:**

Thoroughly wipe off any excess cleaning solutions. Care should be taken to prevent water or cleaning solution to run into connector openings or crevices.

## THE DISPLAY

**CAUTION:**

Use care when cleaning the display. Scratches may occur.

Occasionally, as needed, clean the display window using a soft, lint-free cloth sprayed with an alcohol free glass cleaner. Do not use either isopropyl alcohol or solvent to clean the display. Use of these cleaners can cause damage to the display. The use of paper towels is not recommended as it may scratch the surface.

## CUFFS

Prior to each patient use, inspect the blood pressure cuff and its hose for damage.

### TUFF-CUFF® REUSABLE CUFFS

As necessary, clean the blood pressure cuff using a soft cloth dampened with a 70% Isopropyl Alcohol solution.

**NOTE:**

CAS does not recommend submersion of the cuff. Liquid should not be permitted to enter the cuff bladder because instrument damage may occur. The cuff should be allowed to thoroughly dry before use.

### PEDISPHYG® CUFFS

CAS is aware that, in certain situations, the cuff may become soiled during its use. In these situations a water-based detergent is suitable for wiping the cuff.

As necessary, the preferred method for cleaning the Pedisphyg Cuff is to wipe it down with a damp, soapy cloth. A damp, detergent-free cloth should then be used to rinse the cuff.

**NOTE:**

CAS does not recommend submersion of the cuff. Liquid should not be permitted to enter the cuff bladder because instrument damage may occur. The cuff should be allowed to thoroughly dry before use.

## PNEUMATIC TUBING

Prior to each patient use, inspect the NIBP Inflation Hose for proper connection, cracks and kinks. As necessary, clean the pneumatic tubing using a soft cloth dampened with a germicidal solution.

**PRINTER**

When the printer becomes dirty, wipe with a soft dry cloth. For extreme dirt buildup, soak a cloth with mild detergent, wring well and wipe. Dry by wiping with a soft dry cloth.

**CAUTION:**

Before cleaning the printer, disconnect the AC adapter from the printer.  
Do not use volatile chemicals such as thinner, benzene, etc.  
Never wet the inside of the printer mechanism.

Refer to the printer User's Manual for more information.

**SpO<sub>2</sub> SENSORS**

(Reusable)

As necessary, the sensor may be surface cleaned by wiping it with a 70% isopropyl alcohol pad. Allow the sensor to dry prior to placement on a patient.

**CAUTION:**

Do not soak or immerse the sensor or its cable in any liquid solution. Do not attempt to sterilize.

Refer to the Directions For Use pamphlet enclosed with each sensor for more information.

**TEMPERATURE PROBES**

The temperature probe should periodically be cleaned by wiping with an alcohol dampened cloth or wipe, warm water, or properly diluted non-staining disinfectant. Do not immerse the probes.

**CAUTION:**

Do not soak or immerse the probe or its cable in any liquid solution. Do not attempt to sterilize.





# Section 11

## Maintenance



## 11. MAINTENANCE

### MAINTENANCE INTERVALS

Preventive maintenance of the monitor is an important function that should be performed routinely by the user to ensure safe and efficient monitor operation. CAS Medical Systems recommends that you do the following:

- Perform a pneumatic check every six (6) months.
- Perform a calibration check once per year or when there is doubt about the validity of the pressure readings.
- Replace the battery pack every two (2) years.




If the monitor is in need of repair, it must be referred to the appropriate service personnel. Service performed by unauthorized personnel could be detrimental to the monitor and may void the warranty. For service, contact your dealer or CAS Medical Systems, Inc.

### TEST MODE

The monitor must be in the Configuration Menu in order to perform the following functions:

- All LEDs "ON" Check
- Calibration Check
  - System Pressure
  - Over Pressure
- Pneumatic Pressure Checks
- Temperature Calibration Check
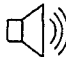
### ENTER THE TEST MODE

To enter the monitor's Configuration Menu, depress and hold the AUDIO  and ALARM LIMITS  pushbutton keys while the monitor is being turned "ON" .

Once in the menu, depress one of the PREVIOUS (CYCLE TIME / ALARM LIMITS) programmed pushbutton keys until the Message Window briefly displays "TestMode" followed by "0 mmHg".

#### NOTE:

While in the Test Mode if no pushbutton is depressed within 15 minutes, the monitor will automatically terminate the Monitor Configuration menu and return to the "Ready" mode.

#### WARNING:

Do not place the monitor in the TEST MODE when a cuff is attached to a patient.

### EXIT THE TEST MODE

When you have completed with the Test Mode, press the CANCEL pushbutton to exit. The Message Window will briefly display "Saving" and return to the "Ready" state.

### LED CHECK

The CAS 740 Monitor incorporates an all lights "ON" check to verify the functionality of the front panel LED displays and indicators.

Enter the Test Mode. Refer to Page 115, ENTER THE TEST MODE. The Message Window will briefly display "TestMode" followed by "0 mmHg".

Depress and hold either the ARROW UP or ARROW DOWN pushbuttons. The monitor will illuminate all appropriate 7-segment displays, bar graph, bell icon and patient mode indicators for as long as the button is depressed.

### CALIBRATION CHECK

Verify the calibration of the monitor once (1) per year.

A Calibration Kit, (product #P9) is included with the monitor. The kit contains a T-connector with a male and a female luer fitting (for a Calibration Check) and a male luer plug (to be used for the Pneumatic Check).

Obtain a mercury manometer whose accuracy meets the AAMI/ANSI Standard for Non-Automated Sphygmomanometers, 2002.

#### **NOTE:**

Monitor must be in the Adult mode prior to performing these pressure checks.

### SYSTEM PRESSURE

Assemble the Calibration Kit according to the diagram provided in the P9 kit.

- 1) Remove the manometer tubing from the inflation bulb. Connect the open ended tubing of the T-connector to the inflation bulb.
- 2) Connect the female luer fitting to the inflation tube leading to the manometer.
- 3) Connect the male luer fitting to the manometer tubing.
- 4) Enter the Test Mode. Refer to Page 115, ENTER THE TEST MODE. The Message Window will briefly display "TestMode" followed by "0 mmHg".
- 5) Use the manometer inflation bulb to slowly inflate the system pausing for 30 seconds at the following points and verify calibration according to the following table:

0 mmHg +/- 1 mmHg

50 mmHg +/- 4 mmHg

100 mmHg +/- 4 mmHg

150 mmHg +/- 4 mmHg

200 mmHg +/- 5 mmHg

**NOTE:**

If the monitor does not display the test pressure for the 30-second period, deflate to zero and verify the proper assembly of the calibration set-up. Re-inflate the system. If the monitor again fails to hold the pressure, it is recommended the monitor be returned to CAS Medical Systems for service.

**OVERPRESSURE**

Inflate the pressure slowly until 290 mmHg +/- 10 mmHg is reached. The Message Window should stop updating, display the message "OverPres" and provide a NIBP Application Error audible tone.

Press the CANCEL pushbutton to exit the Overpressure Test. The monitor returns to the Calibration Check function.

If the monitor does not meet the above specifications, it is recommended the monitor be returned to CAS Medical Systems for service.

**PNEUMATIC PRESSURE CHECKS**

Check the monitor's pneumatic system for air leakage every six (6) months.

**PLUG TUBE**

Obtain the male luer plug found in the Calibration Kit (product #P9) supplied with the monitor.

- 1) Place this plug into the cuff connector at the end of the monitor inflation hose and twist one-quarter turn. The plug must fit securely into the connector for this test to be performed properly.
- 2) Enter the Test Mode. Refer to Page 115, ENTER THE TEST MODE. The Message Window will briefly display "TestMode" followed by "0 mmHg".
- 3) Press the NIBP START pushbutton to begin the Pressure Check.
- 4) The Message Window will display "Chk Prs", will inflate to approximately 180 mmHg and attempt to hold this pressure. The pressure value will be displayed in the SYSTOLIC display window. This test takes about fifteen (15) seconds.

## CAS 740 Monitors

- 5) At the completion of a successful Pressure Check, the Message Window will display "Passed", the monitor will beep two (2) times and will return to the Calibration Check function after five (5) seconds.
- 6) If the monitor fails the Pressure Check, the Message Window will display "Leak", the monitor will beep three (3) times and the return to the Calibration Check function after five (5) seconds.

Due to the volume differences of the hoses offered with the CAS 740 Monitor, the monitor may incorrectly fail the Plug Tube check. Should the monitor fail the Plug Tube Pressure Check, obtain a 500 ml Pressure Cylinder and follow the 500 ml Pressure Check.

### 500 ml PRESSURE CHECK

Obtain a fixed volume 500 ml Pressure Cylinder (CAS p/n 01-02-0248).

- 1) Place the end of the monitor's inflation hose securely onto the luer fitting at the top of the pressure cylinder. The hose must fit securely onto the connector for this test to be performed properly.
- 2) Enter the Test Mode. Refer to Page 115, ENTER THE TEST MODE. The Message Window will briefly display "TestMode" followed by "0 mmHg".
- 3) Press the NIBP START pushbutton to begin the Pressure Check.
- 4) The Message Window will display "Chk Prs", will inflate to approximately 160 mmHg and attempt to hold this pressure. The pressure value will be displayed in the SYSTOLIC display window. This test takes about fifteen (15) seconds.
- 5) At the completion of a successful Pressure Check, the Message Window will display "Passed", the monitor will beep two (2) times and will return to the Calibration Check function after five (5) seconds.
- 6) If the monitor fails the Pressure Check, the Message Window will display "Leak", the monitor will beep three (3) times and the return to the Calibration Check function after five (5) seconds.

Should the monitor fail the 500 ml Pressure Check, it is recommended the monitor be returned to CAS Medical Systems for service.

**TEMPERATURE CALIBRATION CHECK**

Verify the calibration of the monitor's Temperature circuit, once (1) every year.

To perform a Temperature Calibration Check, obtain a Calibration Key. This key can be purchased directly from Welch Allyn or CAS Medical Systems. Refer to Section 13, ACCESSORIES for part number information.

**NOTE:**

The Temperature Calibration Check can be performed, at any time, once the monitor enters the Configuration Menu.

- 1) Turn the monitor "OFF".
- 2) Remove the temperature probe and its connector completely from the monitor and insert the Calibration Key.
- 3) Turn the monitor "ON", and enter the monitor's Configuration Menu. Refer to Page 115, ENTER THE TEST MODE.
- 4) Re-insert and remove the temperature probe from the probe guide to reset the thermometer's electronics. An audible "click" should be heard when the probe is completely placed into its holder.
- 5) Wait for the test to complete, and observe the display reading in the TEMP display window.
- 6) The display value should read 36.3 +/- 0.1 °C or 97.3 +/-0.1°F.

**NOTE:**

The monitor will display the Temperature Calibration Key value using the current temperature units selected.

**NOTE:**

The Temperature Calibration Key will only operate while the monitor is in the Configuration Menu. If the Calibration Key is inserted during normal monitoring, the TEMP display will show three (3) flashing dashes "- - -" and no value will be displayed.

**NOTE:**

While in the Configuration Menu, if no pushbutton is depressed within 60 seconds, the monitor will automatically exit the Monitor Configuration menu. The Message Window will briefly display "Saving" and return to the "Ready" mode.

- 7) Remove the Calibration Key and re-insert the temperature probe connector.
- 8) Install the temperature probe into the probe holder. An audible "click" should be heard when the probe is completely placed into its holder.

## OXIMETRY CALIBRATION CHECK

The oximeter is factory calibrated to determine the percentage of arterial oxygen saturation of functional hemoglobin. No user calibration is required.

## REPLACING THE MONITOR BATTERY

A part number for the battery can be found in the Accessories section of this manual or on the label located on the inside panel of the battery pack. When the battery fails to hold a charge it will need to be replaced.

CAS Medical Systems recommends the battery be changed every two (2) years.

## REMOVING THE BATTERY

- 1) Turn the monitor "OFF" and disconnect the power cord.
- 2) Push down on the battery latch to unlock the battery door from the rear panel of the monitor.
- 3) Carefully remove the battery pack from the rear panel of the monitor. Refer to Figure 24.

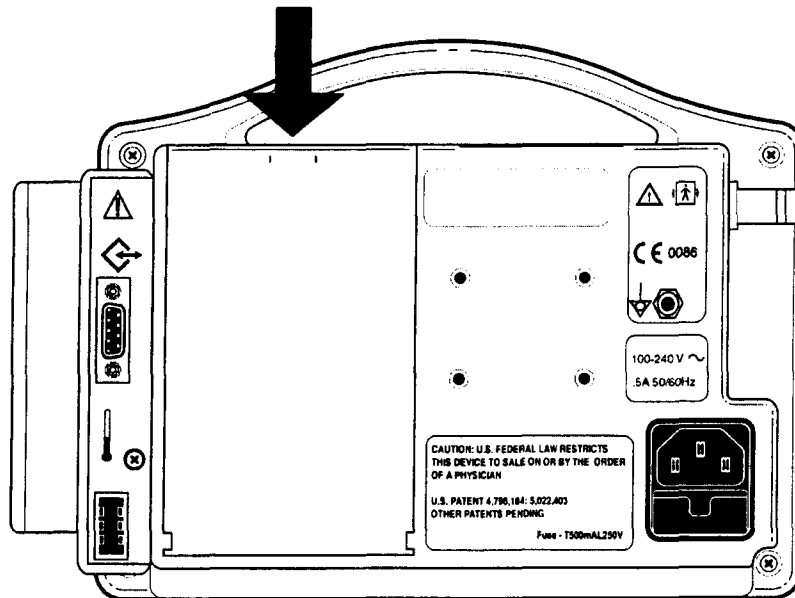


Figure 24: Removing the Monitor Battery Pack



**INSTALLING THE BATTERY**

- 1) Align the Battery Pack guides with the bottom of the monitor.
- 2) Slowly close the battery door to ensure the connector in the monitor and the connector on the battery pack mate together.
- 3) Lock the battery door closed.

Refer to Page 86, BATTERY POWER for additional battery information.

**NOTE:**

When the battery pack is re-inserted, the monitor will automatically turn "ON".

**WARNING:**

Do not disassemble the battery pack or batteries. The batteries contain electrolytes, which can cause injury to eyes, skin and clothing.

**NOTE:**

This product contains a rechargeable battery that is recyclable. Under various state and local laws, it may be illegal to dispose of this battery into the municipal waste stream. Check with your local authorities for instructions on recycling options in your area.

**CHANGING THE FUSES**

The CAS 740 Monitor uses a dual fuse power input receptacle. The receptacle incorporates fuses in the hot and neutral AC input lines that are user serviceable.

The two (2) fuses for the CAS 740 Monitor are each rated at 250V, 500mA, 5 x 20 mm, Slow Blow. Refer to Section 13, ACCESSORIES for part number information.

**CAUTION:**

For continued protection against fire hazard, replace only with identically rated fuses.

A fuse may need to be replaced if the monitor is plugged into an electrical outlet but the Battery Power Visual Indicator is not illuminated the color Green.

**WARNING:**

Before changing the fuse, unplug the power cord.

The fuse holder is incorporated into the power input receptacle and located under the power cord input connector.

## CAS 740 Monitors

To replace the fuses:

- 1) Turn the monitor "OFF" and disconnect the power cord.
- 2) Depress down on the locking tab, which holds the fuse holder in the power input receptacle.
- 3) While holding down on the tab, pull the fuse holder out.
- 4) Remove the fuses.
- 5) Place new fuses directly into the fuse holder.
- 6) Insert the fuse holder into the power input receptacle. There should be an audible "click" when it is secure.

**NOTE:**

The CAS 740M (EMS Monitor) uses a single fuse located on the inside of the monitor. This fuse is not user replaceable. If you suspect a defective fuse, (the Battery Power Visual Indicator is not illuminated), it is recommended to contact CAS Medical for service.

### STORAGE

**WARNING:**

If it becomes necessary to store the monitor for longer than six (6) months, remove the monitor's battery pack and place the monitor in its original packing container if available.

**WARNING:**

Use of unapproved batteries will invalidate the product's warranty and may result in serious safety consequences for the patient and user.

See Section 15, SPECIFICATIONS for storage temperature information.

# Section 12

## External Device Interfacing



## 12. EXTERNAL DEVICE INTERFACING

### OVERVIEW

The CAS 740 Monitor is capable of interfacing to an external Serial printer or have the ability to interface to a Nurse Call System if the DB9 RS232 option is available. Both connections are made through the DB9 connector located on the rear panel of the monitor.

**WARNING:**

The CAS 740 Monitor has been tested with the Citizen CMP-10 Mobile printer to comply with IEC 60601-1-1 and is the only printer that is recommended to be used with the monitor. If another printer is to be used, the user must read the Caution on Page 24 under LEAKAGE CURRENT TEST and follow the guidance given.

### RS232

The CAS 740 monitor uses the DB9 connector to interface to the Citizen CMP-10 Mobile printer using the cable supplied with the printer. The connector information provided in this section is made available to allow the user the ability to print the monitor's History data to an external serial printer. Refer to Figure 25 and Table 10 for connection information.

Refer to Section 15, SPECIFICATIONS for Serial Interface information.

### NURSE CALL INTERFACE

The CAS 740 Monitor provides an isolated relay switch closure output connection between two (2) of the pins on the DB9 RS232 output connector. The output is compatible with most Nurse Call Systems in that there is no polarity to the connection.

When properly connected, the Nurse Call Interface activates the Nurse Call System each time an alarm is activated on the monitor. The delay time for the Nurse Call Interface to activate is less than 0.5 seconds.

The Nurse Call System's relay contacts are rated at 120 VAC at 0.3A; or 30 VDC at 1.0 A.

The Nurse Call Option is available as a normally open (closed on alarm) or normally closed (open on alarm) depending upon how it is wired.

For normally open (N.O.) applications, the Nurse Call system needs to be connected to pins 1 and 9 of the RS232 connector. For normally closed (N.C.) applications, the Nurse Call system needs to be connected to pins 6 and 9. Refer to Figure 25 and Table 10 for connection information.

**WARNING:**

The connection to the Nurse Call Interface should only be installed by a qualified service personnel.

**WARNING:**

The interconnection of auxiliary equipment to the Nurse Call Interface may increase the total leakage current. The user must read the Caution on Page 24 under LEAKAGE CURRENT TEST and follow the guidance given.

**NOTE:**

Even though the Nurse Call Interface allows remote alarm indication, it does not replace appropriate bedside surveillance by trained clinicians.

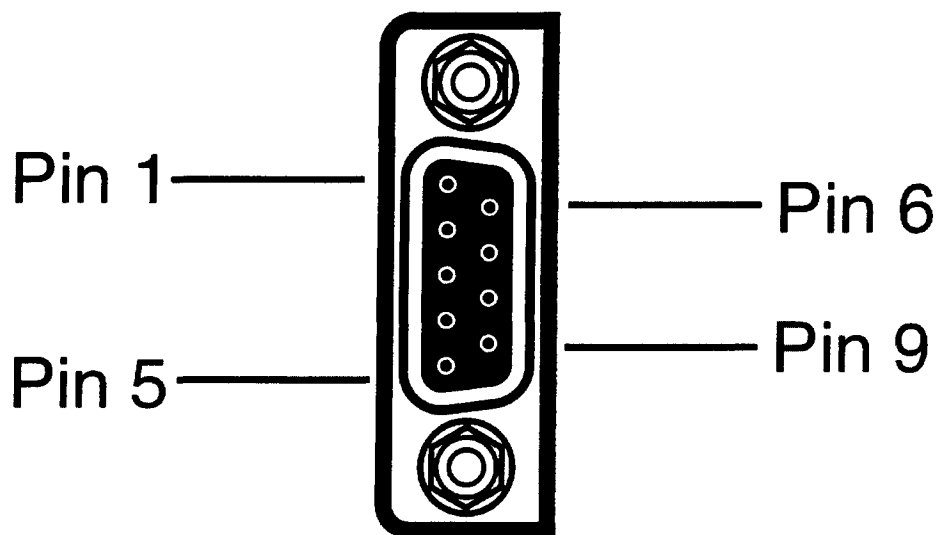


Figure 25: DB9 Male Connector Pin Layout

Pin Number	Signal Description
1	Nurse Call (N.O.)
2	Serial Receive In
3	Serial Transmit Out
4	No Connection
5	Isolated Ground
6	Nurse Call (N.C.)
7	No Connection
8	No Connection
9	Nurse Call (common)

Table 10: DB9 Pin Out

# Section 13

## Accessories





## 13. ACCESSORIES

### BLOOD PRESSURE CUFFS

#### Tuff-Cuff®

Reusable Blood Pressure Cuffs (single tube)

Catalog Number	Description	Size
CR 5216	Large Adult	16 cm x 42 cm
CR5214	Adult	14 cm x 37 cm
CR5212	Small Adult	12 cm x 30 cm
CR5209	Child	9 cm x 27 cm
CR5207	Small Child	7 cm x 21 cm
CR5206	Infant	6 cm x 18 cm

#### Safe-Cuff™

Single – Patient Use Blood Pressure Cuffs (single tube)

Catalog Number	Description	Size
CD2060	X-Large Adult	20 cm x 52 cm
CD1642	Large Adult	16 cm x 41 cm
CD1437	Adult	14 cm x 36 cm
CD1230	Small Adult	12 cm x 31 cm
CD927	Child	9 cm x 25 cm
CD618	Infant	6 cm x 16 cm

#### Pedisphyg®

Single – Patient Use Neonatal Blood Pressure Cuffs (single tube)

Catalog Number	Description	Size
C26	2.5 cm	2.5 cm x 9.0 cm
C39	3.0 cm	3.0 cm x 11.5 cm
C412	4.0 cm	4.0 cm x 14.5 cm
C515	5.0 cm	5.0 cm x 17.75 cm

### INFLATION HOSES

Catalog No.	Description
01-01-0495	Ten (10) FT Coiled NIBP Hose, Adult and Pediatric, Bayonet Cuff
01-02-0185	Six (6) FT Straight NIBP Hose, Neonatal and Infant
01-02-0131	Ten (10) FT Coiled NIBP Hose, Adult and Pediatric

# CAS 740 Monitors

## OXIMETRY

### Masimo

**Note: Use only with Masimo Pulse Oximeter**

Catalog No.	Description	
01-02-0251	LNOP ADT	Adult Single Patient Adhesive Sensor, >30 kg
01-02-0252	LNOP ADT LNG	Adult Single Patient Sensor, Long Cbl, >30 kg
01-02-0253	LNOP PDT	Pediatric Slender Digit Single Patient Adhesive Sensor, >10<50 kg
01-02-0254	LNOP NEO	Neonatal Single Patient Adhesive Sensor, >1<10 kg
01-02-0255	LNOP NEOPT	Neonatal Preterm Single Patient Adhesive Sensor, <1 kg
01-02-0178	LNOP DCI	Adult Reusable Finger Clip Sensor, >30 kg
01-02-0182	PC04	Patient Cable, 4 FT
01-02-0190	LNOP DCIP	Pediatric Sensor
01-02-0191	LNOP YI	Multisite Sensor
01-02-0192	PC08	Patient Cable, 8 FT
01-02-0312	DC-195	Adult Reusable Finger Clip Sensor, >30 kg
01-02-0432	LNOP TC-1	Reusable Ear Sensor, >30 kg

### Nellcor

**Note: Use only with Nellcor Pulse Oximeter**

Catalog No.	Description	
01-02-0179	DS-100A	Adult Reusable Finger Clip Sensor
01-02-0183	DOC-10	OxiMax Patient Cable, 10 FT

### Nonin

**Note: Use only with Nonin Pulse Oximeter**

Catalog No.	Description	
01-02-0100	8500I	SpO <sub>2</sub> Extension Cable, 1 Meter
01-02-0106	8000AP-3	Pediatric Finger Clip Sensor, 9 ft / 3 Meter Cable
01-02-0108	8000Q	Ear Clip Sensor
01-02-0244	8000JFW	Sensor Attachment Tape, 25/case
01-02-0136	8500VI	SpO <sub>2</sub> Extension Cable, 3 Meter
01-02-0250	8000AA-2	Adult Articulated Finger Clip Sensor, 6 ft / 2 Meter Cable
01-03-0117	8000AA-3	Adult Articulated Finger Clip Sensor, 9 ft / 3 Meter Cable
01-03-0120	8000AA-1	Adult Articulated Finger Clip Sensor, 3 ft / 1 Meter Cable
01-03-0121	8000AP	Pediatric Finger Clip Sensor, 3 ft / 1 Meter Cable
01-03-0122	8000K2	Adult Finger Clip Sensor, 3 ft / 1 Meter Cable
01-03-0123	8000R	Reflectance Sensor
01-03-0124	8000J	Adult Flex Sensor (Straight)
01-03-0125	8000H	Reflectance Sensor Holder System
01-03-0126	8000T	Sensor Attachment Tape
01-03-0127	8000TH	Hydro-gel Tape Strips

**TEMPERATURE****Welch Allyn**

Catalog No.	Description	
01-02-0095	06137-000	SureTemp® Calibration Key
01-02-0096	02678-100	Temperature Probe, Oral
01-02-0097	02679-100	Temperature Probe, Rectal
01-02-0103	05031	Temperature Probe Covers, case of 40 boxes (25 Probe covers/box)

**OTHER ACCESSORIES**

Catalog No.	Description
01-01-0047	P9 Calibration Kit (includes T - connector with tubing and male luer plug)
01-02-0172G	Roll Stand with Basket (US Government customers only)
01-02-0173	Swivel Mount Kit
01-02-0174	Carry Bag
01-02-0176	RS232 Interface
01-02-0181G	Printer Bracket Attachment for Roll Stand
01-02-0188	Printer Battery
01-02-0189	Printer, includes Battery, RS232 Cable, Power Supply, One (1) Roll of Paper and Manual
01-02-0243	Universal Mount
01-02-0248	500 ml Fixed Volume Cylinder
01-02-0297	Roll Stand with Quick Release, Steelcraft
01-02-0301	Printer Bracket, used with Steelcraft Roll Stand
01-04-0016	U.S.A. Replacement Power Cord
01-05-0144	International Power Cord IEC320 European Plug
03-08-0386	DC Power Cable
03-08-0450	Monitor Battery Pack (7.2 VDC, 3700 mAh)
09-01-0002	Fuse (250V, 500mA, 5 x 20mm, Slow Blow), 2 per monitor
18-04-0004	Australian Power Cord
18-04-0010	U.K. Power Cord
21-02-0171	CAS 740 User's Manual, English
21-02-0173	CAS 740 User's Manual, Spanish
21-02-0175	CAS 740 User's Manual, German
21-02-0176	CAS 740 User's Manual, French
21-02-0177	CAS 740 User's Manual, Italian
21-02-0178	CAS 740 User's Manual, Dutch
21-02-0179	CAS 740 User's Manual, Portuguese
21-02-0180	CAS 740 User's Manual, Swedish
21-02-0174	CAS 740 Service Manual
28-02-0077	Printer Paper, One (1) Roll

**MONITOR CONFIGURATIONS**

<b>Model</b>	<b>Description</b>
CAS 740-1	MAXNIBP®, 100-240V, 50/60HZ, AC Power Supply and Battery
CAS 740M-1	MAXNIBP®, 12VDC Power input with Battery, Mount included
CAS 740-2MS	MAXNIBP® and Masimo SpO <sub>2</sub> , 100-240V, 50/60HZ, AC Power Supply and Battery
CAS 740M-2MS	MAXNIBP® and Masimo SpO <sub>2</sub> , 12VDC Power input with Battery, Mount included
CAS 740-2NL	MAXNIBP® and Nellcor SpO <sub>2</sub> , 100-240V, 50/60HZ, AC Power Supply and Battery
CAS 740M-2NL	MAXNIBP® and Nellcor SpO <sub>2</sub> , 12VDC Power input with Battery, Mount included
CAS 740-2NN	MAXNIBP® and Nonin SpO <sub>2</sub> , 100-240V, 50/60HZ, AC Power Supply and Battery
CAS 740M-2NN	MAXNIBP® and Nonin SpO <sub>2</sub> , 12VDC Power input with Battery, Mount included
CAS 740-2T	MAXNIBP® and Temperature, 100-240V, 50/60HZ, AC Power Supply and Battery
CAS 740M-2T	MAXNIBP® and Temperature, 12VDC Power input with Battery, Mount included
CAS 740-3MS	MAXNIBP®, Masimo SpO <sub>2</sub> , and Temperature, 100-240V, 50/60HZ, AC Power Supply and Battery
CAS 740M-3MS	MAXNIBP®, Masimo SpO <sub>2</sub> , and Temperature, 12VDC Power input with Battery, Mount included
CAS 740-3NL	MAXNIBP®, Nellcor SpO <sub>2</sub> , and Temperature, 100-240V, 50/60HZ, AC Power Supply and Battery
CAS 740M-3NL	MAXNIBP®, Nellcor SpO <sub>2</sub> , and Temperature, 12VDC Power input with Battery, Mount included
CAS 740-3NN	MAXNIBP®, Nonin SpO <sub>2</sub> , and Temperature, 100-240V, 50/60HZ, AC Power Supply and Battery
CAS 740M-3NN	MAXNIBP®, Nonin SpO <sub>2</sub> , and Temperature, 12VDC Power input with Battery, Mount included

**Table 11: Monitor Configurations**

# Section 14

## Glossary



## 14. GLOSSARY

### Diastolic

The Diastolic pressure is the bottom (lower) number given in a blood pressure reading. This number represents the amount of pressure present in the system between heartbeats.

### Mean Arterial Pressure (MAP)

The MAP pressure represents the mean pressure in the artery during one cardiac cycle.

### NIBP

Non-Invasive Blood Pressure. Abbreviation used when a blood pressure reading is taken externally using a non-invasive procedure.

### SpO<sub>2</sub> (%SpO<sub>2</sub>)

Abbreviation used when oxygen saturation level of the blood is measured with a pulse oximeter.

### Systolic

The Systolic pressure is the upper (top) number given in a blood pressure reading. This number represents the maximum pressure present in the artery.





# Section 15

## Specifications



## 15. SPECIFICATIONS

### NIBP MEASUREMENT

Characteristic	Specification	
Technique:	Oscillometric (MAX NIBP® Technology) Microprocessor software eliminates most ambient noise and motion artifact.	
Patient Range:	Neonate – Adult	
Blood Pressure Range	<u>NEO</u>	<u>ADULT</u>
Systolic:	30 – 135 mmHg	30 – 255 mmHg
Diastolic:	15 – 110 mmHg	15 – 220 mmHg
MAP:	20 – 125 mmHg	20 – 235 mmHg
Pulse Rate Range:	40 – 240 BPM	30 – 240 BPM
Accuracy		
Blood Pressure:	+/-5 mmHg with a standard deviation no greater than 8 mmHg (See Standards)	
Pulse Rate:	+/-2% or +/-2 BPM, whichever is greater	

### OXIMETRY (OPTIONS)

Characteristic	Specification	
Nonin®		
Type:	Functional Oxygen Saturation	
SpO <sub>2</sub> % Range:	0 - 100%	
SpO <sub>2</sub> Accuracy:	<u>Sensor</u> 8000AA 8000AP 8000K2	<u>Accuracy</u> 70 - 100%, +/-2 digits (1 S.D.)
	<u>Sensor</u> 8000J 8000R	<u>Accuracy</u> 70 - 100%, +/-3 digits (1 S.D.)
	<u>Sensor</u> 8000Q	<u>Accuracy</u> 70 - 100%, +/-4 digits (1 S.D.)
Measurement Wavelengths:	Red 660 nanometers Infrared 910 Nanometers	
Power:	3 mW nominal	
Pulse Rate Range:	18 - 240 BPM	
Pulse Rate Accuracy:	+/-3% or +/-1 digit, whichever is greater	
Numerics:	Updated every one (1) second.	

**NOTE:**

For further information on sensors and sensor accuracy, contact Nonin.

# CAS 740 Monitors

Characteristic	Specification																				
Masimo SET®																					
Type:	Functional Oxygen Saturation																				
SpO <sub>2</sub> % Range:	0 - 100%																				
SpO <sub>2</sub> Accuracy:	<table border="0"> <thead> <tr> <th>Sensor</th> <th>Accuracy</th> </tr> </thead> <tbody> <tr> <td>DC-195</td> <td>70 - 100%, +/-2 digits (1 S.D.)</td> </tr> <tr> <td>LNOP® Adt</td> <td></td> </tr> <tr> <td>LNOP DCI</td> <td></td> </tr> <tr> <td>LNOP DCSC</td> <td></td> </tr> <tr> <td>LNOP DC1P</td> <td></td> </tr> <tr> <td>LNOP DC150</td> <td></td> </tr> <tr> <td>LNOP Pdt</td> <td></td> </tr> <tr> <td>LNOP Y-1</td> <td></td> </tr> <tr> <td>LNOPv Ad</td> <td></td> </tr> </tbody> </table>	Sensor	Accuracy	DC-195	70 - 100%, +/-2 digits (1 S.D.)	LNOP® Adt		LNOP DCI		LNOP DCSC		LNOP DC1P		LNOP DC150		LNOP Pdt		LNOP Y-1		LNOPv Ad	
Sensor	Accuracy																				
DC-195	70 - 100%, +/-2 digits (1 S.D.)																				
LNOP® Adt																					
LNOP DCI																					
LNOP DCSC																					
LNOP DC1P																					
LNOP DC150																					
LNOP Pdt																					
LNOP Y-1																					
LNOPv Ad																					
	<table border="0"> <thead> <tr> <th>Sensor</th> <th>Accuracy</th> </tr> </thead> <tbody> <tr> <td>LNOP Neo</td> <td>70 - 100%, +/-3 digits (1 S.D.)</td> </tr> <tr> <td>LNOP Neo Pt</td> <td></td> </tr> </tbody> </table>	Sensor	Accuracy	LNOP Neo	70 - 100%, +/-3 digits (1 S.D.)	LNOP Neo Pt															
Sensor	Accuracy																				
LNOP Neo	70 - 100%, +/-3 digits (1 S.D.)																				
LNOP Neo Pt																					
	<table border="0"> <thead> <tr> <th>Sensor</th> <th>Accuracy</th> </tr> </thead> <tbody> <tr> <td>LNOP EAR</td> <td>70 - 100%, +/-3.5 digits (1 S.D.)</td> </tr> </tbody> </table>	Sensor	Accuracy	LNOP EAR	70 - 100%, +/-3.5 digits (1 S.D.)																
Sensor	Accuracy																				
LNOP EAR	70 - 100%, +/-3.5 digits (1 S.D.)																				
Measurement Wavelengths:	Red 660 Nanometers																				
	Infrared 905 Nanometers																				
Power:	Maximum radiant power at 50 mA pulsed is 0.79mW																				
Pulse Rate Range:	25 - 240 BPM																				
Pulse Rate Accuracy:	+/-3 BPM																				
Numerics:	Updated every one (1) second.																				

**NOTE:**

For further information on sensors and sensor accuracy, contact Masimo.

Characteristic	Specification																
Nelcor® OxiMax®																	
Type:	Functional Oxygen Saturation																
SpO <sub>2</sub> % Range:	1 - 100%																
SpO <sub>2</sub> Accuracy:	<table border="0"> <tr> <td style="text-align: right;"><u>Sensor</u></td> <td style="text-align: right;"><u>Accuracy</u></td> </tr> <tr> <td>MAX-A</td> <td>70 - 100%, +/-2 digits (1 S.D.)</td> </tr> <tr> <td>MAX-AL</td> <td></td> </tr> <tr> <td>MAX-N (Adult)</td> <td></td> </tr> <tr> <td>MAX-P</td> <td></td> </tr> <tr> <td>MAX-I</td> <td></td> </tr> <tr> <td>MAX-FAST</td> <td></td> </tr> <tr> <td>SC-A (Adult)</td> <td></td> </tr> </table>	<u>Sensor</u>	<u>Accuracy</u>	MAX-A	70 - 100%, +/-2 digits (1 S.D.)	MAX-AL		MAX-N (Adult)		MAX-P		MAX-I		MAX-FAST		SC-A (Adult)	
<u>Sensor</u>	<u>Accuracy</u>																
MAX-A	70 - 100%, +/-2 digits (1 S.D.)																
MAX-AL																	
MAX-N (Adult)																	
MAX-P																	
MAX-I																	
MAX-FAST																	
SC-A (Adult)																	
	<table border="0"> <tr> <td style="text-align: right;"><u>Sensor</u></td> <td style="text-align: right;"><u>Accuracy</u></td> </tr> <tr> <td>OxiCliq A</td> <td>70 - 100%, +/-2.5 digits (1 S.D.)</td> </tr> <tr> <td>OxiCliq I</td> <td></td> </tr> <tr> <td>OxiCliq N (Adult)</td> <td></td> </tr> <tr> <td>OxiCliq P</td> <td></td> </tr> </table>	<u>Sensor</u>	<u>Accuracy</u>	OxiCliq A	70 - 100%, +/-2.5 digits (1 S.D.)	OxiCliq I		OxiCliq N (Adult)		OxiCliq P							
<u>Sensor</u>	<u>Accuracy</u>																
OxiCliq A	70 - 100%, +/-2.5 digits (1 S.D.)																
OxiCliq I																	
OxiCliq N (Adult)																	
OxiCliq P																	
	<table border="0"> <tr> <td style="text-align: right;"><u>Sensor</u></td> <td style="text-align: right;"><u>Accuracy</u></td> </tr> <tr> <td>DS-100A</td> <td>70 - 100%, +/-3 digits (1 S.D.)</td> </tr> <tr> <td>D-YS (Infant to Adult)</td> <td></td> </tr> <tr> <td>MAX-N (Neonate)</td> <td></td> </tr> <tr> <td>OXI-A/N (Adult)</td> <td></td> </tr> <tr> <td>OXI-P/I</td> <td></td> </tr> <tr> <td>SC-NEO (Neonate)</td> <td></td> </tr> <tr> <td>SC-PR (Neonate)</td> <td></td> </tr> </table>	<u>Sensor</u>	<u>Accuracy</u>	DS-100A	70 - 100%, +/-3 digits (1 S.D.)	D-YS (Infant to Adult)		MAX-N (Neonate)		OXI-A/N (Adult)		OXI-P/I		SC-NEO (Neonate)		SC-PR (Neonate)	
<u>Sensor</u>	<u>Accuracy</u>																
DS-100A	70 - 100%, +/-3 digits (1 S.D.)																
D-YS (Infant to Adult)																	
MAX-N (Neonate)																	
OXI-A/N (Adult)																	
OXI-P/I																	
SC-NEO (Neonate)																	
SC-PR (Neonate)																	
	<table border="0"> <tr> <td style="text-align: right;"><u>Sensor</u></td> <td style="text-align: right;"><u>Accuracy</u></td> </tr> <tr> <td>D-YS</td> <td>70 - 100%, +/-3.5 digits (1 S.D.)</td> </tr> <tr> <td>D-YSE</td> <td></td> </tr> <tr> <td>D-YSPD</td> <td></td> </tr> <tr> <td>OxiCliq N (Neonate)</td> <td></td> </tr> </table>	<u>Sensor</u>	<u>Accuracy</u>	D-YS	70 - 100%, +/-3.5 digits (1 S.D.)	D-YSE		D-YSPD		OxiCliq N (Neonate)							
<u>Sensor</u>	<u>Accuracy</u>																
D-YS	70 - 100%, +/-3.5 digits (1 S.D.)																
D-YSE																	
D-YSPD																	
OxiCliq N (Neonate)																	
	<table border="0"> <tr> <td style="text-align: right;"><u>Sensor</u></td> <td style="text-align: right;"><u>Accuracy</u></td> </tr> <tr> <td>MAX-R</td> <td>80 - 100%, +/-3.5 digits (1 S.D.)</td> </tr> </table>	<u>Sensor</u>	<u>Accuracy</u>	MAX-R	80 - 100%, +/-3.5 digits (1 S.D.)												
<u>Sensor</u>	<u>Accuracy</u>																
MAX-R	80 - 100%, +/-3.5 digits (1 S.D.)																
	<table border="0"> <tr> <td style="text-align: right;"><u>Sensor</u></td> <td style="text-align: right;"><u>Accuracy</u></td> </tr> <tr> <td>D-YS (Neonate)</td> <td>70 - 100%, +/-4.0 digits (1 S.D.)</td> </tr> <tr> <td>OXI-A/N (Neonate)</td> <td></td> </tr> </table>	<u>Sensor</u>	<u>Accuracy</u>	D-YS (Neonate)	70 - 100%, +/-4.0 digits (1 S.D.)	OXI-A/N (Neonate)											
<u>Sensor</u>	<u>Accuracy</u>																
D-YS (Neonate)	70 - 100%, +/-4.0 digits (1 S.D.)																
OXI-A/N (Neonate)																	
Measurement Wavelengths:	Red 660 Nanometers																
	Infrared 890 Nanometers																
Power:	Not exceeding 15 mW																
Pulse Rate Range:	20 - 240 BPM																
Pulse Rate Accuracy:	+/-3 digits																
Numerics:	Updated every one (1) second.																

**NOTE:**  
For further information on sensors and sensor accuracy, contact Nellcor.

# CAS 740 Monitors

## TEMPERATURE (OPTIONAL)

Characteristic	Specification
Temperature Range:	28.9 to 42.2 °C (84.0 to 108.0 °F)
Accuracy:	+/-0.1°C (+/-0.2°F), Meets or exceeds ASTM Standards

## PATIENT ALARMS

CAS 740  
with NIBP

Patient Parameter	Neonatal Limit Range		Adult Limit Range	
	Low	High	Low	High
SYS	35 – 130	35 – 130	35 – 250	35 – 250
DIA	20 – 105	20 – 105	20 – 215	20 – 215
MAP	25 – 120	25 – 120	25 – 230	25 – 230
Pulse	45 – 235	45 – 235	35 – 235	35 – 235

CAS 740  
with NIBP  
and SpO<sub>2</sub>

Patient Parameter	Neonatal Limit Range		Adult Limit Range	
	Low	High	Low	High
SYS	35 – 130	35 – 130	35 – 250	35 – 250
DIA	20 – 105	20 – 105	20 – 215	20 – 215
MAP	25 – 120	25 – 120	25 – 230	25 – 230
%SpO <sub>2</sub>	70 – 95	80 – 99	70 – 95	80 – 99
Pulse	25 – 235	25 – 235	25 – 235	25 – 235

### NOTE:

Each alarm limit may also be selected "OFF" individually or as a whole.  
Low Limits cannot be set above the associated High Limit.  
High Limits cannot be set lower than the associated Low Limit.

## CONTROL PANEL

Characteristic	Specification
Display:	LED display of measurement results, instructions, troubleshooting messages and signal strength bar.
Parameters Displayed:	Systolic Pressure, Diastolic Pressure and Mean Arterial Pressure (MAP) Pulse Rate %SpO <sub>2</sub> Temperature (in Fahrenheit or Celsius)

## SAFETY LIMITS

Characteristic	Specification
Automatic Cuff Deflation:	If cuff pressure exceeds 290 mmHg (Adult); 145 mmHg (Neonate) If measurement time exceeds 120 seconds (Adult), 90 seconds (Neonate) If safety timer detects microprocessor failure

**OPERATING MODES**

Characteristic	Specification
Patient:	NIBP function in Neonatal or Adult
NIBP:	Manual, STAT or Automatic (at preset intervals)
History:	Review of previous measurements
%SpO <sub>2</sub> :	Continuous Monitoring
Temperature:	Predictive or Continuous Monitoring

**POWER**

Characteristic	Specification
Source:	External line or internal battery
AC Power:	100 - 240 VAC, 50/60 Hz, 0.5A; Fuse Rating – T500mAL250V (two provided)
DC Power (EMS Option):	+12 VDC; 7W; Fuse Rating – Wickman Type TE5, Time Lag, 3.15A, 125VAC or approved equivalent (one provided)
Battery:	Nickel Metal Hydride (NiMH) battery pack (user removable) Charge Time: 4 hours Operation on battery: 100 NIBP readings when set in the 5-minute Automatic Mode
Leakage Current:	100 microamp (maximum)

**FEATURES**

Characteristic	Specification
Self Test:	System self test is performed each time power is turned on.
Auto Zero:	Zero pressure reference is automatically established after every reading.
Inflation:	Initial inflation to 150 mmHg (Adult) or 85 mmHg (Neonatal) or user selectable. (100, 120, 140, 160, 180, 200) - Adult ; (60, 80, 100, 120) - Neonatal. Subsequent inflation to approximately 30 mmHg greater than previous Systolic pressure.
Deflation:	Automatic
Max Measurement Time:	Limited to 120 seconds (Adult), 90 seconds (Neonate)

# CAS 740 Monitors

## OPERATING ENVIRONMENT

Characteristic	Specification
Operating Temperature:	0°C to 50°C (32°F to 122°F)
Humidity:	15 to 95%, non-condensing
Altitude:	10,000 to -1,000 ft (690 – 1050 hPa)

Monitors may not meet performance specifications if stored or used outside temperature and humidity ranges. When moving the monitor from a storage location, wait at least one-hour prior to use to allow the monitor to adjust to room temperature.

## STORAGE/TRANSPORT ENVIRONMENT

Characteristic	Specification
Storage / Transport Temperature:	-20°C to 60°C (-4°F to 140°F)
Humidity:	15 to 95%, non-condensing
Altitude:	10,000 to -1,000 ft (690 – 1050 hPa)

## PHYSICAL DIMENSIONS & WEIGHT

Characteristic	Specification
Base Unit	
H x W x D:	6.75 in x 8.5 in x 3.0 in (17 cm x 21.5 cm x 7.5 cm)
Weight:	3 lbs approx. (1.4 kg)



**SERIAL INTERFACE**

Characteristic	Specification
Interface:	Bi-directional serial communication
Speed:	9600 – Printer 115200 – CAS Serial Protocol
Signal Level:	RS232C
Data Length:	8 bits
Start Bit:	1 bit
Stop Bit:	1 bit
Parity:	None
Flow Control:	None

**OPTIONAL ACCESSORIES**

Infrared Printer  
 Swiveled Hard Mount (for ambulance applications)  
 Roll Stand and Basket  
 Protective Carrying Case

## STANDARDS

Accuracy complies with that given in American National Standard for Electronic or Automated Sphygmomanometers, ANSI/AAMI SP10, 2002. Adult blood pressure measurements determined with this device are equivalent to those obtained by an auscultatory blood pressure measurement device and neonatal ones are equivalent to those obtained by an intra-arterial blood pressure device, within the limits prescribed by the American National Standard for Electronic or Automated Sphygmomanometers. The 4th Korotkoff sound was used to determine Diastolic pressure. Study findings are available.

Units comply with the following requirements:

- EN 60601-1
- EN 60601-1-2
- EN 865
- EN 60601-2-30
- EN 60601-2-49
- ETL Listed - UL 2601, CAN/CSA C22.2 No.601.1
- CE marking according to Directive 93/42/EEC

Masimo SET® is a registered trademark of Masimo, Inc.

Possession or purchase of this device does not convey any express or implied license to use the device with replacement parts which would, alone, or in combination with this device, fall within the scope of one or more of the patents relating to this device.


**MAXNIBP®**, Tuff-Cuff® and Pedisphyg® are registered trademarks of CAS Medical Systems, Inc.

Nellcor® and OxiMax® are registered trademarks of Nellcor Puritan Bennett, Inc.

NONIN®, NONIN® Finger Clip Sensor and Flexi-Form Sensors are registered trademarks of Nonin Medical, Inc.

Safe-Cuff™ is a trademark of CAS Medical Systems, Inc.

SureTemp® is a registered trademark of Welch Allyn, Inc.

 **CAS** is a registered trademark of CAS Medical Systems, Inc. All units covered by U.S. patent 4,796,184 and 5,022,403. Other patents pending.

Monitors are  marked.

## 16. PURCHASING RECORD

CAS 740 VITAL SIGNS MONITOR

Installed Options:

NIBP ( ) SpO<sub>2</sub> ( ) Temperature ( )

RS232 ( )

Model: \_\_\_\_\_

Serial Number: \_\_\_\_\_

Date of Purchase: \_\_\_\_\_

Dealer Name: \_\_\_\_\_

Representative: \_\_\_\_\_

Phone Number: \_\_\_\_\_

Fax Number: \_\_\_\_\_

Email: \_\_\_\_\_









## TECHNICAL BULLETIN

September 7, 2006

CAS Medical Systems, Inc. has increased the rating for the fuses set found at the Power Entry Module. Prior to this change the fuses were rated 0.5A. As a result of the change the monitor accompanying this User Manual now utilizes a 1.25A fuse set.

Please note that the prior fuse set called out replacement part number "**T500mAL250V**". The replacement fuse calls for "**T1.25AL250V**". You will find the label on the monitor has been changed, however the User Manual, until it is changed, still refers to the older rating in the following places:

- Rear case graphical depictions, numerous places throughout manual;
- Paragraph on "Fuses" section 11;
- Table of "Other Accessories", section 13;
- "Power" specifications, section 15.

Until the User Manual becomes updated please refer to the label adjacent to the power entry module for fuse rating and replacement information.

